

(19) World Intellectual Property Organization
International Bureau



(43) International Publication Date
14 February 2008 (14.02.2008)

PCT

(10) International Publication Number
WO 2008/019013 A2

(51) International Patent Classification:
A61M 5/172 (2006.01) *G06F 19/00* (2006.01)

(21) International Application Number:
PCT/US2007/017120

(22) International Filing Date: 31 July 2007 (31.07.2007)

(25) Filing Language: English

(26) Publication Language: English

(30) Priority Data:
11/499,248 3 August 2006 (03.08.2006) US

(71) Applicant (for all designated States except US): SMITHS
MEDICAL MD, INC. [US/US]; 1265 Grey Fox Road, St.
Paul, MN 55112 (US).

(72) Inventor; and

(75) Inventor/Applicant (for US only): BLOMQUIST,
Michael, L. [US/US]; 11050 Amen Circle, Blaine, MN
55449 (US).

(74) Agent: BRUESS, Steven, C.; Merchant & Gould P.C.,
P.O. Box 2903, Minneapolis, MN 55402-0903 (US).



(81) Designated States (unless otherwise indicated, for every
kind of national protection available): AE, AG, AL, AM,
AT, AU, AZ, BA, BB, BG, BH, BR, BW, BY, BZ, CA, CH,
CN, CO, CR, CU, CZ, DE, DK, DM, DO, DZ, EC, EE, EG,
ES, FI, GB, GD, GE, GH, GM, GT, HN, HR, HU, ID, IL,
IN, IS, JP, KE, KG, KM, KN, KP, KR, KZ, LA, LC, LK,
LR, LS, LT, LU, LY, MA, MD, ME, MG, MK, MN, MW,
MX, MY, MZ, NA, NG, NI, NO, NZ, OM, PG, PH, PL,
PT, RO, RS, RU, SC, SD, SE, SG, SK, SL, SM, SV, SY,
TJ, TM, TN, TR, TT, TZ, UA, UG, US, UZ, VC, VN, ZA,
ZM, ZW.

(84) Designated States (unless otherwise indicated, for every
kind of regional protection available): ARIPO (BW, GH,
GM, KE, LS, MW, MZ, NA, SD, SL, SZ, TZ, UG, ZM,
ZW), Eurasian (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM),
European (AT, BE, BG, CH, CY, CZ, DE, DK, EE, ES, FI,
FR, GB, GR, HU, IE, IS, IT, LT, LU, LV, MC, MT, NL, PL,
PT, RO, SE, SI, SK, TR), OAPI (BF, BJ, CF, CG, CI, CM,
GA, GN, GQ, GW, ML, MR, NE, SN, TD, TG).

Published:

— without international search report and to be republished
upon receipt of that report

(54) Title: INTERFACE FOR MEDICAL INFUSION PUMP

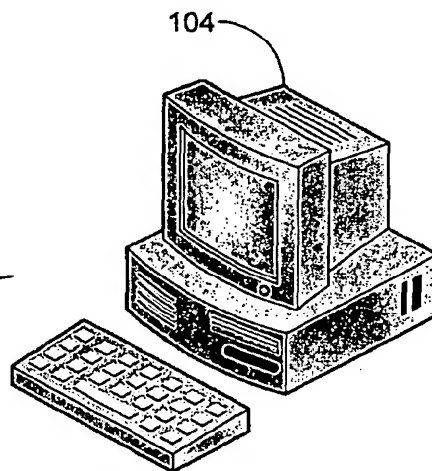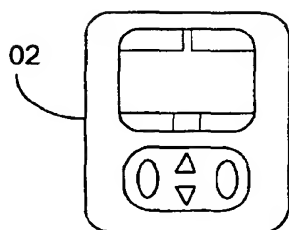

(57) Abstract: Methods of programming and use of a medical infusion pump are disclosed. One method includes displaying a first menu containing a plurality of therapy menu items. The therapy menu items each correspond to a therapy. The method also includes displaying a second menu containing at least one qualifier menu item corresponding to a qualifier. The method further includes displaying a third menu containing at least one drug menu item corresponding to a drug. The second menu is displayed upon selection of a therapy menu item, and the third menu is displayed upon selection of a qualifier menu item.

WO 2008/019013 A2

INTERFACE FOR MEDICAL INFUSION PUMP

This application is being filed on 31 July 2007, as a PCT International Patent application in the name of Smiths Medical MD, Inc., a U.S. national corporation, applicant for the designation of all countries except the US, and
5 Michael L. Blomquist, a citizen of the U.S., applicant for the designation of the US only, and claims priority to U.S. Utility Patent Application Serial No. 11/499,248, filed August 3, 2006.

Background

Patients at hospitals and other care centers regularly require
10 controlled drug intake as a part of the patient's prescribed therapy. One form of controlled drug intake is accomplished by infusing fluidic drugs with a medical infusion pump.

Medical infusion pumps, in general, provide regulated drug delivery to a patient. These pumps are used to deliver a selected drug or other therapeutic
15 agent to a patient at a predetermined rate that is programmed into the pump. However, programming and managing such pumps can be difficult and cumbersome. Programming typically includes preloading a pump program into a pump and then entering pump parameters or data into the pump through a keypad that is directly in the pump. Each time the pump is programmed, the data must be
20 reentered by hand.

Managing the status and locations of pumps also can be difficult. A single pump can be us programmed for delivering different fluids in different therapies and in different locations within a hospital. Similarly, the status of a pump and alarms can be difficult to monitor because the pumps are often in locations other
25 than where the caregiver is located and have small displays on which information can be difficult to see.

Summary

A method of programming a medical infusion pump with a protocol is disclosed. The method includes displaying a first menu containing a plurality of
30 therapy menu items, where each therapy menu item corresponds to a therapy. The method further includes displaying a second menu containing at least one qualifier menu item, where each qualifier menu item corresponds to a qualifier. The second

menu displays upon selecting a therapy menu item from the first menu. The method also includes displaying a third menu containing at least one drug menu item, where each drug menu item corresponds to a drug. The third menu displays upon selection of a qualifier menu item from the second menu.

5 According to a second aspect, an apparatus for programming a medical infusion pump with a protocol is disclosed. The apparatus includes a memory configured to store a plurality of therapy menu items, qualifier menu items, and drug menu items. The apparatus further includes a programmable circuit in electrical communication with the memory. The programmable circuit is
10 programmed to generate a user interface configured to display a first menu containing a plurality of therapy menu items, where each therapy menu item corresponds to a therapy. The programmable circuit is further programmed to display a second menu containing at least one qualifier menu item, where each qualifier menu item corresponds to a qualifier. The programmable circuit displays
15 the second menu upon selection of a therapy menu item from the first menu. The programmable circuit is further programmed to display a third menu containing at least one drug menu item, each drug menu item corresponding to a drug. The programmable circuit displays the third menu upon selection of a qualifier menu item from the second menu.

20 According to a third aspect, an apparatus for programming a medical infusion pump with a protocol is disclosed. The apparatus includes a memory configured to store a plurality of therapy menu items, qualifier menu items, and drug menu items. The apparatus also includes a display. The apparatus further includes a programmable circuit in electrical communication with the memory and the display.
25 The programmable circuit is programmed to generate a user interface. The user interface displays a first menu containing a plurality of therapy menu items, where each therapy menu item corresponds to a therapy administrable by the medical infusion pump. The user interface also displays, upon selecting a therapy menu item from the first menu, a second menu containing at least one qualifier menu item,
30 where each qualifier menu item corresponds to a qualifier associated with the therapy. The user interface further displays, upon selecting a qualifier menu item from the second menu, a third menu containing at least one drug menu item, where each drug menu item corresponds to a drug associated with the therapy and the

qualifier. The programmable circuit is programmed to, upon selecting a drug menu item from the third menu, identify a set of pump data corresponding to the selected therapy menu item, qualifier menu item, and drug menu item. The programmable circuit is further programmed to select a program corresponding to the set of pump data. The programmable circuit is also programmed to upload the program to the medical infusion pump.

According to a further aspect, an apparatus for programming a medical infusion pump with a protocol is disclosed. The apparatus includes a memory configured to store a plurality of therapy menu items, qualifier menu items, and drug menu items. The apparatus further includes a programmable circuit in electrical communication with the memory and programmed to generate a graphical user interface. The graphical user interface has therapy, qualifier, and drug window objects. The graphical user interface is configured to present a list of therapy menu items in the therapy window object, where each therapy menu item corresponds to a therapy. The graphical user interface is also configured to, upon selecting a therapy menu item in the therapy window object, present a list of qualifier menu items in the qualifier window object, where each qualifier menu item corresponds to a qualifier. The graphical user interface is also configured to, upon selecting a qualifier menu item in the qualifier window object, present a list of drug menu items in the drug window object, where each drug menu item corresponds to a drug.

According to yet another aspect, a method of programming a medical infusion pump with a protocol is disclosed. The method includes presenting on a monitor a list of therapy menu items in a therapy window object, where each therapy menu item corresponds to a therapy. The method also includes, upon selecting a therapy menu item in the therapy window object, presenting on the monitor a list of qualifier menu items in a qualifier window object, where each qualifier menu item corresponds to a qualifier. The method further includes, upon selecting a qualifier menu item in a qualifier window object, presenting on the monitor a list of drug menu items in the drug window object, where each drug menu item corresponds to a drug.

According to a further aspect, a method of programming a medical infusion pump with a protocol is disclosed. The method includes displaying a first menu containing a plurality of first index items. The method also includes, upon

selecting a first index item, displaying a second menu containing at least one second index item. The method further includes, upon selecting a second index item from the second menu, displaying a third menu containing at least one third index item. The first index item, the second index item, and the third index item in combination
5 refer to a protocol for programming into the medical infusion pump.

According to another aspect, an apparatus for programming a medical infusion pump with a protocol is disclosed. The apparatus includes a memory configured to store a plurality of first menu items, second menu items, and third menu items. The apparatus also includes a display. The apparatus further includes a
10 programmable circuit in electrical communication with the memory and the display. The programmable circuit is programmed to display a first menu containing a plurality of first index items. The programmable circuit is further programmed to, upon selecting a first index item, display a second menu containing at least one second index item. The programmable circuit is also programmed to, upon selecting
15 a second index item from the second menu, display a third menu containing at least one third index item. The first index item, the second index item, and the third index item in combination refer to a protocol for programming into the medical infusion pump.

Brief Description of the Drawings

20 FIG. 1 illustrates a pump-computer communication system according to a possible embodiment of the present disclosure;

FIG. 2 illustrates an infusion pump network according to a possible embodiment of the present disclosure;

25 FIG. 3 illustrates the architecture of a computing system that can be used to implement aspects of the present disclosure;

FIG. 4 illustrates the architecture of a pump that can be used to implement aspects of the present disclosure;

FIG. 5 is an exemplary infusion pump network according to a possible embodiment of the present disclosure;

30 FIG. 6A is an exemplary data structure for a pump protocol library according to a possible embodiment of the present disclosure;

FIG. 6B is an exemplary data structure for a pump protocol library according to a possible embodiment of the present disclosure;

FIG. 6C is an exemplary data structure for pump protocols according to a possible embodiment of the present disclosure;

5 FIG. 7 is an exemplary architecture of administrative software for setting global pump protocols according to a possible embodiment of the present disclosure;

FIG. 8 is one example of a computer user interface library import screen in accordance with the present disclosure;

10 FIG. 9 is one example of a computer user interface for administrative software in accordance with the present disclosure;

FIG. 10 is one example of a computer user interface location tab in accordance with the present disclosure;

15 FIG. 11 is one example of a computer user interface therapy tab in accordance with the present disclosure;

FIG. 12 is one example of a computer user interface protocol tab in accordance with the present disclosure;

FIG. 13 is one example of a computer user interface drug bar code display screen in accordance with the present disclosure;

20 FIG. 14 is one example of a computer user interface prescription order form display screen in accordance with the present disclosure;

FIG. 15 is one example of a computer user interface therapy selection screen in accordance with the present disclosure;

25 FIG. 16 is one example of a computer user interface qualifier selection screen in accordance with the present disclosure;

FIG. 17 is one example of a computer user interface drug selection screen in accordance with the present disclosure;

FIG. 18 is one example of a computer user interface drug delivery tab in accordance with the present disclosure;

30 FIG. 19 is one example of a computer user interface weight-based drug delivery tab in accordance with the present disclosure;

FIG. 20 is one example of a computer user interface secondary drug delivery tab in accordance with the present disclosure;

FIG. 21 is one example of a computer user interface alarm tab in accordance with the present disclosure;

FIG. 22 is one example of a computer user interface security tab in accordance with the present disclosure;

5 FIG. 23 is one example of a computer user interface appearance tab in accordance with the present disclosure;

FIG. 24 is one example of a computer user interface report tab in accordance with the present disclosure;

10 FIG. 25 is one example of a computer user interface library export screen in accordance with the present disclosure;

FIG. 26 is a flow diagram of methods and systems for custom programming of a medical infusion pump according to a possible embodiment of the present disclosure;

15 FIG. 27 is one example of a computer user interface library import screen in accordance with the present disclosure;

FIG. 28 is a flow diagram of methods and systems for editing and loading a protocol for a medical infusion pump according to a possible embodiment of the present disclosure;

20 FIG. 29 is one example of a computer user interface protocol selection screen in accordance with the present disclosure;

FIG. 30 is one example of a computer user interface therapy selection screen in accordance with the present disclosure;

FIG. 31 is one example of a computer user interface qualifier selection screen in accordance with the present disclosure;

25 FIG. 32 is one example of a computer user interface drug selection screen in accordance with the present disclosure;

FIG. 33 is a flow diagram of methods and systems for custom programming of a medical infusion pump according to a possible embodiment of the present disclosure;

30 FIG. 34 is one example of a computer user interface drug delivery customization screen in accordance with the present disclosure;

FIG. 35 is one example of a computer user interface drug delivery customization screen in accordance with the present disclosure;

FIG. 36 is one example of a computer user interface medical infusion pump programming screen in accordance with the present disclosure;

FIG. 37 is a flow diagram of methods and systems for displaying medical infusion pump customizations according to a possible embodiment of the present disclosure;

FIG. 38 is one example of a computer user interface pump comparison screen in accordance with the present disclosure;

FIG. 39 is a schematic front view of a medical infusion pump displaying a change bar according to a possible embodiment of the present disclosure;

FIG. 40 is one example of a computer user interface report generation screen in accordance with the present disclosure; and

FIG. 41 is one example of a computer user interface report screen in accordance with the present disclosure.

Detailed Description

Various embodiments will be described in detail with reference to the drawings, wherein like reference numerals represent like parts and assemblies throughout the several views. Reference to various embodiments does not limit the scope of the claims attached hereto. Additionally, any examples set forth in this specification are not intended to be limiting and merely set forth some of the many possible embodiments for the appended claims.

The following discussion is intended to provide a brief, general description of a suitable computing environment in which the invention may be implemented. Although not required, the invention will be described in the general context of computer-executable instructions being executed by a computer, for example, a hand held computer, a personal computing system, or a medical infusion pump. The structure, creation, and use of a message store hierarchical folder structure are described after the discussion of an exemplary operating environment.

Additionally, the logical operations of the various embodiments of the invention described herein are implemented as: (1) a sequence of computer implemented operations running on a computing system; and/or (2) interconnected machine modules within the computing system. Modules represent functions

executed by program code such as commonly available programming languages or as the code found in a dynamic-link library (DLL). The implementation used is a matter of choice dependent on the performance requirements of the pump and the computing systems with which it interfaces. Accordingly, the logical operations
5 making up the embodiments of the invention described herein can be referred to alternatively as operations, modules, and the like.

FIG. 1 illustrates an exemplary embodiment of an infusion pump network 100 having a medical infusion pump 102, a computing system 104, and a communications link 106. The medical infusion pump 102 is configured to deliver
10 therapeutic fluids, such as drugs, saline, or nutrition to a patient. Examples, of medical infusion pumps 102 include ambulatory pumps, stationary pumps, and pole mounted pumps.

The computing system 104 is configured to execute computer-readable instructions, such as computer software. The computing system 104 can be
15 located in a variety of locations such as the point of care (POC) where a patient is being treated, in a healthcare facility at a location remote from the POC, or even at an off-site location remote from the healthcare facility itself. In further embodiments, the medical infusion pump 102 acts as the computing system 104.

In the exemplary embodiment, the computing system 104 is
20 programmed to generate and store pump protocols for execution in the context of a pump application program. Each pump protocol includes a series of pump parameters. Pump parameters refer to settings that define an operational aspect of a medical infusion pump. The pump parameters dictate the control of the pump.

Pump protocols are collections of these pump parameters defining the
25 variable operational characteristics of a medical infusion pump during application of a specific therapy, qualifier, and drug. The pump protocol includes a listing of operational parameters to be included in the pump, and correlates to an index for referring to a specific protocol containing a specific set of pump parameters. The index can be associated with a therapy, qualifier, and drug, and is either contained
30 within the protocol or associated with a specific protocol. The pump protocol includes patient specific pump parameters and non-patient specific pump parameters. Patient specific pump parameters refer to those parameters which are set on a patient-by-patient basis, and for example include the basal delivery rate or

bolus amount. Non-patient specific pump parameters refer to those parameters which are set for the pump to perform specific tasks, and do not account for the specific patient to which they are applied. These parameters are generally related to the pump, the infusion pump network, or the medical care to be provided by the pump and/or pump network. Non-patient specific pump parameters can include, for example, a range of permissible values for basal delivery, a range of values and patterns for basal delivery, a range of permissible values for boluses, a range of values and patterns for extended boluses, a starting value within a particular range of values, alarm values, protocols for data communication, and various flag settings.

A pump application program is a program having instructions (e.g., executable code, rules, and/or data) that control operation of the pump for a specific therapy or type of delivery (e.g., continuous delivery, intermittent delivery, pain control, chemotherapy, total parenteral nutrition, etc.). For example, a pump application program might contain instructions that define operation of a pump to accomplish various of the pump parameters. Pump application programs include, for example, pump protocols including both patient specific and non-patient specific pump parameters, and instructions for allocating memory, user interfaces, or algorithms for monitoring various sensors and driving a motor for the pump mechanism.

The communications link 106 connects the pump 102 and computing system 104. In various embodiments, the communications link 106 can include serial or parallel connections, wired or wireless connections, and a direct or networked connection to a computer. Additionally, the pump 102 and the computing system 104 can communicate using any protocol appropriate for data communication. Examples of network connections to a computer include Intranet, Internet, and LAN (e.g., Ethernet). Examples of wired connections to a computer include USB, RS-232, Firewire, and power-line modem connection. Examples of wireless connections include bluetooth, 802.11a/b/g, infrared (IR), and radio frequency (RF).

FIG. 2 illustrates an exemplary embodiment of an infusion pump network 200 having a server 206 networked with a plurality of computing systems 104₁-104_n. The network 200 can be any wired or wireless network that enables data communication between the server, computing systems, and medical infusion

pumps. Examples of networks include the Internet, Intranets, and LANs. Each computing system 104 can communicate with a medical infusion pump 102₁-102_n through a communication link 106.

In the exemplary embodiment, the individual computing systems
5 104_n execute software for generating and managing pump application programs and sets of pump operating parameters. The pump application programs and sets of pump operating parameters are stored on the server 206 so they can be accessed by other individual computing systems 104_n. The individual computing systems 104_n are also programmed to retrieve previously created pump application programs and
10 sets of pump operating parameters that are stored on the server 206 for viewing, editing, and downloading to medical infusion pumps 102_n.

In alternative embodiments, the medical infusion pumps 102_n can directly access the server to retrieve pump application programs and sets of pump operating parameters. For example, the medical infusion pumps 102_n can be loaded
15 with client software such as a web browser and communicate directly with the network 200, either through a wired or wireless connection as described herein.

In other alternative embodiments, one or more of the computing systems is not configured to communicate directly with a medical infusion pump 102_n, but rather provides administrative access to the server 206 for generating,
20 viewing, and editing pump application programs and sets of pump operating parameters. Additionally, servers, workstations, and other computing systems unaffiliated with the medical infusion pumps 102_n can be included in the network 200.

In yet other alternative embodiments, the software is executed in the
25 server 206. For example, the server functions as an application service provider that communicates user interface and other data entries in mark-up language such as HTML or some other language or protocol that allows a user to execute software from a remote location. In these embodiments, the server 206 can function as an application service provider in which the server provides access to the software for
30 generating and storing pump application programs and pump protocols that a user can create and download to a medical infusion pump. For example, the server 206 could be located at a pump manufacture, pharmaceutical manufacture, pharmacist, or some other third party separate from the user. The server 206 in such an

embodiment can be accessed either from an individual computing system 104 or by a medical infusion pump 102 that has networking capabilities and client software.

Example embodiments of a server 206 and a medical infusion pump 102 having a web browser are disclosed in United States patent application serial no. 11/066,425, which was filed on February 22, 2005 and is entitled Server for Medical Device, the entire disclosure of which is hereby incorporated by reference.

FIG. 3 illustrates an exemplary architecture that can be used to implement aspects of the present disclosure, including the computing systems 104 and the server 206. The computing system architecture includes a general purpose computing device in the form of a computing system 300. The computing system 300 can be used, for example, as the computing system or server of FIG. 2, and can execute program modules included in the administrative software or user software disclosed below.

The computing system 300 including at least one processing system 302. A variety of processing units are available from a variety of manufacturers, for example, Intel or Advanced Micro Devices. The computing system 300 also includes a system memory 304, and a system bus 306 that couples various system components including the system memory 304 to the processing unit 302. The system bus 306 may be any of a number of types of bus structures including a memory bus, or memory controller; a peripheral bus; and a local bus using any of a variety of bus architectures.

The system memory 304 can include read only memory (ROM) 308 and random access memory (RAM) 310. A basic input/output system 312 (BIOS), containing the basic routines that help transfer information between elements within the computing system 300, such as during start up, is typically stored in the ROM 308.

The computing system 300 can also include a secondary storage device 313, such as a hard disk drive, for reading from and writing to a hard disk (not shown), and/or a compact flash card 314.

The hard disk drive 313 and compact flash card 314 are connected to the system bus 306 by a hard disk drive interface 320 and a compact flash card interface 322, respectively. The drives and cards and their associated computer

readable media provide nonvolatile storage of computer readable instructions, data structures, program modules and other data for the computing system 300.

Although the exemplary environment described herein employs a hard disk drive 313 and a compact flash card 314, other types of computer-readable media, capable of storing data, can be used in the exemplary system. Examples of
5 these other types of computer-readable mediums include magnetic cassettes, flash memory cards, digital video disks, Bernoulli cartridges, CD ROMs, DVD ROMs, random access memories (RAMs), or read only memories (ROMs).

A number of program modules may be stored on the hard disk 313,
10 compact flash card 314, ROM 308, or RAM 310, including an operating system 326, one or more application programs 328, other program modules 330, and program data 332. A user may enter commands and information into the computing system 300 through an input device 334. Examples of input devices might include a keyboard, mouse, microphone, joystick, game pad, satellite dish, scanner, digital
15 camera, touch screen, and a telephone. These and other input devices are often connected to the processing unit 302 through an interface 340 that is coupled to the system bus 306. These input devices also might be connected by any number of interfaces, such as a parallel port, serial port, game port, or a universal serial bus (USB). Wireless communication between input devices and interfaces 340 is
20 possible as well, and can include infrared, bluetooth, 802.11a/b/g, cellular, or other radio frequency communication systems. A display device 342, such as a monitor or touch screen LCD panel, is also connected to the system bus 306 via an interface, such as a video adapter 344. The display device 342 might be internal or external. In addition to the display device 342, computing systems, in general, typically
25 include other peripheral devices (not shown), such as speakers, printers, and palm devices.

When used in a LAN networking environment, the computing system 300 is connected to the local network through a network interface or adapter 352. When used in a WAN networking environment, such as the Internet, the computing
30 system 300 typically includes a modem 354 or other communications type, such as a direct connection, for establishing communications over the wide area network. The modem 354, which can be internal or external, is connected to the system bus 306 via the interface 340. In a networked environment, program modules depicted

relative to the computing system 300, or portions thereof, may be stored in a remote memory storage device. It will be appreciated that the network connections shown are exemplary and other methods of establishing a communications link between the computing systems may be used.

5 The computing system 300 might also include a recorder 360 connected to the memory 304. The recorder 360 includes a microphone for receiving sound input and is in communication with the memory 304 for buffering and storing the sound input. The recorder 360 also can include a record button 361 for activating the microphone and communicating the sound input to the memory
10 304.

- A computing device, such as computing system 300, typically includes at least some form of computer-readable media. Computer readable media can be any available media that can be accessed by the computing system 300. By way of example, and not limitation, computer-readable media might comprise
15 computer storage media and communication media.

 Computer storage media includes volatile and nonvolatile, removable and non-removable media implemented in any method or technology for storage of information such as computer readable instructions, data structures, program modules or other data. Computer storage media includes, but is not limited to,
20 RAM, ROM, EEPROM, flash memory or other memory technology, CD-ROM, digital versatile disks (DVD) or other optical storage, magnetic cassettes, magnetic tape, magnetic disk storage or other magnetic storage devices, or any other medium that can be used to store the desired information and that can be accessed by the computing system 300.

25 Communication media typically embodies computer-readable instructions, data structures, program modules or other data in a modulated data signal such as a carrier wave or other transport mechanism and includes any information delivery media. The term "modulated data signal" refers to a signal that has one or more of its characteristics set or changed in such a manner as to encode
30 information in the signal. By way of example, and not limitation, communication media includes wired media such as a wired network or direct-wired connection, and wireless media such as acoustic, RF, infrared, and other wireless media. Combinations of any of the above should also be included within the scope of

computer-readable media. Computer-readable media may also be referred to as computer program product.

FIG. 4 illustrates the architecture of a medical infusion pump 400 that can be used to implement aspects of the present disclosure. A microprocessor 402 is
5 in electrical communication with and controls a pump motor 404, a screen 406, an audible alarm 408, and a vibratory alarm 410. Other embodiments can use a microcomputer, or any other type of programmable circuit, in place of the microprocessor.

The pump motor 404 drives a drive mechanism 412. The drive
10 mechanism 412 delivers the therapeutic fluid to a patient. The drive mechanism can be connected to a plunger system, a peristaltic drive mechanism, or another type of fluid delivery system.

The screen 406 can have many different configurations such as an LCD screen. The screen 406 displays a user interface that presents various items of
15 information useful to a patient or caregiver. The audible alarm 408 is a beeper, and an alarm provides actual alarms, warnings, and reminders. Similar to other portable electronic devices such as a cellular telephone, the vibratory alarm 410 provides an alarm to either supplement the audio alarms or replace the audio alarm when an audible beep would be disruptive or not heard. A user can selectively enable or
20 disable the audible 408 and vibratory 410 alarms. In one possible embodiment, however, both the audible 408 and vibratory 410 alarms cannot be disabled at the same time.

The microprocessor 402 is in electrical communication with both a random access memory (RAM) 416 and a read only memory (ROM) 418, which are
25 onboard the pump 400 but external to the microprocessor 402 itself. In one possible embodiment, the microprocessor 402 includes internal memory as well. The RAM 416 is a static RAM stores that data that can change over time such as pump settings and a historical log of events experienced by the medical infusion pump 400. The ROM 418 stores code for the operating system and the application programs. The
30 ROM 418 can be any type of programmable ROM such as an EPROM. In one possible embodiment, the RAM 416 has 500 kilobytes of memory capacity and the ROM 418 has 2 megabytes of memory capacity.

An infrared (IR) port 420 is in electrical communication with the microprocessor. As explained in more detail below, the IR port 420 provides data communication with an external device such as a computer for programming an application program, programming pump settings, and downloading historical data logs. The medical infusion pump 400 can include other types of communication ports in place of or in addition to the IR port 420. Examples of other possible communication ports include a radio frequency (RF) port or a port that provides a hard-wired data communication link such as an RS-232 port, a USB port, or the like.

A real-time clock 422 provides a clock signal to the microprocessor 402. An advantage of having a real-time clock 422 is that it provides the program with the actual time in real-time so that the programs executed by the medical infusion pump can track and control the actual time of day that drug delivery and other events occur. Various durations described here are used for alerts, alarms, reminders, and other functions. In one possible embodiment, the timers are formed by the real-time clock 422 and software executed by the microprocessor 402.

A keypad 424 also provides input to the microprocessor 402. Although other possible types of keypads are possible, one type of keypad has four buttons and is a membrane-type of keypad, which provides resistance to water and other environmental conditions. The keypad 424 contains soft keys for which the function of the keys can change as a user executes different menu selections and commands.

An audio bolus button 425 optionally provides input to the microprocessor 402. The audio bolus button 425 can program the pump 400 to audibly administer a bolus of drugs or other therapeutic fluids without requiring visual confirmation using the pump. In an example embodiment, the audio bolus button 425 can be pressed a series of times to trigger bolus delivery of a selected volume, based on a preprogrammed trigger granularity. A single button press can represent a bolus of 5 grams, as selected by a user, and subsequent presses of the audio bolus button can represent multiples thereof.

Other inputs into the microprocessor 402 can include an occlusion sensor 426, which is sensitive to occlusions in the therapeutic fluid delivery line; a cartridge sensor 428, which is sensitive to the presence of a therapeutic fluid cartridge; and a motion detector 430, which detects motion of a gear (not shown) in

the drive mechanism 412. In an exemplary embodiment, the cartridge sensor 428 includes one or more sensors configured to detect insertion of a therapeutic fluid cartridge. The pump 400 can detect the type of cartridge present via a mechanical interface, and can include in the pump software instructions regarding operation in
5 conjunction with the cartridge. Examples of cassette sensing features are described, for example, in United States Patent No. 5,531,697, filed on April 15, 1994, issued on July 2, 1996, and entitled Systems and Methods for Cassette Identification for Drug Pumps.

FIG. 5 illustrates a schematic architecture of a medical infusion pump
10 network 500 according to an exemplary embodiment. The medical infusion pump network 500 includes an administrator computer 502 communicatively connected to the server 206 of FIG. 2, which includes a database 504. The medical infusion pump network 500 also includes one or more medical infusion pumps 102 and computing systems 104.

15 The administrator computer 502 and computing systems 104 are systems such as those described above in conjunction with FIG. 3. The administrator computer 502 includes administrative software installed on or accessible to the computer for generating one or more libraries 508 of pump protocols 510. An exemplary embodiment of the administrative software is
20 described below in FIGS. 7-25.

In the present disclosure, libraries refer to collections of pump protocols generated using the administrative software described herein. Libraries can be stored in files, databases, or other data structures. Libraries contain pump protocols as well as indices pointing to the protocols, and are loaded in user software
25 to select a specific pump protocol for operation of a medical infusion pump.

The computing systems 104 include user software for accessing one or more libraries 508 of protocols 510 and programming a medical infusion pump 102 with a protocol 510 or a library 508. In one possible embodiment, the computing systems 104 are optional in that the user software resides directly on the
30 medical infusion pumps 102. An exemplary embodiment of the user software is described below in FIGS. 26-41.

The medical infusion pumps 102 connect either to a computing system 104 or directly to the server 206, and are described above in conjunction with

FIG. 4. In a first embodiment, the medical infusion pumps 102 are configured to accept a pump protocol from the server 206 or the computing system 104. In a second embodiment, the medical infusion pumps 102 are configured to accept a library 508 of pump protocols 510 directly from the server 206 or from the computing system 104.

The database 504 contains pump protocol data 506 and log files 516. The pump protocol data 506 forms a plurality of libraries 508 which in turn each include a number of protocols 510. Each protocol 510 is stored as a data record, and includes a set of parameters, including patient specific pump parameters 512a and non-patient specific pump parameters 512b, as described above. Each library 508 can contain one or more pump protocols 510.

The log files 516 include log data regarding access and usage of the libraries 508, and can include information related to the administrator computer 502, the medical infusion pumps 102, or the computing systems 104 authorized to connect to the server 206. In one possible embodiment, the log files include access records, which record instances in which medical infusion pumps access a library 508 on the server 206.

FIGS. 6A-6C illustrate exemplary data sets accepted by medical infusion pumps 102. The specific data set accepted by a medical infusion pump 102 is dependent upon that pump, but in various embodiments, the pumps accept data sets representing pump libraries, pump protocols, or pump programs. FIG. 6A shows a library 508 in greater detail. The library 508 can be loaded into the memory of a medical infusion pump 102, allowing a user of the pump to select a protocol for operation. The library 508 includes a number of protocols 510, which include patient specific pump parameters 512a, non-patient specific pump parameters 512b, and an index 514. The number of protocols 510 within a given library can vary, and depends upon the number defined using the administrative software.

The total number of pump parameters 512a-512b remains constant for each particular model of pump, but can vary between types of pumps. Additionally, the pump parameters 512a-512b can be configured in a number of formats within each protocol 510 for the same type of pump. For example, the number of patient specific pump parameters 512a can vary between protocols due to the specific type of drug and therapy applied. For example, a protocol defining a

continuous drug delivery may only require a single patient specific protocol, namely, the drug delivery rate. In another example, a protocol defining intermittent drug delivery may require additional patient specific pump parameters, such as the time between drug delivery phases, a bolus amount, patient bolus amounts, and other parameters. The number of non-patient specific pump parameters 512b represents the difference between the total number of parameters programmable into a pump and the number of patient specific pump parameters 512a as dictated by the therapy and drug applied.

The index 514 can be any generic index referencing a specific location within the library. Each index is unique within the library, although another library may contain the same index and relate that index to a different set of pump parameters contained within that library. In the exemplary embodiment shown, the index 514 includes therapy, qualifier, and drug regions. By selecting a combination of a therapy, a qualifier, and a drug, a user of the system can select one of the protocols 510 from the library 508. Therapies, as referred to herein, are the methods of patient treatment for diseases or generalized rehabilitation. For example, a therapy can be an epidural treatment or patient-controlled analgesia. Qualifiers include factors affecting the administration of a therapy, such as weight, age, or sensitivity of a patient to a specific therapy. Drugs refer to any therapeutic fluids deliverable by a medical infusion pump.

Each of the protocol entries on the server can be assigned an identification code in order to ensure that the medical infusion pumps access a correct library and/or protocol, and that the protocols on the medical infusion pumps and computing systems associated with the server 206 of FIG. 5 are up to date. These identification codes are generated by the server and stored in conjunction with the protocol in the server 206, as well as in the medical infusion pump 102 and/or its associated computing system 104. The identification codes can be generated using globally unique identifiers (GUID), and are used to track the specific protocol and/or library accessed by each pump 102 or computing system 104 in the database 504. A GUID is a 128-bit pseudo random number used to provide a statistically unique identifier for corresponding the protocol on the pump 102 to the protocol as stored in the server 206. The GUID can be generated by the server 206 and transmitted alongside the protocol and/or library when transmitted to the computing system 104

or infusion pump 102. The server copy of the GUID can be stored alongside the protocol in the database 504, and the pump copy can be stored in RAM in the medical infusion pump 102 or computing system 104. Each time the medical infusion pump 102 or computing system 104 accesses the protocol on the server
5 206, the GUID assigned to that protocol as stored in the pump 102 or computing system 104 can be matched to the protocol as stored in the database 504 to ensure that the correct protocol is accessed. Different protocols in different libraries can use the same index of therapy, qualifier, and drug, and can be stored in the same database, but have a unique GUID and are therefore uniquely identifiable. In a
10 possible embodiment, the GUID system can be used in conjunction with user access control systems, such as are disclosed in conjunction with both the user software and administrative software, below. In a further embodiment, a new GUID can be generated and associated with each protocol when first created, or optionally each time the protocol is edited using administrative software, such as is described below.

15 FIG. 6B illustrates a second possible embodiment of a library 508. The library 508 contains a number of protocols 510, which, as in FIG. 6A, include patient specific pump parameters 512a, non-patient specific pump parameters 512b, and an index 514. The number of protocols 510 included in each library 508 may vary, so the number of protocols 510 in the library 508 of FIG. 6B may be different
20 from the number of protocols in the library shown in 6A. In the embodiment shown, the library 508 contains one protocol 510. Additionally, one or more of the therapy, qualifier, and drug regions in the index 514 can be replaced by other index criteria, such as locations, pump programs, doctor identification, or other indexable criteria capable of referring to a unique pump protocol within the library. In the example
25 shown, the "therapy" region is used to select a pump program, such as a continuous delivery program, an intermittent delivery program, or a specific type of program such as a pain management program. The "qualifier" region is used to select a doctor, and may be the name of a doctor using the infusion pump network of FIG. 5.

FIG. 6C illustrates a series of possible pump protocols 510a-510c.
30 The pump protocols 510a-510c incorporate patient specific pump parameters 512a and non-patient specific pump parameters 512b. Protocols 510 are specific to a pump 102, in that the pump has a specific number and type of parameters that are programmable. Therefore, the total number of pump parameters remains constant.

However, the number of patient specific pump parameters 512a can vary depending upon the protocol selected for programming into the pump, which in turn dictates that the number of non-patient specific pump parameters 512b varies as well. In a possible embodiment, one of the protocols 510 is selected using a computing system
5 104 or infusion pump 102. If selected using the computing system, the protocol is then programmed into the medical infusion pump 102.

In another possible embodiment, the pump protocol 510 is selected using the infusion pump 102 or the computing system 104. The pump protocol 510 is then incorporated into a pump program to provide a set of instructions dictating
10 the operation of a medical infusion pump 102 according to the protocol 510. The complete pump program is then downloaded into the pump 102. In yet another possible embodiment, the pump program is downloaded to the pump 102 at a different time from the pump protocol 510. In still a further embodiment, multiple pump programs reside within the pump 102, and the pump protocol 510 contains a
15 parameter which dictates which pump program is to be used. In a further embodiment, the pump program within the pump is altered based on one or more of the pump parameters included in the pump protocol 510.

FIG. 7 illustrates exemplary architecture of administrative software 700 for generating one or more libraries of pump protocols. The software 700 can
20 operate within the server 206, pump 102, computing system 104, or a combination thereof.

The administrative software 700 allows a user, for example a doctor, nurse, pharmacist, or other caregiver, to create, define, and edit pump application programs and protocols for execution in and control of medical infusion pumps 102.
25 For example, the administrative software 700 can generate protocols and programs that can be loaded using the user software described in FIGS. 26-41, below.

The administrative software 700 provides protocol-based programming of medical infusion pumps in which the user creates a pump application program by designating a particular therapy and other criteria such as a
30 location and qualifiers (e.g., patient age, weight, skin surface area). Once criteria are selected, the administrative software 700 applies rules and other logic that assembles sets of pump parameters into a pump protocol. For example, the administrative software 700 might be used to select one delivery pattern and enable bolus delivery

if the selected therapy is for delivering pain medication and another delivery pattern and not enable bolus delivery if the selected therapy is for parenteral nutrition. In another example, the administrative software 700 might be used to select one range of permissible delivery rates if one of the criteria indicates the patient is an adolescent and different range of permissible delivery rates if the patient is an adult. Other embodiments permit programming a medical infusion pump 102 without using therapy-based programming. Additional embodiments of protocol- or therapy-based programming is discussed in more detail in United States patent application serial no. 11/003,147, filed on December 3, 2004 and entitled
5
10 Programming Medical Pumps with Electronic Standing Order Template, the entire disclosure of which is hereby incorporated by reference.

Operation of the software 700 begins at a start module 702. The start module 702 corresponds to initial execution of the administrative software by clicking on an icon on the computer or by some other mechanism for executing
15 software. Upon startup, the software 700 connects to a library loaded in the database 504 of FIG. 5.

Following the start module, operational flow optionally proceeds to a load library module 704, which allows a user to access a listing of library files available to the administrative software 700. The library files contain one or more
20 libraries, which in turn contain a collection of pump protocols as described above. The collection of library files can be stored in the server 206 or in one or more individual computing systems 102. The load library module 704 allows the user to select a library file containing one or more libraries for viewing, editing, and downloading to a medical infusion pump 102. If a user does not want to download
25 or otherwise access an existing library, it can selectively bypass the load library module 704. An example of when a user bypasses the load library module 704 is if the user plans to only create a new library or edit one or more protocols within the currently loaded library. In an alternative embodiment, the software always executes the load library module 704 and the user then selectively chooses whether
30 to load any previously created libraries via a stored library file.

Following the start module 702 and optional load library module 704, operational flow proceeds to a login module 706. The login module 706 regulates user rights in the software 700. User rights define access levels to the currently

connected library in user software, and are configurable for users such as doctors, nurses, or other caregivers. Based on the user rights assigned to a caregiver, that user will have a set access level allowing the user to view, add, or edit pump libraries within the user software, described in detail below. Access levels can be
5 set according to a variety of criteria. Examples include the type of caregiver (e.g., physician, nurse, pharmacist), location (e.g., hospital, clinic, pharmacy, manufacturer), or a particular department within a location.

In possible embodiments, different access levels also provide different rights with respect to a particular pump protocol or pump operational
10 parameters. For example, one access level might give a user a right to edit, create, and download pump protocols and/or pump application programs. One access level might permit a user the right to edit, create, and download only specified pump parameters, such as the patient specific pump parameters described in conjunction with FIGS. 5-6. One access level might permit a user the right to only edit or
15 download pump parameters. One access level might permit a user the right to only view and download pump parameters. Different embodiments can include the ability to provide an access level for a user any combination of rights to create, edit, view, and/or download pump application programs and/or pump operational parameters. An example of lock levels are disclosed in United States Patent
20 6,475,180, issued on November 5, 2002 and entitled Drug Pump Systems and Methods, the entire disclosure of which is hereby incorporated by reference.

Once the user is logged in and the library is optionally loaded, the user selectively executes three different modules, a library module 708, a therapy module 710, and a protocol module 712.

25 The library module 708 assigns a label that identifies an entity and user attributes for the selected entity. Entity attributes can be properties specific to the library, such as a name of a doctor, a name of a healthcare regimen, or a location of the medical infusion pump or pump network, for example the hospital or department at which the pump is located. User attributes define the users allowed to
30 access and modify pump parameters for protocols associated with a particular library by using the medical infusion pump, and can also define users allowed to modify pump protocols using user software, as described below. The library module 708

contains a library definition module 714 and a user rights module 716, which are configured to perform these tasks, respectively.

The therapy module 710 adds and modifies therapies, qualifiers associated with the therapies, and drugs. The therapy module includes a therapy definition module 718, a qualifier definition module 720, and a drug definition module 722. The therapy definition module 718 controls addition and editing of therapies, which are the methods of patient treatment for diseases or generalized rehabilitation as previously described. The qualifier module 720 defines qualifiers and associates the qualifiers with one or more therapies. The drug definition module 722 defines one or more drugs that can be used in the medical infusion pump.

The protocol module 712 adds, edits, and defines protocols by associating therapies, qualifiers, and drugs with pump parameters to form libraries of pump protocols for a medical infusion pump. The protocol module 712 allows a user to select a therapy defined in the therapy definition module 718. The protocol module further allows the user to associate a qualifier defined in the qualifier definition module 720 with the selected therapy. For example, one or both of "adults" and "children 5-10 years" qualifiers can be associated with an epidural therapy. The protocol module 712 also associates one or more drugs with each therapy and qualifier combination, indicating that use of the drug is appropriate for that therapy and qualifier. The protocol module 712 guides the user in defining a protocol by assigning default pump parameters to be associated with the selected therapy, qualifier, and drug.

The protocol module 712 allows a user to associate more than one qualifier to each therapy, and also allows a user to associate more than one drug to each therapy and qualifier combination. For example, a protocol used in an epidural therapy for an adult can include a higher basal delivery rate parameter than a protocol used with a child for the same therapy. Likewise, usage of one drug can require a higher or lower dosage than another drug for the same therapy and qualifier because of concentration, reaction, or other factors.

Operational flow proceeds from the modules 708, 710, 712 to an optional export library module 724. The export library module 724 saves the defined or edited pump application programs and parameters in a file or other data structure that can be loaded by the administrative software 700 at another time or

location, or can be loaded by user software such as described below in FIGS. 26-41. If the export library module 724 does not execute, the library remains within the normally connected database 504 of FIG. 5, but is not extracted into a file for portability to a pump not connected to the network or to another medical infusion pump network. One or more libraries can be exported into a single library file.

Operation of the software 700 terminates at an end module 726. The end module 726 corresponds to termination of the administrative software 700 by clicking on a close window button on the computer or by some other mechanism for terminating execution of software.

In one embodiment of the administrative software 700, a user of the software 700 defines each protocol included in a library. In defining each protocol, the user assigns the index to the protocol, such as the therapy, qualifier, and drug defined in the modules 708, 710, 712 above. In a second possible embodiment, the administrative software includes a number of default settings or pump parameter modifications used when specific therapies, qualifiers, or drugs are selected. The user selects a therapy, qualifier, and drug to associate with a pump protocol. The administrative software 700 can include instructions dictating that selection of one or more of the therapies, qualifiers, and drugs sets or modifies one or more of the patient specific pump parameters or non-patient specific pump parameters. In one example of this second embodiment, a user setting a drug having a maximum safe consumption rate will trigger the administrative software 700 to preset an acceptable range of programmable delivery rates and a default delivery rate in the protocol, as well as alarms or other non-patient specific pump parameters. In another possible example of this second embodiment, a user setting a qualifier indicating a low age, such as "Children 5-10 years old", will set or adjust the protocol to result in a low delivery rate and demand dose being incorporated into the protocol, and will set one or more parameters related to alarms for use in a medical infusion pump.

Referring now to FIG. 8, a library import screen 800 is shown. The library import screen corresponds to the load library module 704, and is optionally used to load a library file containing one or more libraries of pump programs in the administrative software 700 of FIG. 7, above. The library import screen 800 includes a file selection field 802 and a password field 804. The term field as used

herein can include the window, or screen, generally, and can also include menus, selectable lists, or buttons within the window.

The file selection field 802 displays one or more library files that are available to be selected. The file selection field 802 allows a user to select one or
5 more of the protocol library files, in conjunction with selection control buttons 806a, 806b. The password field 804 controls access to the selected library file by requiring a user to input a correct password associated with the selected file, or library within the file.

A location field 808 presents the directory path of the selected library
10 file and a browse button 810 provides browsing capabilities to allow a user to find a library file other than those displayed in the library selection field 802.

In use, the library import screen 800 initially presents a listing of library files in the file selection field 802 available to the administrative software. The user selects one of the displayed library files, or presses the browse button 810
15 to search for additional library files. The user selects a library file by clicking on the displayed library file, or by other selection method. The directory path of the selected library file appears in the location field 808, and the user then enters a password in the password field 804 corresponding to the selected library file. The user confirms the choice using the selection control buttons 806a, 806b. The
20 administrative software confirms that the user-entered password is correct, and accesses the library file.

Referring now to FIG. 9, a user interface 900 for the administrative software 700 of FIG. 7 is shown. The user interface 900 includes features associated with the modules incorporated in the administrative software 700, and corresponds
25 to the screen displayed following the start operation 702 of FIG. 7. The user interface 900 includes a number of features providing global control and status information to a user.

Global control features included in the user interface 900 relate to library and pump access settings. A library field 902 provides a listing of libraries
30 currently loaded by the software. The library field 902 contains the libraries which have been loaded using the load library module 704 of FIG. 7. A login button 904 generates a login screen that checks whether a user has the right to modify pump parameters and protocols.

The user interface 900 can display information related to the status of the network of medical infusion pumps. A database identification field 906 identifies the database 504 currently connected to the administrator computer 502, as shown in FIG. 5. In additional exemplary embodiments, additional information can be displayed, such as the location of the server and/or medical infusion pumps connected to the infusion pump network or a serial number or other identification number of the pumps or libraries.

The user interface 900 also includes a location tab 920, a therapy tab 940, and a protocol tab 960. Referring back to FIG. 5, the location tab 920 corresponds to the location module 708, the therapy tab 940 corresponds to the therapy module 710, and the protocol tab 960 corresponds to the protocol module 712. By clicking on or otherwise selecting a tab, a user transfers operation to the module corresponding to that tab.

Referring now to FIGS. 10-12, the user interface 900 is shown, and location tab 920, therapy tab 940, and protocol tab 960 are described in greater detail. FIG. 10 shows the user interface 900 after selection of the location tab 920, or after initial execution of the login module 706 and optional load library module 704.

A region of the location tab 920 supersedes the therapy tab 940 and protocol tab 960. Within the region exposed by the location tab 920, a library selected in the library field 902 populates a library description field 1002 and a user accounts field 1004. The library description field 1002 describes the library that is currently loaded, and includes attributes of the location in which the library is used, or other information about the library. The location attributes can include the name of the hospital or the department in the hospital associated with the library. The additional information can include information related to the users of the library, the contents of the library, or other information. The library description field 1002 corresponds to the location attributes module 714 of FIG. 7. The user accounts field 1004 displays a listing of user accounts associated with the library. Each of the user accounts generally correspond to a user, but can also correspond to a location of a computing system 104 or medical infusion pump 102. The user accounts define the persons allowed to access and modify pump protocols, and/or parameters associated with the library at the point of care, such as a medical infusion pump 102 or

computing system 104 executing user software. The user accounts field 1004 enables a user to add, edit, or delete users from the listing using the buttons 1006a-1006c provided. The user accounts field corresponds to the user rights module 716 of FIG. 7.

5 In a possible embodiment, a log report field 1008 can optionally direct the server 206 to generate a record for various events occurring in the database 504. For example, a log report can be created each time a protocol is sent to a medical infusion pump 102. Alternately, a log report is created each time a library is sent to a medical infusion pump 102. In further embodiments, a log report
10 is created each time a library is edited or accessed. In still further embodiments, globally unique identifiers are entered into the log report related to instances where computing systems 104 and/or infusion pumps 102 access a library in the database 504. In the embodiment shown, the log report field 1008 is a selectable check box, but can be implemented as any other type of selectable field.

15 A password field 1010 sets a password for the currently selected library file. The password protects access to the library from the perspective of user software residing on either a computing system 104 or an infusion pump 102 such that only users with knowledge of the password associated with the library file can load the library using the load library module 704 of FIG. 7 and user interface 800.

20 Referring now to FIG. 11, the user interface 900 is shown with the therapy tab 940 selected. A region of the therapy tab 940 supersedes the location tab 920 and the protocol tab 960. The therapy tab 940 corresponds to the therapies module 710, which defines therapies, qualifiers, and drugs used by medical infusion pumps associated with the administrative software 700 of FIG. 7. The therapy tab
25 940 includes a therapy definition field 1102, a qualifier definition field 1104, and a drug definition field 1106.

 The therapy definition field 1102 corresponds to the therapy definition module 718 of FIG. 7, and includes a therapy listing 1108, a therapy notes field 1110, and control buttons 1112a-1112c. The therapy listing 1108 displays the
30 therapies currently defined in the library displayed in the library listing 802. Two example therapies, "Patient Controlled Analgesia" and "Epidural", are shown in the therapy listing 1108, and are associated with a library named "NeoNatal Intensive Care Unit". The therapy notes field 1110 contains administrative user-defined notes

related to the therapy selected in the therapy listing 1108. The notes relate to a therapy, and include information related to administration of the therapy, such as messages related to dosage, bolus amount, or administration. The therapy notes field 1110 presents administrative user-created notes to convey therapy-specific information to caregivers using or programming the pump. Control buttons 1112a-1112c allow a user to add, edit, and delete therapies from the therapy listing 1108 and therapy notes field 1110.

The qualifier definition field 1104 corresponds to the qualifier definition module 720 of FIG. 7, and includes a qualifier listing 1114, a qualifier notes field 1116, and control buttons 1118a-1118c. The qualifier definition field 1104 lists the qualifiers associated with the therapy selected from the list shown in the therapy listing 1108. For example, an administrative user might add a "Children 5-10 yrs" entry to the qualifier listing 1114. Each therapy relates to one or more qualifiers, such as a general qualifier, a weight based qualifier, or an age based qualifier. The qualifier notes field 1116 contains notes describing the difference in application of the therapy based on the qualifier selected. In the case of the "Children 5-10 years" qualifier, the qualifier notes field 1116 can indicate a lower than normal dosage to be administered. Notes associated with the qualifier display in the qualifier notes field 1116 only when the qualifier is selected.

The drug definition field 1106 corresponds to the drug definition module 722 of FIG. 7, and includes a drug listing 1120 and control buttons 1122a-1122d. The drug listing 1120 includes information related to the therapeutic fluid used in the medical infusion pump. This information can include the drug name, identification code, units, concentration, and usage. The control buttons 1122a-1122c allow a user to add, edit, or delete drugs in the drug listing 1120. A bar code control button 1122d generates a bar code screen including information related to a drug.

Referring now to FIG. 12, the user interface 900 is shown with protocol tab 960 selected. The protocol tab 960 supersedes the location tab 920 and the therapy tab 940. The protocol tab 960 corresponds to the protocol module 712 of FIG. 7, and associates therapies, qualifiers, and drugs to allow a user to define protocols by setting pump parameters for inclusion within pump programs. The protocol tab 960 includes a protocol field 1202 and control buttons 1204a-1204e.

The protocol field 1202 lists the combinations of therapies, qualifiers, and drugs for which a protocol is defined. The protocol field 1202 also lists a pump type 1203 for which each protocol is defined. By specifying a pump type for which each protocol is defined, it is possible to enable pump-specific protocol programming, while still
5 using the administrative software 700 and user software, described below, for all pump types configurable using the protocol-based programming scheme described herein.

The control buttons 1204a-1204e operate to add, edit, view, and/or delete the protocol for the combinations of therapies, qualifiers, and drugs. The
10 control buttons 1204a-1204d allow a user to set pump parameters so as to define the protocol. Control button 1204e generates a prescription form screen representing the protocol for the selected therapy, qualifier, and drug.

FIG. 13 shows a bar code screen 1300 displaying a bar code for a drug defined in the administrative software 700 of FIG. 7. A drug identification
15 code is associated with each drug in the drug listing 1120 of FIG. 11. A drug is selected from among those in the drug listing 1120 of FIG. 11, and the bar code screen 1300 displays upon clicking on the control button 1122d.

The bar code screen 1300 includes a print preview field 1302, a printer drop down menu 1304, a label information menu 1306, a copies drop down
20 menu 1308, and control buttons 1310a-1310c. The print preview field 1302 displays a bar code associated with the selected drug. The printer drop down menu 1304 lists available printers configured to print the bar code. The label information drop down menu 1306 defines the label configuration to which the bar code is directed. The label configuration includes, for example, the size and layout of the label paper. The
25 copies drop down menu 1308 dictates the number of copies of the bar code that are printed. Control buttons 1310a-1310c provide printing, cancellation, and help procedures to a user.

In use, the bar code corresponds to a drug identification code associated with the drug. In a possible embodiment, a pharmacist can print the
30 barcode associated with the drug using the administrative software 700 via the bar code screen 1300 and affix the printed bar code label to the drug container. The labels indicate the drug and dosage being delivered by the pump. This provides an

easily accessible and visually prominent indication of the drug being delivered by a medical infusion pump.

A caregiver connecting a drug to a medical infusion pump 102 will scan the bar code on the drug container, which will correspond to the drug
5 identification code associated with the drug in the administrative software 700. This ensures that the caregiver affixes the drug to the pump which corresponds to the protocol selected using administrative software 700.

FIG. 14 shows a prescription form screen 1400 displaying prescription information related to a specific therapy, qualifier, and drug. The
10 prescription form screen displays upon selection of control button 1204e on protocol tab 960. The prescription form screen 1400 corresponds to the currently selected protocol, and therefore incorporates information specific to the selected protocol. Prescription information includes, for example, directions for application of the drug according to the defined protocol, and parameters such as drug delivery rates or
15 bolus amounts. The prescription form also includes usage notes describing operation or application of the selected therapy, qualifier, and drug.

A doctor using the prescription form dictates the protocol used in the pump by using a prescription form analogous to the prescription form screen. The prescription form 1400 generated by the administrative software will correspond to
20 the prescription form completed by the doctor. In the embodiment shown, the prescription form screen 1400 that is generated includes information specific to the drug and therapy administered, which may be notes related to administration of the drug and therapy as dictated by the doctor. For example, the prescription form will include drug information, such as the name, type, concentration, and notes regarding
25 the drug, and will also include information related to patient specific pump parameters associated with a selected therapy, qualifier, and drug selected from a library.

Referring now to FIG. 15-17, a protocol definition user interface 1500 defines the relationships between the therapies, qualifiers, and drugs added to
30 the library by the therapy module 710 and using the therapy tab 960. A user defines a protocol by setting an index by sequentially selecting a therapy, a qualifier, and a drug, and then by associating pump parameters with that index. In setting the index using the user interface 1500, a therapy drop down menu 1502 connects to a therapy

notes field 1504, a qualifier drop down menu 1506 connects to a qualifier notes field 1508, and a drug drop down menu 1510 connects to a drug notes field 1512.

FIG. 15 shows the initial state of the protocol definition user interface 1500. A list of therapies appears in the therapy drop down menu 1502, while the
5 remaining menus 1506, 1510 and fields 1504, 1508, 1512 remain blank.

FIG. 16 shows the state of the protocol definition user interface 1500 after a therapy is selected in the therapy drop down menu 1502. Notes related to the selected therapy appear in the therapy notes field 1504, and a list of qualifiers associated with the selected therapy appears in the qualifier drop down menu 1506.
10 The notes correspond to the therapy notes entered in the therapy notes field 1110. The list of qualifiers corresponds to the qualifiers associated to the therapy by listing in the qualifier listing 1114. For example, two qualifiers, "Adults & Children over 5 yrs" and "Children 3-5 yrs" can be associated with the "Patient Controlled Analgesia" therapy. When that therapy is selected in the therapy drop down menu
15 1502, the two associated qualifiers populate the qualifier drop down menu 1506.

FIG. 17 shows the state of the protocol definition interface 1500 after a qualifier is selected in the qualifier drop down menu 1506. Notes related to the qualifier appear in the qualifier notes field 1508, and a list of drugs available within the library or database appears in the drug drop down menu 1510. The notes
20 correspond to changes in the therapy due to the qualifier selected. Continuing the example from FIG. 14, the "Adults & Children over 5 yrs" qualifier is selected. Notes related to customization of the "Patient Controlled Analgesia" therapy based on the selected qualifier are displayed in the qualifier notes field 1508. Furthermore, a list of four drugs/drug concentrations appear in the drug drop down menu 1510,
25 which are the drugs defined and available within the library or database.

Analogously, upon selection of a drug from the drug drop down menu 1510, drug notes (not shown) appear in the drug notes field 1512.

Referring now to FIGS. 18-24, a parameter user interface 1800 provides additional tabbed screens allowing a user to define pump parameters. The
30 pump parameters complete the protocol definition associated with the therapy, qualifier, and drug combination selected in FIGS. 15-17. As described in conjunction with FIG. 7, selection of the index defined by a therapy, qualifier, and drug using the protocol definition interface 1500 may trigger the administrative

software 700 to set or modify one or more of the pump parameters. Alternately, none of the pump parameters are set during the definition and selection of the index formed by the therapy, qualifier, and drug using the protocol definition interface 1500. In such an embodiment, the parameter user interface 1800 is used to define the pump parameters that are programmable into a medical infusion pump.

5 The parameter user interface 1800 includes a status region 1802, a protocol activation field 1804, and control buttons 1806a-1806c. The status region 1802 displays the therapy, qualifier, and drug associated with the assignable pump parameters. The protocol activation field 1804 publishes the protocol within the library such that the protocol is visible to user software accessing the library when it resides within the database 504 of FIG. 5. The protocol activation field 1804 allows a user of the administrative software 700 to control when protocols become available for use by user software. In some circumstances, it can be advantageous to prevent user software from accessing data, particularly while that data is being edited in the administrative software 700. One example of user software used to access pump protocols is illustrated below in conjunction with FIGS. 26-41. The control buttons 1806a-1806c provide save, cancel, and help options to a user of the administrative software.

20 The parameter user interface 1800 further includes a number of tabs, including a drug delivery tab 1810, a secondary drug delivery tab 1820, an alarm tab 1830, a security tab 1840, a display/sound tab 1850, and a report tab 1860. Parameters set within each of the tabs are discussed in FIGS. 18-24.

FIG. 18 shows the parameter user interface 1800 with the drug delivery tab 1810 selected. The drug delivery tab 1810 allows a user to set drug delivery rate parameters, and includes a general settings region 1812 and a patient specific parameters region 1814.

25 The general settings region 1812 includes verification settings 1816 and weight based settings 1818. The verification settings 1816 includes drug verification and caregiver verification settings. Specifically, the verification settings require that a caregiver verifies that the correct drug is provided to the medical infusion pump. The verification settings also require a second caregiver to verify the settings of the medical infusion pump. The weight based settings 1818 set a weight based protocol at a programmable, variable weight limit. By weight based

protocol, it is intended that dosage delivery rates, boluses, thresholds, and other delivery parameters change from a "dosage per hour" basis to a "dosage per weight factor" rate, where the weight factor can be on a per unit measure weight basis for the user of the medical infusion pump 102 (i.e. "per kilogram" or other), or based on
5 the user's body surface area, a weight based therapy, or other options.

The patient specific parameters region 1814 includes a continuous rate region 1822, a demand dose region 1824, and a demand dose lockout region 1826. Continuous rate refers to the constant drug delivery rate of the medical infusion pump, also referred to as the basal rate. Demand dose refers to an added
10 drug delivery bolus amount delivered by the pump upon a demand by a patient. Demand dose lockout refers to the time interval after a demand dose is delivered, during which another demand dose will not be delivered by the pump.

The continuous rate region 1822 includes a meter, shown as a slider bar 1828 and an indicator 1829. The meter generally has two or more locations,
15 each corresponding to a parameter value that can be programmed in the medical infusion pump. Generally, the positional relationship of the meter indicates the setting of the meter. In a possible embodiment of the slider bar 1828 shown, the indicator 1829 is movable relative to the slider bar 1828 to set a default value, or "initial value" continuous drug delivery rate parameter. In a second possible
20 embodiment, the default value is set using an initial value gauge 1832.

The continuous rate region also includes hard limit gauges 1834, soft limit gauges 1836, and user interface options 1838a-1838c. The initial value gauge 1832, hard limit gauges 1834, and soft limit gauges 1836 include values, which may include numerical ranges. The hard limit gauges 1834 set a hard maximum and hard
25 minimum which form an acceptable pump programming range. The range of acceptable pump activity represents the absolute maximum and minimum values programmable into the pump by user software, as described below. This configuration allows for control of the range of values visible to a user of the medical infusion pump or associated computing system.

30 The limits set by the soft limit gauges 1836 represent a manually exceedable threshold value. The soft limit can be overridden by a user of a medical infusion pump on a pump-by-pump basis. Pump activity outside the range defined by soft limits can trigger an alarm or otherwise alert a caregiver that a pump is

functioning outside of the usual operational range of the pump. A variety of alarm levels or alerts can be set by the soft limit gauges 1836. For example, the alert can be a flag set in the software. The alert could additionally be an audible alarm, or a visual indicator displayed on at least a portion of the medical infusion pump. The
5 visual indicator could be a flashing indicator or changed/changing color on the display of the medical infusion pump.

In a second possible embodiment, the hard limit gauges set a non-limiting range, and the user software described below can be programmed within its full operational range. In such an embodiment, pump activity outside the range set
10 by the hard limit gauges 1834 can trigger an alarm or otherwise alert a caregiver that a pump is functioning outside of the usual operational range of the pump. In this embodiment, the soft limit gauges 1836 set a narrower range, operation outside of which can trigger a warning or second alarm indicating pump activity outside of an expected range of pump operation. This warning or second alarm indicates a pump
15 condition less serious than the alarm triggered by the hard limits.

The user interface options 1838a-1838c enable the user software to display and edit the delivery rate, editing of the delivery rate, and requiring comments by users of a pump who wish to exceed the soft limits when setting the drug delivery rate. Selection of the display option 1838a publishes the pump
20 parameter so that the value is visible to a user of the pump or computing system.

The demand dose region 1824 includes a slider bar 1842 and an indicator 1843. The slider bar and indicator operate in a similar manner to the slider bar 1828 and indicator 1829 in the continuous rate region 1822, but control demand dose settings. Likewise, the demand dose region 1824 includes an initial value
25 gauge 1844, as well as hard limit gauges 1846 and soft limit gauges 1848 setting visible thresholds and triggering alarms as in the continuous rate region 1822. Demand dose options 1852a-1852c provide analogous display, editing, and comment options to the user interface options 1838a-1838c.

The demand dose lockout region 1826 includes a slider bar 1854 and
30 an indicator 1855, and also includes an initial value gauge 1856, hard limit gauges 1858, and soft limit gauges 1862. Each of these features functions analogously to those discussed above in conjunction with the continuous rate region 1822. The

demand dose lockout region also includes lockout options 1864a-1864c analogous to the user interface options 1838a-1838c.

FIG. 19 shows the parameter user interface 1800 with the drug delivery tab 1810 modified to provide a weight based drug delivery protocol. The modification of the user interface 1800 occurs in the drug delivery tab 1810 upon user selection of the weight based settings 1818 discussed in FIG. 18. The continuous rate region 1822 and demand dose region 1824 are modified to reflect dosage rates on a "per kilogram" or other weight measure basis.

FIG. 20 shows the parameter user interface 1800 with the secondary drug delivery tab 1820 selected. The secondary drug delivery tab 1820 provides additional medical infusion pump programming options for assigning pump parameters in a specific protocol. The secondary drug delivery tab 1820 includes a dosing limit region 2002, a patient specific parameter region 2004, a titration region 2006, a maximum delivery rate gauge 2008, and a maximum clinician bolus gauge 2010.

The dosing limit region 2002 displays limits for total drug delivery within a specified amount of time. The dosing limit region 2002 includes options for setting a limit on doses per hour, a timed medication delivery limit, or other limits. A user selects one of the options for setting the dosage limit.

The patient specific parameter region 2004 sets parameters related to the dosage limits coordinated to the options in the dosing limit region 2002. The patient specific parameter region 2004 includes a timed delivery limit region 2012, an hourly demand doses region 2014, and a reservoir region 2016.

The timed delivery limit region 2012 sets the delivery limit on a per assigned time period when the option for setting a limit on doses per hour is selected in the dosing limit region 2002. The timed delivery limit region 2012 includes a meter, shown as a slider bar 2018 and indicator 2019. The timed delivery limit region 2012 also includes an initial value gauge 2020, hard limit gauges 2022, soft limit gauges 2024, and control options 2026a-2026c. Operation of the slider bar 2018, indicator 2019, gauges 2020-2024, and control options 2026a-2026c is analogous to operation of the slider bar features as discussed in conjunction with FIG. 18, above.

The hourly demand doses region 2014 sets the delivery limit on a per hour basis when the doses per hour limit is selected in the dosing limit region 2002. The maximum demand doses region 2014 includes a slider bar 2028, indicator 2029, gauges 2030-2034, and control options 2036a-2036c, operation of which is likewise analogous to operation of the slider bar features discussed in FIG. 18.

The timed delivery limit region 2012 and hourly demand doses region 2014 are operated in the alternative, in that only one of the two regions is active at one time. The region that is active depends upon the option selected in the dosing limit region 2002. Selection of a timed delivery limit in the dosing limit region 2002 activates the timed delivery limit region 2012 and deactivates the hourly demand doses region 2014. Selection of an hourly delivery limit activates the hourly demand doses region 2014 and deactivates the timed delivery limit region 2012.

The reservoir region 2016 sets the initial volume settings and display settings for tracking the volume of fluid in the reservoir attached to the medical infusion pump. The reservoir region 2016 includes a meter, shown as a slider bar 2040 and indicator 2041, operation of which is analogous to the slider bar and indicator discussed in conjunction with FIG. 18. The reservoir region further includes an initial value gauge 2042, a disable option 2044, and control options 2046a-2046b. The initial value gauge 2042 provides an alternate method for setting the initial value of the reservoir volume to the slider bar 2040 and indicator 2041. The disable option 2044 disables the reservoir volume monitor. The control options 2046a-2046b provide display and edit capabilities to a user of the pump.

The remaining regions, i.e. the titration region 2006, the maximum deliver rate gauge 2008, and the maximum clinician bolus gauge 2010, set pump specific settings related to drug delivery limits. The titration region 2006 enables or disables titration in the medical infusion pump, and sets an optional titration limit in the pump. The maximum delivery rate gauge 2008 sets a maximum delivery rate for the infusion pump. The maximum delivery rate is measured in milliliters per hour, and includes both the basal delivery rate and bolus delivery. The maximum clinician bolus gauge 2010 sets the maximum bolus which can be delivered by a caregiver. The maximum clinician bolus may be a larger bolus than the standard

patient-controlled bolus, but must be administered under the supervision of a caregiver. Other regions can be included in the tab 1820 as well.

FIG. 21 shows the parameter user interface 1800 with the alarm tab 1830 selected. The alarm tab includes a number of regions used for maintenance and hardware-related alarms, as opposed to the drug delivery threshold alarms discussed above in conjunction with the slider bars. The alarm tab 1830 allows the user to enable and disable alarms in the medical infusion pump. The alarm tab 1830 includes a reservoir low alarm region 2102, a reservoir empty alarm region 2104, a pump stop alarm 2106, a maintenance alarm 2108, and a hardware alarm 2110. The reservoir low alarm region 2102 provides an alarm indicator when a drug supply volume falls below a threshold level. The threshold level can be a standard level or an administrative user-set volume level. The reservoir empty alarm region 2104 sets a single occurrence or repeating alarm when a drug supply volume falls below a threshold level. The pump stop alarm 2106 sets an alarm which occurs when the medical infusion pump stops operating. The maintenance alarm 2108 enables a maintenance alarm, which alerts a user when maintenance is needed. The hardware alarm 2110 provides options for detection of optional components used with the infusion pump. For example, the hardware alarm 2110 can trigger upon detection of an air detector or other component. The pump hardware detector region 2110 can provide the option of enabling an alarm sent to a caregiver.

FIG. 22 shows the parameter user interface 1800 with the security tab 1840 selected. The security tab 1840 includes options related to security of the medical infusion pumps. The security tab 1840 includes a security region 2202 and an epidural region 2204. The security region 2202 presents a number of security options for controlling access to each medical infusion pump. Security options can include an automatic lock, a lock code, clinician code, and an initial lock level. The epidural region 2204 selectably places the pump into a mode configured for epidural therapy.

FIG. 23 shows the parameter user interface 1800 with the display/sound tab 1850 selected. The display/sound tab 1850 includes options related to the display and sound settings for each medical infusion pump. The display/sound tab 1850 includes an auto-review option 2302, a main display region 2304, a units option 2306, a date format option 2308, a send protocol region 2310,

and a sound region 2312. The auto-review option 2302 enables a power-up display of options in the medical infusion pump. The main display region 2304 presents a number of programmable fields that will, by default, be displayed on the medical infusion pump. For example, the main display region 2304 includes programmable text display entry, power source information, and drug delivery rate. Other display options related to properties of the medical infusion pump are available as well. The units option 2306 selectably assigns a units location for display on the medical infusion pump. The date format option 2308 assigns date formats for displaying on the medical infusion pump. The send protocol region 2310 sets one or more patient markers or date/time stamps in the infusion pump network upon distribution of programming instructions to the medical infusion pump. The sound region 2312 enables sound in the medical infusion pump, such as beeping sounds when one or more keys/key sequences are depressed.

FIG. 24 shows the parameter user interface 1800 with the report tab 1860 selected. The report tab 1860 presents options for displaying reports by a medical infusion pump related to events occurring in the pump. The report tab 1860 includes a new patient marker region 2402 and a custom report region 2404. The new patient marker region 2402 includes a list of options to perform related to report generation upon association of a medical infusion pump with a new patient. For example, the new patient marker region 2402 includes a number of report-clearing options and power-on options to be performed when a medical infusion pump is assigned to a new patient. The custom report region 2404 provides a number of selectable options for display in a custom report regarding events tracked in the medical infusion pump network. The custom report region 2404 provides all of the options tracked in the infusion pump network 500, which can be stored in the log files 516. A user selects one or more of the options to generate a report. For example, the custom report region 2404 can include dose counters, doses per hour, pain scale information, drug delivery information, an event log, and a patient marker. The report generated by the custom report region 2404 options displays on the screen of the medical infusion pump 102 or computing system, 104. Both the new patient marker region 2402 and the custom report region 2404 can include additional options as well.

Referring now to FIG. 25, a library export screen 2500 is optionally used to save the pump parameters and protocols for exporting from the database 504 of FIG. 5. One or more libraries are exported using the library export screen 2500, including all of the pump protocols, and containing all of the defined therapies, 5 qualifiers, drugs, and combinations thereof. The library export screen 2500 corresponds to the export library module 724 of FIG. 7, and provides file extraction, in that one or more protocol libraries built using the administrative software 700 are exported to a data file or database. The export library screen 2500 includes a library field 2502 and an export option region 2504. The library field 2502 displays a 10 number of libraries available to the system that can be exported to a portable file. The export option region 2504 provides export options to the administrative user seeking to export the library data to a file. For example, export options can include creating a read-only file, or a read-write template for use in another instantiation of the software system disclosed herein. The read-only file might be selected in the 15 case that the library is to be loaded onto a medical infusion pump that is disconnected from the database 504 of FIG. 5. The read-write file might be selected if the library is to be transferred to a separate infusion pump network 500 altogether, in which administrative software on that subsequent infusion pump network may be used to subsequently edit that library. Additionally, the export option region 2504 20 includes an assignable password to add security to the library file exported such that a user attempting to access the protocols contained in the library file must know and enter the correct password.

The above description and figures corresponding to the administrative software 700 provides a therapy-centric programming schema for a 25 medical infusion pump. For example, a certain drug used in conjunction with a medical infusion pump may be appropriate for use with a specific therapy for an adult, but may not be appropriate for the same therapy for a child. Certain drugs may only be appropriate in certain therapies, and under certain qualifying conditions. Pump parameters are initially set according to the protocols defined in 30 the administrative software 700, but are customizable on a pump-by-pump basis using user software associated with a specific pump and/or patient.

FIG. 26 illustrates exemplary architecture of user software 2600 for accessing a pump application program and programming a medical infusion pump.

The software 2600 can operate within the pump 102, computing system 104, or a combination thereof. The user software 2600 allows a user, for example a doctor, nurse, pharmacist, or other caregiver, to select and customize pump application protocols, and parameters for execution in and control of a medical infusion pump
5 102. Depending upon the pump configuration, the user may select a protocol or a library for loading into a medical infusion pump. Additional data structures could be loaded into the medical infusion pump as well. Although the user software 2600 is discussed in conjunction with the administrative software 700 previously described in FIGS. 7-25, it is understood that the user software systems described
10 herein are operable in conjunction with additional hardware/software embodiments.

The medical infusion pumps as described store pump data in memory, such as the memory shown above in FIG. 3. The pump data can include pump parameters, parameter values, programs, and other functional and data systems configured to operate the medical infusion pump. As referred to herein, a
15 set of pump parameters can include the entire memory contents of the pump, or can include a subset of the memory contents, such as selected data values, that can be altered to change operation of the pump.

The user software 2600 is instantiated by a start module 2602. The start module 2602 corresponds to initial execution of the user software 2600 by
20 clicking on an icon on the computer or by some other mechanism for executing software. Upon execution of the start module 2602, the user software 2600 connects to a library on a server 206 containing one or more pump protocols.

Following the start module 2602, operational flow optionally proceeds to a library import module 2604. The library import module 2606 provides
25 the ability to import one or more libraries into the software 2600. This feature can be used by a computing system 104 or medical infusion pump 102 that is not connected to the same medical infusion pump network 500 as the database storing the library 504. If the computing system 104 or medical infusion pump 102 is connected to the medical infusion pump network 500, each component is by default
30 connected to the server 206 and database 504.

The libraries available to be imported include pump protocols and parameters, and may have been created using the administrative software of FIGS. 7-25. The collection of pump protocols can be accessed from the server 206 or in

one or more individual computing systems 102. The library import module 2604 allows the user to select one or more pump application programs for downloading to a medical infusion pump 102.

Once connected to the desired library either by default or via the library import module 2604, operational flow proceeds to a login module 2606. The login module 2606 regulates user rights in the software 2600 by controlling access to the libraries 508 in the database 504 of FIG. 5. User rights define access levels in the software for users such as doctors, nurses, other caregivers, or patients. A user will have a set access level allowing the user to view or edit pump application protocols and parameters within the user software 2600. Access levels are set using the user rights module 716 of the administrative software 700, described above in conjunction with FIG. 7. Access levels can be set by a user of the administrative software 700 according to a variety of criteria, such as the type of caregiver (e.g., physician, nurse, or pharmacist).

Different access levels also can provide different rights with respect to pump operational parameters. For example, one access level might give a user a right to edit patient specific pump parameters. One access level might permit a user the right to only view and download the patient specific pump operational parameters. Different embodiments can include the ability to provide an access level for a user any combination of rights to edit, view, and/or download pump operational parameters, protocols, or libraries.

Once the user is logged in, the user selectively executes three different modules, a protocol selection module 2608, a task module 2610, and a report module 2612.

The protocol selection module 2608 selects a protocol for use with a medical infusion pump from the protocols loaded in the user software 2600. The protocol selection module 2608 guides the user through selection of a therapy, qualifier, and drug combination defined to be a protocol by the administrative software. The protocol selection module 2608 includes a therapy selection module 2614, a qualifier selection module 2616, and a drug selection module 2618 for this purpose. The therapy selection module 2616 selects a therapy to be administered by the drug infusion pump 102. The therapy is one of the therapies included in the library selected in the library import module 2606. The qualifier selection module

2616 selects a qualifier from those associated with the therapy in the library. The drug selection module 2618 selects a drug associated with the therapy and drug. The protocol selection module 2608 further allows customization of the protocol by allowing a user to modify pump parameters, such as the drug delivery rate, the demand dose, the demand dose lockout, drug delivery limits, and reservoir volume.

The task module 2610 guides a user through maintenance and monitoring tasks that are required for each medical infusion pump 102. These maintenance and monitoring tasks can include pump settings comparison and testing, as well as changing the reservoir holding the drug delivered by the medical infusion pump 102. The task module 2610 includes a pump settings module 2620, a comparison module 2622, and a reservoir module 2624. The pump settings module 2620 compares the local pump settings to a standing order set for the pump, for example by a caregiver who programmed or customized the medical infusion pump. During operation of the pump settings module 2620, the user software 2600 receives a GUID specific to a protocol and generated by the server 206, and stores the GUID on the pump or computing system. The GUID generated by the server and stored on the pump or computing system is made available to the server when the pump or computing system accesses the library to look up and verify the future access of the correct protocol and/or library. This ensures that the pump settings are compared to the correct protocol stored on the server. The comparison module 2622 compares the local pump settings to a protocol, such as the protocols defined using the administrative software 700. The reservoir module 2624 determines if a drug reservoir is nearly empty and guides a patient or caregiver through the drug cartridge changing process occasionally required during use of a medical infusion pump.

The report module 2612 generates a report from preexisting logged information for a selected medical infusion pump 102. The report module includes a location module 2626 and a type module 2628. The location module 2626 requests the location of the pump from which a report is generated, such as a specific pump or a previously saved report stored on the pump or computing system. The type module 2628 presents a number of types of reports which can be generated from the logged information, such as a drug delivery report or event log, and can display the report responsive to a user request.

Operation of the software terminates at an end module 2630. The end module 2630 corresponds to termination of the administrative software 2600 by clicking on a close window button on the computer or by some other mechanism for terminating execution of software.

5 FIG. 27 shows a library import screen 2700 used to optionally load a library of pump protocols into the user software 2600. The library import screen 2700 corresponds to the login module 2604, and presents a location field 2702 and a library field 2704 used to browse to and select a library from a library database or file. An optional password field 2706 accepts user input to allow the user to enter a
10 password for accessing a password protected library file, such as one created using the export library screen shown in FIG. 25. Control buttons 2708a-2708b provide confirmation and cancellation options to the user, allowing the user to complete the library access process.

 FIGS. 28-32 illustrate an exemplary process and user interface
15 through which the protocol selection module 2608 leads a user to select a protocol for use with a medical infusion pump. The methods, systems and user interfaces described in conjunction with FIGS. 28-32 can be performed either on a computing system, such as a computing system 104 associated with a medical infusion pump 102, or on a medical fusion pump 102 configured to accept a library of pump
20 protocols directly. Medical infusion pumps 102 accepting a single pump protocol use a corresponding computing system 104.

 Referring to FIG. 28, operation of the protocol selection module 2608 is instantiated by a start module 2802. The start module 2802 corresponds to initial execution of the user software 2600, or initial selection of the protocol selection
25 module 2608 by clicking on an icon or tab on the computer or by some other mechanism for executing software.

 Following the start module 2802, operational flow proceeds to a load protocols module 2804. The load protocols module 2804 populates the user software 2600 with the protocols from the loaded library. For example, the load
30 protocols module 2804 can populate a listing of protocols for selection using the user software 2600.

 Following the load protocols module 2804, operational flow proceeds to a select protocol module 2806. The select protocol module 2806 selects a

protocol from among the protocols loaded into the user software 2600 by guiding a user through selection of a therapy, qualifier, and drug defining one of the protocols loaded in the software 2600. The select protocol module 2806 corresponds to the therapy selection module 2614, qualifier selection module 2616, and drug selection
5 module 2618 shown in FIG. 26.

Following the select protocol module 2806, operational flow proceeds to a settings module 2808. The settings module 2808 provides editing and customization of the pump parameters assigned to a medical infusion pump as dictated by the protocol selected by the user. The parameters include, for example,
10 the drug delivery rate, demand dose rate, or demand dose lockout.

Following the settings module 2808, operational flow proceeds to an optional pump programming module 2810. The pump programming module 2810 programs a medical infusion pump with the settings both as defined by the protocol and selected by the select protocol module 2806, and as customized by the settings
15 module 2808. The pump programming module 2810 executes if the user software 2600 resides on a computing system 104 connected to a medical infusion pump 102. The pump programming module may not execute if the software 2600 resides on the medical infusion pump 102 itself, because the protocols are already loaded into the pump alongside the library accessed by the user software 2600.

20 After the pump programming module 2810 completes, operation of the protocol selection module 2608 terminates at an end module 2812. The end module 2812 corresponds to successful programming of the medical infusion pump.

FIGS. 29-32 illustrate a user interface 2900 used to guide a caregiver through the pump programming process. The user interface 2900 operates on a
25 computing system associated to a medical infusion pump. The user interface 2900 includes a login button 2902, a connection status indicator 2904, a library indicator 2906, and a protocol tab 2920, a tasks tab 2940, and a reports tab 2960.

The login button 2902, when selected, generates a login screen that checks whether a user has the right to access pump programs, protocols, or
30 parameters. The connection status indicator 2904 displays the connection status of the user software. Connection status can include a connection to a medical infusion pump or connection to a networked server. The library indicator 2906 displays the current library loaded using the library import screen 2700 of FIG. 27.

Referring now to FIG. 29, the user interface 2900 is shown with the protocol tab 2920 selected. FIG. 29 corresponds to an initial state of the protocol selection module 2608 following the load protocols module 2804 of FIG. 28. The protocol tab 2920 guides a user through the process of selecting a protocol, customizing one or more parameters in the protocol, and programming a medical infusion pump with the customized pump program. The protocol tab includes a therapy selection field 2908, a therapy notes field 2910, a qualifier selection field 2912, a qualifier notes field 2914, a drug selection field, 2916, a drug notes field 2918, and a continue button 2922.

The therapy selection field 2908 lists the therapies included in the currently loaded library. For example, the two therapies shown are "Epidural" and "Patient Controlled Analgesia". The therapy notes field displays the notes associated with the selected therapy. In the initial state, the therapy selection field 2908 and therapy notes field 2910 are active, and the qualifier fields 2912, 2914, drug fields 2916, 2918, and continue button 2922 are inactive. No therapy is initially selected in the therapy listing field 2908, so the therapy notes field 2910 remains empty.

FIG. 30 shows the user interface 2900 with the protocol tab 2920 selected and a therapy selected from the therapy selection field 2908. Notes related to the selected therapy appear in the therapy notes field 2910, and the qualifier selection field 2912 and qualifier notes field 2914 activate. A listing of qualifiers associated with the selected therapy appears in the qualifier selection field 2912. The therapy notes shown recite "Epidural Patient Controlled Analgesia" corresponding to the selected therapy, but could contain particular information related to the therapy, such as warnings, descriptions, or other information about application of the therapy. The qualifiers, which appear once the therapy is selected, are shown to include "Adult and Child over 5" and "Child 5 years and under". No qualifier is initially selected, so the qualifier notes field 2914 remains empty. The drug selection field 2916, drug notes field 2918, and the continue button 2922 remain inactive.

FIG. 31 shows the user interface 2900 with the protocol tab 2920 selected and both a therapy selected from the therapy selection field 2908 and a qualifier selected in the qualifier selection field 2912. Notes related to the qualifier

appear in the qualifier notes field 2914, and the drug selection field 2916 and drug notes field 2918 are active. For example, "Adult and Child over 5" is shown to be, and the qualifier notes field 2914 displays specific notes applicable to those patients. A listing of drugs associated with the therapy and qualifier appears in the drug selection field 2916. Three exemplary drug menu listings including "Fentanyl 10 mcg/ml", "HYDRO Morphine 1 mg/ml" and "Morphine 1 mg/ml" are shown. No drug is initially selected, so the drug notes field 2918 remains empty. The continue button 2922 remains inactive.

FIG. 32 shows the user interface 2900 with the protocol tab 2920 selected and a therapy, qualifier, and drug selected in each of the respective selection fields 2908, 2912, and 2916. Once a therapy, qualifier, and drug are selected a specific protocol is designated from among the protocols defined in the administrative software 700. Each notes field 2910, 2914, and 2918 displays information related to the selected therapy, qualifier, and drug, respectively. The continue button 2922 activates, allowing the user to continue with customization of the selected protocol by changing one or more pump parameters related to the therapy, qualifier, and drug.

FIG. 33 illustrates an exemplary process by which the administrative software 2600 customizes one or more parameters of a selected pump protocol, and corresponds to the settings module 2608 of FIG. 26. The settings module 2608 allows customization of one or more pump parameters while monitoring whether the customized value is within an acceptably safe dosage or drug delivery range.

The settings module 2608 is instantiated by a start module 3302. The start module 3302 corresponds to selection of the confirmation button 2922 of FIGS. 29-32.

Following the start module 3302, operational flow proceeds to a display module 3304. The display module 3304 displays the default protocol settings of the protocol loaded onto a medical infusion pump screen or a computing system associated with the pump. The display module 3304 presents a number of meters to a user related to the drug delivery rates and other parameters controlled by the pump. The meters provide user controls for modifying one or more pump parameters.

In various embodiments, the display module 3304 can be configured to display a variety of coloring and image features. In one possible embodiment, the meters are slider bars that can include graphical thresholds set by the administrative software. A cautionary color change of the user interface (i.e. green or gray to
5 yellow or red) can represent a warning to the user that the current setting is outside of the administratively set thresholds.

In another possible embodiment, the overall background color of the user interface is color-coded to correspond to hospital coding procedures, and can represent one or more location-specific warning or status conditions. Additionally,
10 the color coding can be located behind an image displayed on the pump screen, and can be keyed to a location of the medical infusion pump, the drug administered, or a warning condition within the pump.

Following the display module 3304, operational flow proceeds to a custom settings module 3306. The custom settings module 3306 receives the
15 current customized pump settings based on user-customization of one or more pump parameters. The custom settings module 3306 can provide a user customization interface for setting pump parameters to values other than the initial or default values set in the administrative software 700.

Following the custom settings module 3306, operational flow
20 depends upon the implementation of the hard limits and soft limits in the administrative software 700. If the hard limit gauges in FIGS. 18-20 provide a warning but do not dictate an absolute maximum/minimum for the range of programmable values within the user software 2600, operational flow proceeds to a hard limit determination operation 3308. If the hard limit gauges dictate an absolute
25 maximum/minimum for the range of programmable values within the user software 2600, the hard limit will likely not be exceeded, so operational flow proceeds directly to the soft limit determination operation 3312.

The optional hard limit determination operation 3308 determines if the pump settings are outside the "hard limits" set in the administrative software
30 700. If the pump settings exceed the hard limit (i.e. above the maximum or below the minimum value), operational flow branches "yes" to a hard limit indicator module 3310. If the pump settings do not exceed the hard limit, operational flow branches "no" to a soft limit determination operation 3312.

The optional hard limit indicator module 3310 executes in conjunction with the hard limit determination operation 3308, and generates an indicator to a user that the hard limit set in the administrative software is exceeded by the current settings of the medical infusion pump. If the hard limit determination operation 3308 is bypassed or otherwise absent from the user software 2600, the hard limit indicator module 3310 can be absent/bypassed as well. The hard limit indicator module 3310 creates an alert indicator on the display of the pump or associated computing system, or sends an alert to the server or other computing system to alert a caregiver that an alert condition has been reached by the pump due to exceeding the hard limit. Operational flow proceeds to the display module 3704 to update the display and to allow additional user modification of the pump settings.

The soft limit determination operation 3312 determines if the pump settings are outside the "soft limits" set in the administrative software 700. If the pump settings exceed the soft limit, operational flow branches "yes" to a soft limit indicator module 3314. If the pump settings do not exceed the soft limit, operational flow branches "no" to return to the display module 3304.

The soft limit indicator module 3314 generates an indicator to a user that the soft limit set in the administrative software is exceeded by the current parameter settings. The soft limit indicator module 3314 creates an alert indicator different from the hard limit indicator module 3310 if the hard limit indicator module 3310 exists or executes within the software 2600. For example, the soft limit alert indicator can be a different color, display a different message, or send a different alert to a remote medical care provider.

Following the soft limit indicator module 3314, operational flow proceeds to the display module 3304 to update the display and allow additional user modifications of the pump settings. Upon termination of operation of the medical infusion pump, operational flow terminates at the end module 3316.

Referring now to FIGS. 34-35, an exemplary user interface 3400 for customizing pump parameters is shown. The user interface 3400 corresponds to the settings module 2808 of FIG. 28, and operates generally as described in FIG. 33. The user interface 3400 includes a status indicator 3402, a continuous rate region, 3404, a demand dose region 3406, a demand dose lockout region 3408, a timed limit

region 3410, and a reservoir region 3412. The user interface 3400 further includes control buttons 3414a-3414c.

The regions 3404-3412 correspond to the patient specific pump parameters 512a of FIG. 5, and can include one or more of continuous rate, demand dose, demand dose lockout, timed limits, reservoir volume, and other patient-specific parameters. Only those patient specific pump parameters that are associated with the selected index of therapy, qualifier, and drug appear in the user interface 3400, and can be as few as one parameter and can incorporate as many parameters as are programmable within a medical infusion pump.

The user interface 3400 presents a standardized interface to a patient or caregiver using the user software 2600, such as at the point of care of a patient or other location on the infusion pump network. The user interface 2600 corresponds to any of a number of types of user software 2600 and administrative software 700. The user interface 3400 also can be configured to be used with various types of medical infusion pumps 102. Most pumps require programming with both patient specific pump parameters and non-patient specific pump parameters, but vary as to the data structure in which this data is passed to the pump. The user interface 3400 reflects from the user software 2600 the number of regions to create corresponding to the number of patient specific pump parameters, which are generic to various types of pumps. Therefore, the user interface 3400 can be used with any of a number of pumps having various brands, interfaces, data structures, or other variances.

The status indicator 3402 displays the library, therapy, qualifier, and drug which define the protocol loaded by the user software 2600. In the exemplary user interface, the library is an "ICU" library, the therapy selected is "Patient Controlled Analgesia", the qualifier is "Adult and Child over 5", and the drug is "Fentanyl 10 mcg/ml".

The continuous rate region 3404 defines the continuous, or basal, rate of drug delivery in the specific medical infusion pump. The continuous rate region 3404 includes a numerical reading 3416 and a meter, shown as a slider bar 3418 and an indicator 3419. The meter generally has two or more locations, each corresponding to a parameter value that can be programmed in the medical infusion pump. Generally, the positional relationship of the meter indicates the setting of the

meter. The numerical reading 3416 reflects the current value for the continuous drug delivery rate.

In the embodiment of the meter shown as the slider bar 3418, the indicator 3419 slides along the slider bar 3418, and the positional relationship
5 between the slider bar 3418 and indicator 3419 dictates the continuous drug delivery rate. Threshold indicators 3420 determine the safe limits within which the continuous rate can be set, and can represent either the hard limits or the soft limits set for the parameter by the administrative software 700. In the embodiment of the user software 2600 whose absolute threshold levels are limited by the hard limits set
10 in the administrative software 700, the ends of each bar 3418 represent the hard limits and the threshold indicators 3420 represent the soft limits for the continuous rate. The thresholds are tested using the method described in FIG. 33.

The demand dose region 3406 customizes the demand dose, or bolus, delivered by the medical infusion pump. The demand dose region includes a
15 numerical reading 3422 and a meter, shown as a slider bar 3424 and indicator 3425. The slider bar 3424 includes threshold indicators 3426 displaying either the hard or soft limit defined in the administrative software. The numerical reading 3422, slider bar 3424, indicator 3425, and threshold indicators 3426 operate analogously to those in the continuous rate region 3404, but set the bolus level parameter rather than the
20 continuous rate parameter.

The demand dose lockout region 3408 customizes the time period after a bolus is delivered in which no additional bolus can be provided. The demand dose lockout region 3408 includes a numerical reading 3428 and a meter, shown as a slider bar 3430 and indicator 3431. The slider bar 3430 includes threshold
25 indicators 3432 displaying either the hard or soft limit defined in the administrative software. The numerical reading 3428, slider bar 3430, indicator 3431, and threshold indicators 3432 operate analogously to those in the continuous rate region 3404, but set the demand dose lockout period parameter rather than the continuous rate parameter.

30 The timed limit region 3410 customizes the amount of the selected drug deliverable by a medical infusion pump within a specified timeframe. The timed limit region 3410 also includes a numerical reading 3434 and a meter, shown as a slider bar 3436 and indicator 3437. The slider bar 3436 includes threshold

indicators 3438 displaying either the hard or soft limit defined in the administrative software. The numerical reading 3434, slider bar 3436, indicator 3437, and threshold indicators 3438 operate analogously to those in the continuous rate region 3404, but set the timed drug delivery threshold parameter rather than the continuous rate parameter.

The reservoir region 3412 defines the size of the reservoir used in conjunction with the medical infusion pump. The size of the reservoir is relevant to computing drug delivery volumes for the purpose of setting alarms and other indicators for replacing or refilling the reservoir. The reservoir region 3412, like the other regions, includes a numerical reading 3440 and a meter, shown as a slider bar 3442 and indicator 3443. The slider bar 3442 includes threshold indicators 3444 displaying either the hard or soft limit defined in the administrative software. The numerical reading 3440, slider bar 3442, and indicator 3443 operate analogously to those in the continuous rate region 3404, but set the reservoir volume parameter rather than the continuous rate parameter. In the embodiment shown, no threshold indicators are included in the reservoir region. This is because the reservoir region 3412 is allowed to use the entire operational range of the reservoir, since no hard limits are set in the reservoir volume region 2016 of FIG. 20 in the administrative software. However, in additional embodiments, one or more threshold indicators can be incorporated into the reservoir region to trigger a warning when the drug reservoir associated with the medical infusion pump contains less than the volume of drugs set by the threshold volume.

The meters in each region 3404-3412 may be adjustable in that each of the patient specific pump parameters are adjustable using a meter. The administrative software 700 enables the adjustment of one or more of the meters by allowing adjustment of those patient specific pump parameters using the options displayed on the user interface 1800 of FIGS. 18-20.

The control buttons 3414a-3414c allow a user to send the currently set parameters to the associated medical infusion pump, cancel the parameter customization, or receive help in the process.

FIG. 35 shows the user interface 3400 wherein the indicator 3437 in the timed limit region 3410 resides along the slider bar 3436 at a location outside the range defined by the threshold indicators 3438. A colored region 3502 appears

around the slider bar and indicator, providing a visual warning to the user of the user software 2600 that abnormal or unadvisable pump settings exist. The colored region 3502 can change color (i.e. green or gray to yellow or red) indicating that the value is outside a threshold level.

5 It is noted that additional screen coloration or textual messages can be used to graphically send messages to a user or programmer of the medical infusion pump or associated computing system. For example, a color code system can be used to reflect a variety of conditions of the medical infusion pump. For example, a color could represent the current coding of the hospital or other health care facility at
10 which the pump may be located. Additionally, the color code can represent a warning condition, a location at which the pump is used, a drug being administered by the pump, or an alert condition. Of course, the color code could represent additional characteristics of the medical infusion pump as well.

 The color code can display on the computing system associated with
15 the medical infusion pump, or can be reflected on a monitor associated to the medical infusion pump itself. Text messages can be sent from the server to be displayed on the monitor of the pump or computing system, such as warnings regarding medication, usage tips for the medical infusion pump, or other medical advice. Additionally, the color code can be placed behind images displayed on the
20 pump which can also represent a region of the hospital, an image of the drug being administered, or other background images. Additionally, the screen coloration described can be represented as a flashing screen, a color changing (cross-fading or otherwise) screen, or various other color patterns.

 FIG. 36 shows a pump programming user interface 3600 for guiding
25 a user through the process of sending a pump program to a medical infusion pump. The programming screen can display properties of a pump to be programmed, and can include the settings 3602a-3602e that are to be sent to the pump, a protocol indicator field 3604, and a confirmation button 3606.

 The settings 3602a-3602e reflect the customized pump parameters set
30 using the user interface 3600 of FIGS. 34-35. The customized pump parameters correspond to the patient specific pump parameters 512a of FIG. 5, and can include one or more of continuous rate, demand dose, demand dose lockout, timed limits, reservoir volume, and other patient-specific parameters.

The protocol indicator field 3604 displays the current protocol selected, in this case shown as "Patient Controlled Analgesia", "Adult and Child over 5", and "Fentanyl 10 mcg/ml". The confirmation button 3606 sends the pump program, including the settings 3602a-3602e for the pump parameters, to the pump.

5 Once the pump is programmed, the pump program executes according to the protocol selected in conjunction with the customized parameters. Referring now to FIG. 37, an exemplary software process 3700 for displaying medical infusion pump customizations is shown. The process occurs within the pump 102, computing system 104, or a combination thereof, and can be part of a
10 pump program sent to a medical infusion pump.

 Customizations in protocol programming refer to differences between the actual pump operation and a standing order (i.e. settings programmed by an administrative user). The standing order can be an original pump parameter or initial value, and the pump operation for comparison can be either a current pump
15 parameter or simply a non-original pump parameter.

 A start module 3702 initiates the process 3700. Operational flow proceeds to a comparison module, shown as a compare pump settings to standing order module 3704. The pump settings stored on the pump or computing system as shown above in FIG. 5 can be compared against a standing order stored on a server
20 as a pump protocol. To ensure that the correct standing order is accessed for comparison, a GUID assigned to the loaded pump settings corresponds to the protocol stored on the server. A display change bar module 3706 presents an indicator on either the medical infusion pump or associated computing system. In the display change bar module, the original pump parameter, or portion of a standing
25 order, can be juxtaposed against the non-original pump parameter, such as in a table format. One or both of the original and non-original pump parameters can be a standard or customized pump parameter. Pump parameters that are displayed in the display change bar module 3706 include drug delivery rate, drug capacity, remaining capacity of the medical infusion pump, bolus levels or occurrences, alarm
30 occurrences, threshold, and/or frequency, drug delivery periods, or other parameters.

 In one possible embodiment, a legend can indicate the meaning of the pump parameters being displayed, and a time/date stamp can display the time at which the original and/or non-original pump parameter was measured. In additional

embodiments, differences between the original and non-original pump parameters can be highlighted when displayed, such as by using a color change or other indicator.

Operational flow terminates at an end module 3708.

5 FIG. 38 is an exemplary schematic illustration of a change bar 3800 displayed on a medical infusion pump 102. The change bar 3800 as shown is displayed on the screen 406 of the pump described above in FIG. 4. Alternately, the change bar 3800 can be displayed on a computing system 104 associated with the medical infusion pump 102.

10 The change bar 3800 includes a plurality of change bar entries 3802. The change bar entries correspond to pump parameters, and in the figure shown the entries 3802a-3802b correspond to patient specific pump parameters. The change bar can display non patient specific pump parameters as well.

 The change bar 3800 can compare any of a number of original and
15 non-original pump parameters. In one embodiment, the change bar 3800 compares the current operation of the pump to the originally programmed operation of the pump. In another embodiment, the change bar 3800 compares the operation of the pump as initially programmed to the suggested programming of the pump based on the original pump protocol. In a further embodiment, the change bar 3800 compares
20 historical activity of the pump to the current pump protocol.

 In one embodiment of the change bar 3800, the change bar entries 3802a-3802b change color when the difference between the original and the non-original pump parameters exceeds a threshold amount. The threshold amount can be, for example, the soft limits set in the administrative software 700. In a further
25 embodiment, the text can change color when a difference greater than the threshold is detected. The change bar can incorporate additional graphics and images on the display.

 Referring now to FIG. 39, user software 2600 is again described in which the user interface 2900 of FIGS. 29-32 is shown with the tasks tab 2940
30 selected. The tasks tab 2940 includes a tasks region 3902 for selection of one or more task options, of which a comparison is displayed when the user selects the "continue" option 3904 shown. Tasks include operations related to maintenance of the pump once in operation. The tasks region 3902 includes a number of pump

comparison options, such as a compare pump settings option 3906, a compare pump settings to protocol option 3908, and a change reservoir option 3910. The compare pump settings option 3906 compares the pump settings to the original protocol from which the pump parameters were based. The compare pump settings to protocol
5 option 3908 compares all pump settings for the protocol selected.

The user software 2600 accesses the protocol loaded on the server 206 to compare the current pump settings to the original or current protocol using the options 3904, 3906. To accomplish this, it is necessary for the user software 2600 to clarify to the server 206 which protocol is being compared within the
10 database 504 of FIG. 5. The user software 2600, in conjunction with the server 206 uses the globally unique identifier (GUID) described above in FIG. 5 to provide the identifier for corresponding the protocol on the pump 102 to the protocol as stored in the server 206. The GUID can be generated by the server 206 and transmitted alongside the protocol and/or library when transmitted to the computing system 104
15 or infusion pump 102, as described above.

The change reservoir option 3910 guides a user of the software 2600 through changing a drug reservoir used in conjunction with the medical infusion pump.

FIG. 40 shows the user interface 2900 with the reports tab 2960
20 selected. During or after pump operation, reports including drug delivery or event logs are used to detect the condition of a patient or of the medical infusion pump. The events which can be tracked using reports tab 2960 are those which are available due to being automatically tracked by the medical infusion pump. The reports tab 2960 presents a number of selectable options for generating reports of
25 pump activity, including time/date information, drug delivery information, and other event information. The reports tab 2960 includes source fields 4002 and option regions 4004. The source fields 4002 present a variety of sources from which reports can be drawn. The source fields 4002 can include the medical infusion pump and stored reports saved within the network. The option regions 4004 present a
30 number of options related to the selected source. For example, a report generated directly from a medical infusion pump can be produced based on prescription settings, and event log, a patient history, drug delivery, or a reported pain scale. A report generated from a saved report can be produced by indicating the patient

identification, the pump identification, or the report type and date. A view field 4006, upon selection by a user, generates the report based on the source and options selected.

FIG. 41 shows a report user interface 4100 for displaying operation of a medical infusion pump. The report user interface 4100 shows the report generated using the options selected in the report tab 2960. The report shown in the report screen is a drug delivery report, and can be printed, saved, or discarded by the user. The drug delivery report can include the volume of the drug delivered, as well as the timing of delivery of the drug. Additional attributes of the medical infusion pump can be reported in the drug delivery report as well.

Aspects of the invention described as being carried out by a computing system or otherwise described as a method of control or manipulation of data may be implemented in one or a combination of hardware, firmware, and software. Embodiments of the invention may also be implemented as instructions stored on a machine-readable medium, which may be read and executed by at least one processor to perform the operations described herein. A machine-readable medium may include any mechanism for storing or transmitting information in a form readable by a machine (e.g., a computer). For example, a machine-readable medium may include read-only memory (ROM), random-access memory (RAM), magnetic disc storage media, optical storage media, flash-memory devices, electrical, optical, acoustical or other form of propagated signals (e.g., carrier waves, infrared signals, digital signals, etc.), and others.

In the foregoing detailed description, various features are occasionally grouped together in a single embodiment for the purpose of streamlining the disclosure. This method of disclosure is not to be interpreted as reflecting an intention that the claimed embodiments of the subject matter require more features than are expressly recited in each claim. Rather, as the following claims reflect, inventive subject matter lies in less than all features of a single disclosed embodiment. Thus, the following claims are hereby incorporated into the detailed description, with each claim standing on its own as a separate preferred embodiment. Therefore, the spirit and scope of the appended claims should not be limited to the description of the preferred versions contained herein.

The Claimed Invention Is:

1. A method of programming a medical infusion pump with a protocol,
the method comprising:
5 displaying a first menu containing a plurality of therapy menu items, each
therapy menu item corresponding to a therapy;
upon selecting a therapy menu item from the first menu, displaying a second
menu containing at least one qualifier menu item, each qualifier menu
item corresponding to a qualifier;
10 upon selecting a qualifier menu item from the second menu, displaying a
third menu containing at least one drug menu item, each drug menu
item corresponding to a drug.
2. The method of claim 1, further comprising:
selecting a drug menu item; and
15 identifying a set of pump data corresponding to the selected therapy menu
item, qualifier menu item, and drug menu item.
3. The method of claim 2, wherein the set of pump data includes the
protocol.
4. The method of claim 2, further comprising:
20 selecting a program corresponding to the set of pump data corresponding to
the selected therapy menu item, qualifier menu item, and drug menu
item.
5. The method of claim 4, further comprising uploading the protocol to
the medical infusion pump, the protocol including the therapy, the qualifier, and the
25 drug.
6. The method of claim 4, wherein the program includes therapy-
specific parameters.

7. The method of claim 4, wherein the program includes a drug delivery rate parameter.
8. The method of claim 4, wherein the program includes a bolus parameter.
- 5 9. The method of claim 4, wherein the protocol includes patient-specific pump parameters.
- 10 10. The method of claim 4, wherein the protocol includes a patient-specific drug delivery rate parameter.
- 10 11. The method of claim 4, wherein the protocol includes a patient-specific bolus parameter.
12. The method of claim 1, wherein the list of qualifiers includes age.
13. The method of claim 1, wherein the list of qualifiers includes weight.
14. The method of claim 1, further comprising displaying notes associated with the protocol.
- 15 15. The method of claim 1, further comprising:
upon selecting a therapy menu item from the first menu, displaying notes associated with the therapy.
- 20 16. The method of claim 1, further comprising:
upon selecting a qualifier menu item from the second menu,
displaying notes associated with the qualifier.
17. The method of claim 1, further comprising:
upon selecting a drug menu item from the third menu, displaying notes associated with the drug.

18. The method of claim 1, further comprising displaying a drug delivery report.

19. The method of claim 18, wherein the drug delivery report displays the volume and timing of delivery of the drug.

5 20. An apparatus for programming a medical infusion pump with a protocol, the apparatus comprising:
a memory configured to store a plurality of therapy menu items, qualifier menu items, and drug menu items;
a programmable circuit in electrical communication with the memory, the
10 programmable circuit programmed to generate a user interface configured to:
display a first menu containing a plurality of therapy menu items, each therapy menu item corresponding to a therapy;
upon selecting a therapy menu item from the first menu, display a
15 second menu containing at least one qualifier menu item, each qualifier menu item corresponding to a qualifier;
upon selecting a qualifier menu item from the second menu, display a
third menu containing at least one drug menu item, each drug menu item corresponding to a drug.

20 21. The apparatus of claim 20, wherein the graphical user interface is further configured to:
select a drug menu item from the third menu; and
identify a set of pump data corresponding to the selected therapy menu item, qualifier menu item, and drug menu item.

25 22. The apparatus of claim 21, wherein the graphical user interface is further configured to:
select a program corresponding to the set of pump data corresponding to the selected therapy menu item, qualifier menu item, and drug menu item.

23. The apparatus of claim 22, further comprising a monitor in electrical communication with the programmable circuit, the monitor configured to display the user interface.

5 24. The apparatus of claim 20, wherein the qualifier menu items include age.

25. The apparatus of claim 20, wherein the qualifier menu items include weight.

26. The apparatus of claim 20, wherein the drug menu items include at least one drug associated with each therapy.

10 27. The apparatus of claim 20, wherein the programmable circuit is further programmed to display a drug delivery report.

28. The method of claim 27, wherein the drug delivery report includes the volume and timing of delivery of the drug.

15 29. An apparatus for programming a medical infusion pump with a protocol, the apparatus comprising:
a memory configured to store a plurality of therapy menu items, qualifier menu items, and drug menu items;
a display;
a programmable circuit in electrical communication with the memory and the
20 display, the programmable circuit programmed to generate a user interface configured to:
display a first menu containing a plurality of therapy menu items, each therapy menu item corresponding to a therapy administrable by the medical infusion pump;
25 upon selecting a therapy menu item from the first menu, display a second menu containing at least one qualifier menu item, each qualifier menu item corresponding to a qualifier associated with the therapy;

- upon selecting a qualifier menu item from the second menu, display a third menu containing at least one drug menu item, each drug menu item corresponding to a drug associated with the therapy and the qualifier;
- 5 upon selecting a drug menu item from the third menu, identify a set of pump data corresponding to the selected therapy menu item, qualifier menu item, and drug menu item;
- select a program corresponding to the set of pump data; and
upload the program to the medical infusion pump.
- 10 30. An apparatus for programming a medical infusion pump with a protocol, the apparatus comprising:
- a memory configured to store a plurality of therapy menu items, qualifier menu items, and drug menu items;
- a programmable circuit in electrical communication with the memory, the
- 15 programmable circuit programmed to generate a graphical user interface
- having therapy, qualifier, and drug window objects, the graphical user interface configured to;
- present a list of therapy menu items in the therapy window object,
- 20 each therapy menu item corresponding to a therapy;
- upon selecting a therapy menu item in the therapy window object,
- present a list of qualifier menu items in the qualifier window object, each qualifier menu item corresponding to a qualifier;
- upon selecting a qualifier menu item in the qualifier window object,
- 25 present a list of drug menu items in the drug window object, each drug menu item corresponding to a drug.
31. The apparatus of claim 30, wherein the programmable circuit is further programmed to generate a graphical user interface configured to:
- select a drug menu item in the drug window object; and
- 30 identify a set of pump data corresponding to the selected therapy menu item, qualifier menu item, and drug menu item.

32. The apparatus of claim 31, wherein the programmable circuit is further programmed to:

select a program corresponding to the set of pump data corresponding to the selected therapy menu item, qualifier menu item, and drug menu item.

33. The apparatus of claim 31, wherein the memory is further configured to store therapy notes, qualifier notes, and drug notes.

34. The apparatus of claim 33, wherein the programmable circuit is further programmed to generate a graphical user interface having therapy notes, qualifier notes, and drug notes window objects, the graphical user interface configured to;

upon selecting a therapy menu item in the therapy window object, present therapy notes in the therapy notes window object;

upon selecting a qualifier menu item in the qualifier window object, present qualifier notes in the qualifier notes window object;

upon selecting a drug menu item in the drug window object, present drug notes in the drug notes window object.

35. The apparatus of claim 34, wherein the therapy notes include information related to application of the therapy.

36. The apparatus of claim 34, wherein the qualifier notes include information related to modification of the therapy by the qualifier.

37. The apparatus of claim 34, wherein the drug notes include information related to administration of the drug.

38. A method of programming a medical infusion pump with a protocol, the method comprising:

presenting on a monitor a list of therapy menu items in a therapy window object, each therapy menu item corresponding to a therapy;

upon selecting a therapy menu item in the therapy window object, presenting
on the monitor a list of qualifier menu items in a qualifier window
object, each qualifier menu item corresponding to a qualifier;
upon selecting a qualifier menu item in a qualifier window object, presenting
5 on the monitor a list of drug menu items in the drug window object,
each drug menu item corresponding to a drug.

39. The method of claim 38, further comprising:
selecting a drug menu item; and
identifying a set of pump data corresponding to the selected therapy menu
10 item, qualifier menu item, and drug menu item.

40. The method of claim 38, further comprising:
select a program corresponding to the set of pump data corresponding to the
selected therapy menu item, qualifier menu item, and drug menu
item.

41. The method of claim 38, further comprising:
upon selecting the therapy menu item in the therapy window object,
presenting on the monitor therapy notes in a therapy notes window
object;
upon selecting the qualifier menu item in the qualifier window object,
20 presenting on the monitor qualifier notes in a qualifier notes window
object;
upon selecting the drug menu item in the drug window object, presenting on
the monitor drug notes in a drug notes window object.

42. A method of programming a medical infusion pump with a protocol,
25 the method comprising:
displaying a first menu containing a plurality of first index items;
upon selecting a first index item, displaying a second menu containing at
least one second index item; and
upon selecting a second index item from the second menu, displaying a third
30 menu containing at least one third index item;

wherein the first index item, the second index item, and the third index item in combination refer to a protocol for programming into the medical infusion pump.

43. The method of claim 42, wherein the plurality of first index items are
5 selected from the group of indices consisting of:
therapies;
qualifiers;
drugs;
pump programs;
10 caregivers; and
locations.

44. The method of claim 42, wherein the plurality of second index items
are selected from the group of indices consisting of:
therapies;
15 qualifiers;
drugs;
pump programs;
caregivers; and
locations.

20 45. The method of claim 42, wherein the plurality of third index items are
selected from the group of indices consisting of:
therapies;
qualifiers;
drugs;
25 pump programs;
caregivers; and
locations.

46. The method of claim 42, wherein the protocol includes patient
specific pump parameters.

47. The method of claim 42, wherein the protocol includes non patient specific pump parameters.

48. An apparatus for programming a medical infusion pump with a protocol, the apparatus comprising:

5 a memory configured to store a plurality of first menu items, second menu items, and third menu items;

a display;

a programmable circuit in electrical communication with the memory and the display, the programmable circuit programmed to:

10 display a first menu containing a plurality of first index items;

upon selecting a first index item, display a second menu containing at least one second index item; and

upon selecting a second index item from the second menu, display a third menu containing at least one third index item;

15 wherein the first index item, the second index item, and the third index item in combination refer to a protocol for programming into the medical infusion pump.

49. The apparatus of claim 48, wherein the protocol includes patient specific pump parameters.

20 50. The apparatus of claim 48, wherein the plurality of first index items are selected from the group of indices consisting of:

therapies;

qualifiers;

drugs;

25 pump programs;

caregivers; and

locations.

51. The apparatus of claim 48, wherein the plurality of second index items are selected from the group of indices consisting of:

30 therapies;

5 qualifiers;
 drugs;
 pump programs;
 caregivers; and
 locations.

10 52. The apparatus of claim 48, wherein the plurality of third index items
 are selected from the group of indices consisting of:
 therapies;
 qualifiers;
 drugs;
 pump programs;
 caregivers; and
 locations.

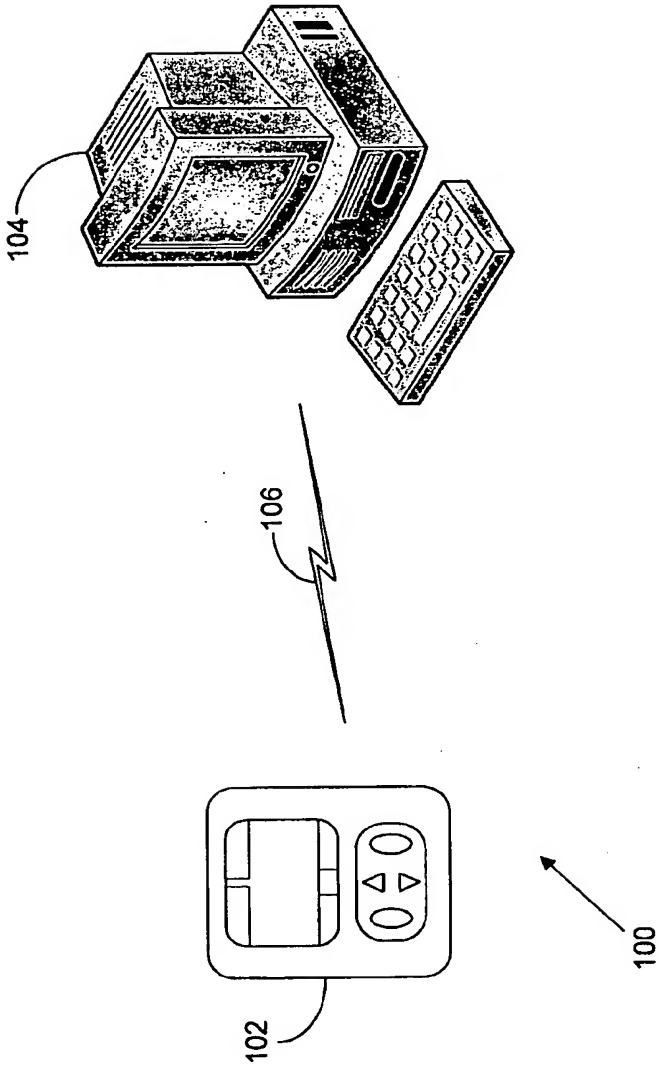


FIG. 1

FIG. 2

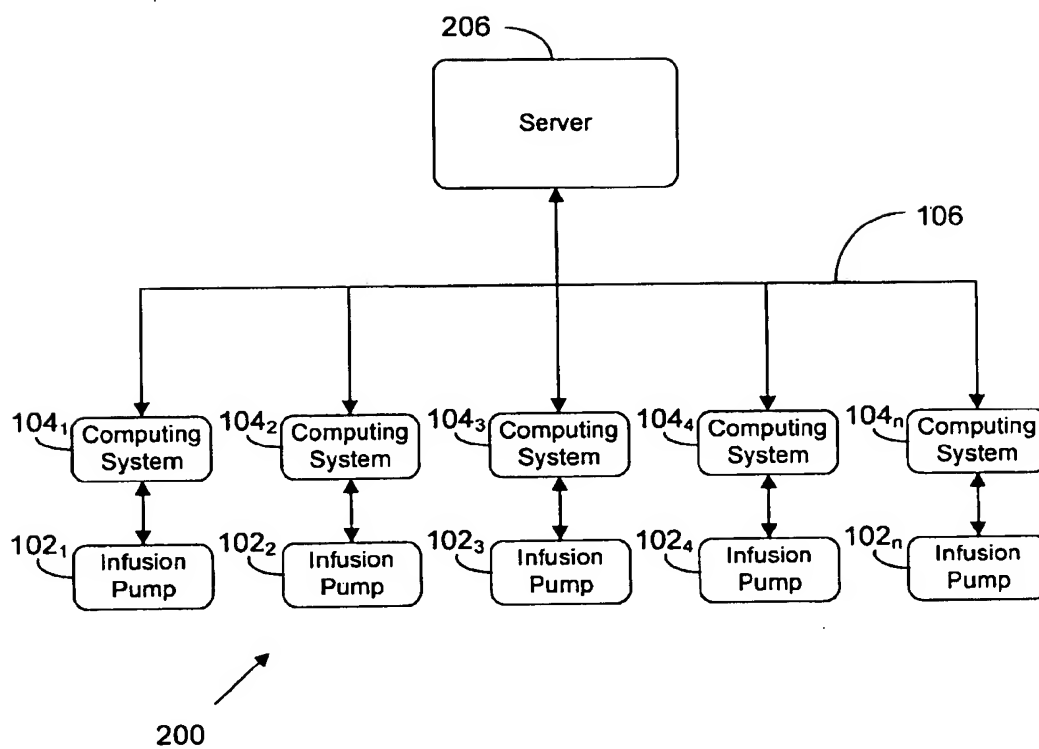


FIG. 3

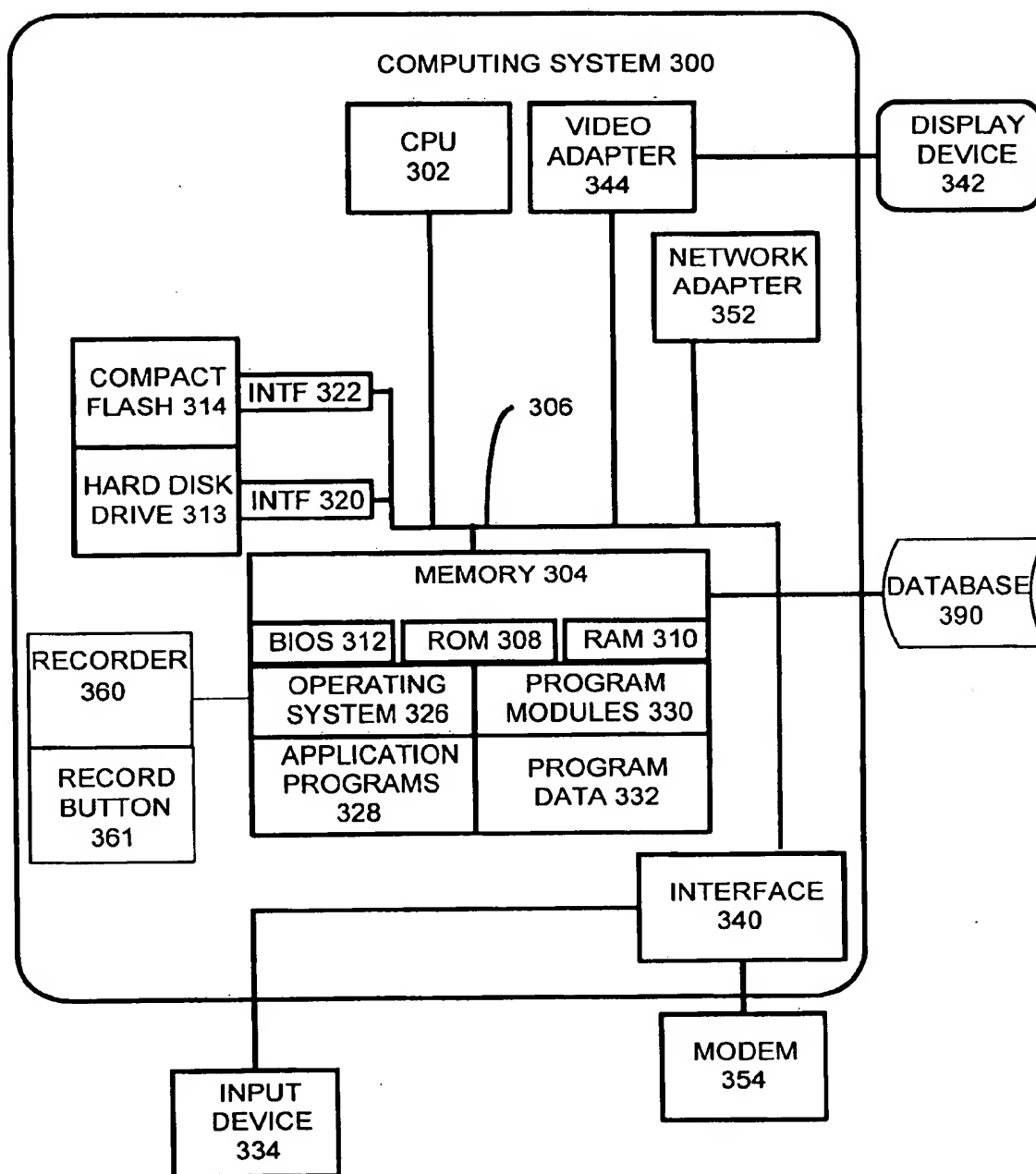


FIG. 4

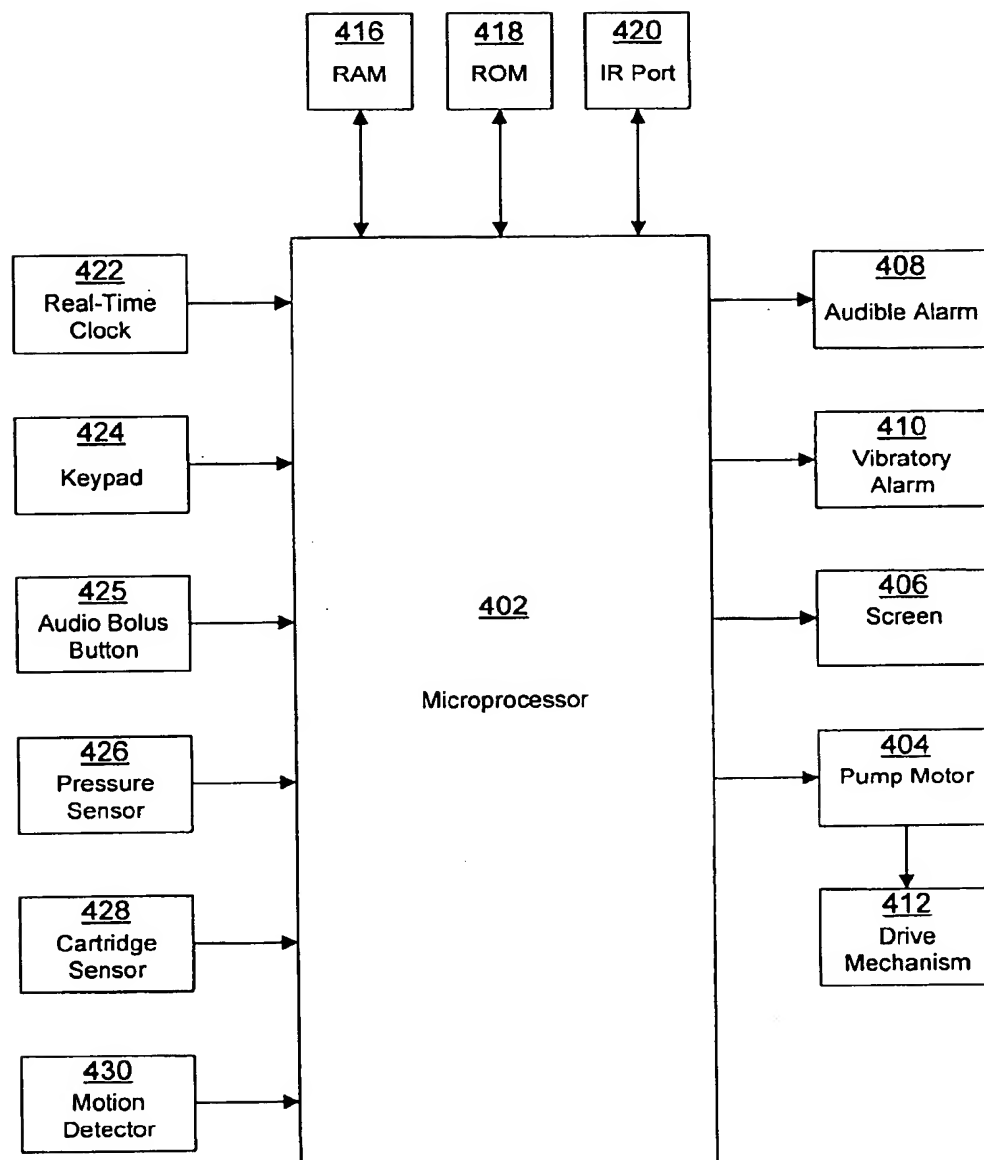
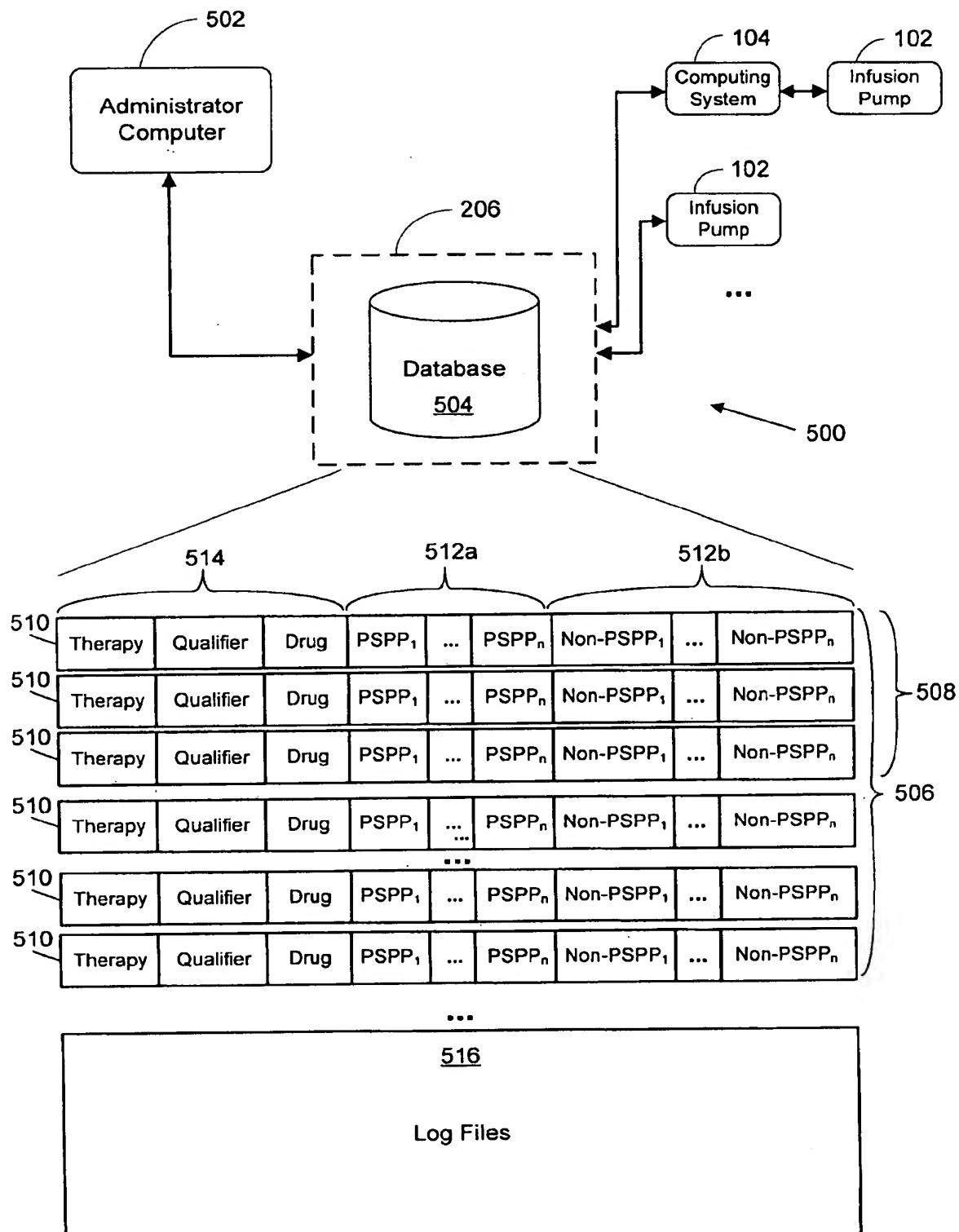


FIG. 5

6/41

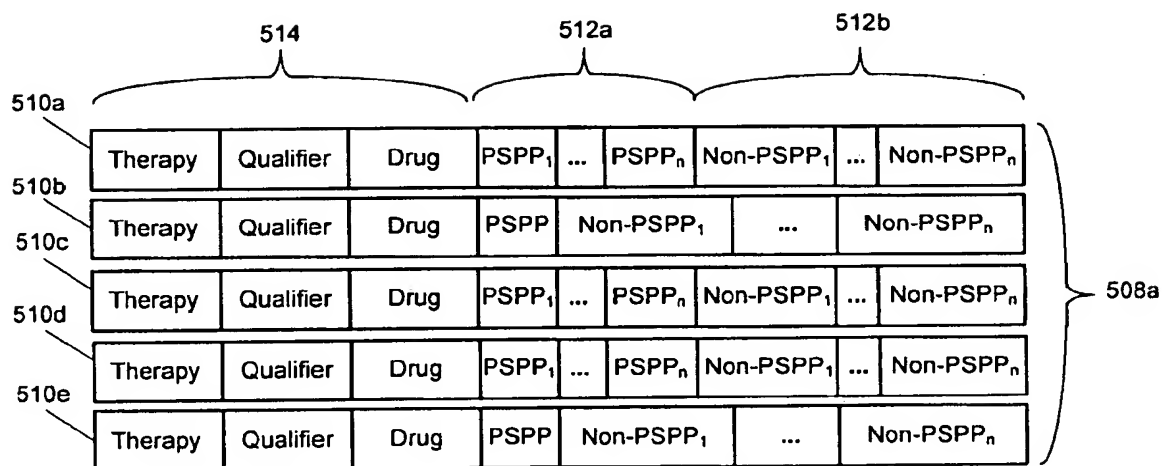


FIG. 6A

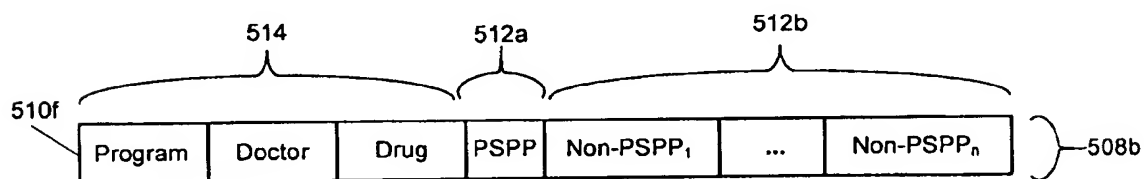


FIG. 6B

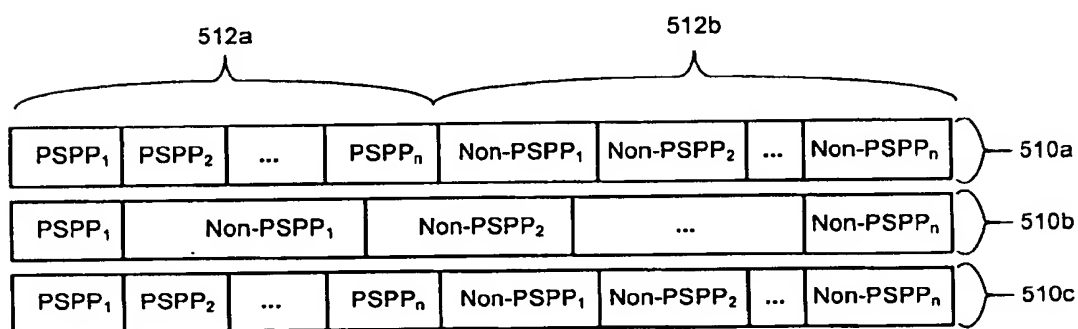


FIG. 6C

7/41

FIG. 7

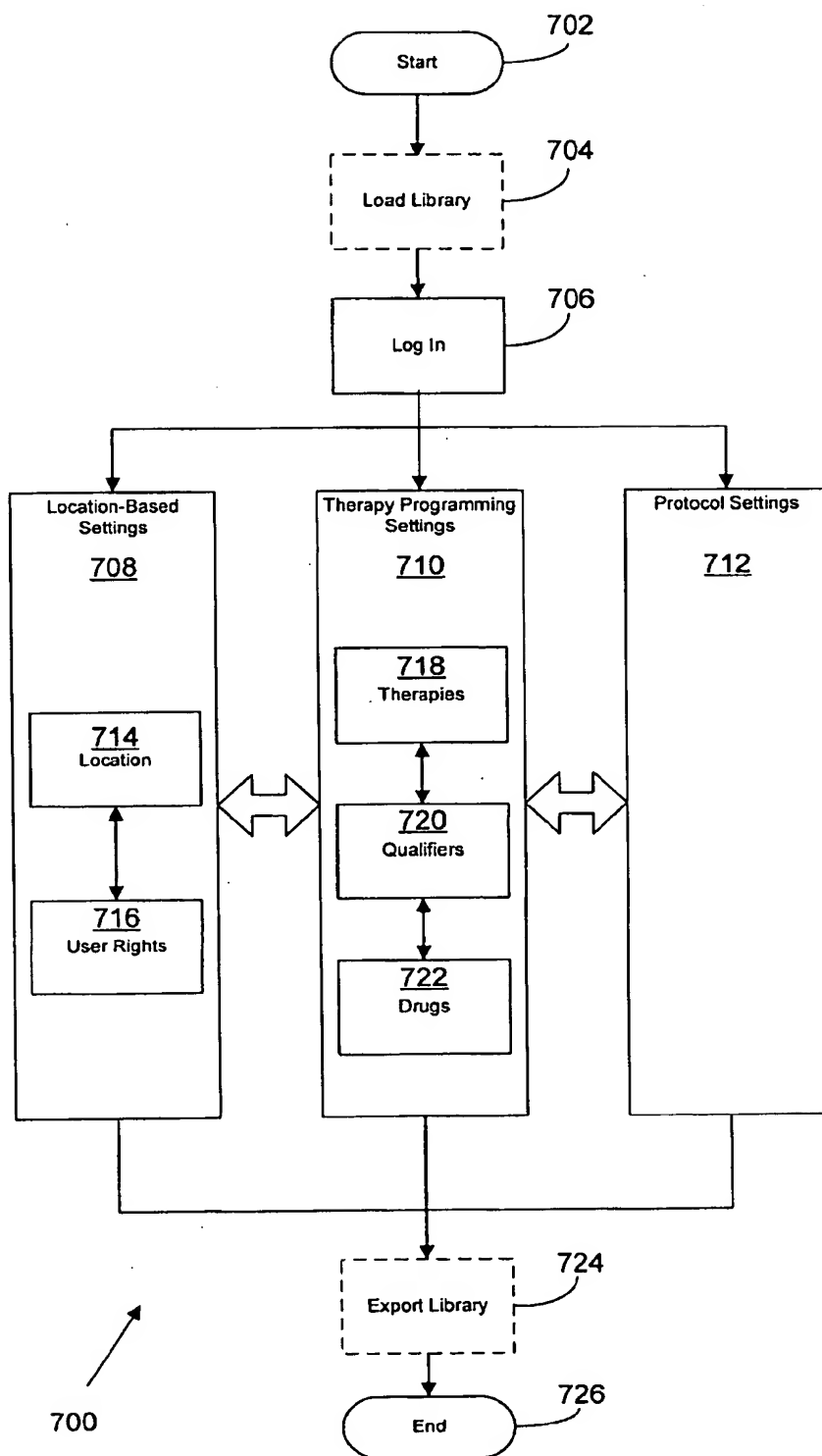


FIG. 8

800

Import Library

Import from Location: 808

Browse 810

Select Protocol Libraries:

- ☐ Neonatal Intensive Care Unit (NICU)
- ☐ Post Anesthesia Care Unit (PACU)
- ☐ Intensive Care Unit (ICU)
- ☐ Our Library Base Template
- ☐ Library One
- ☐ Library Two

802

Select All 806a Clear All 806b

Enter Import Database Password

804

OK Cancel Help

FIG. 9

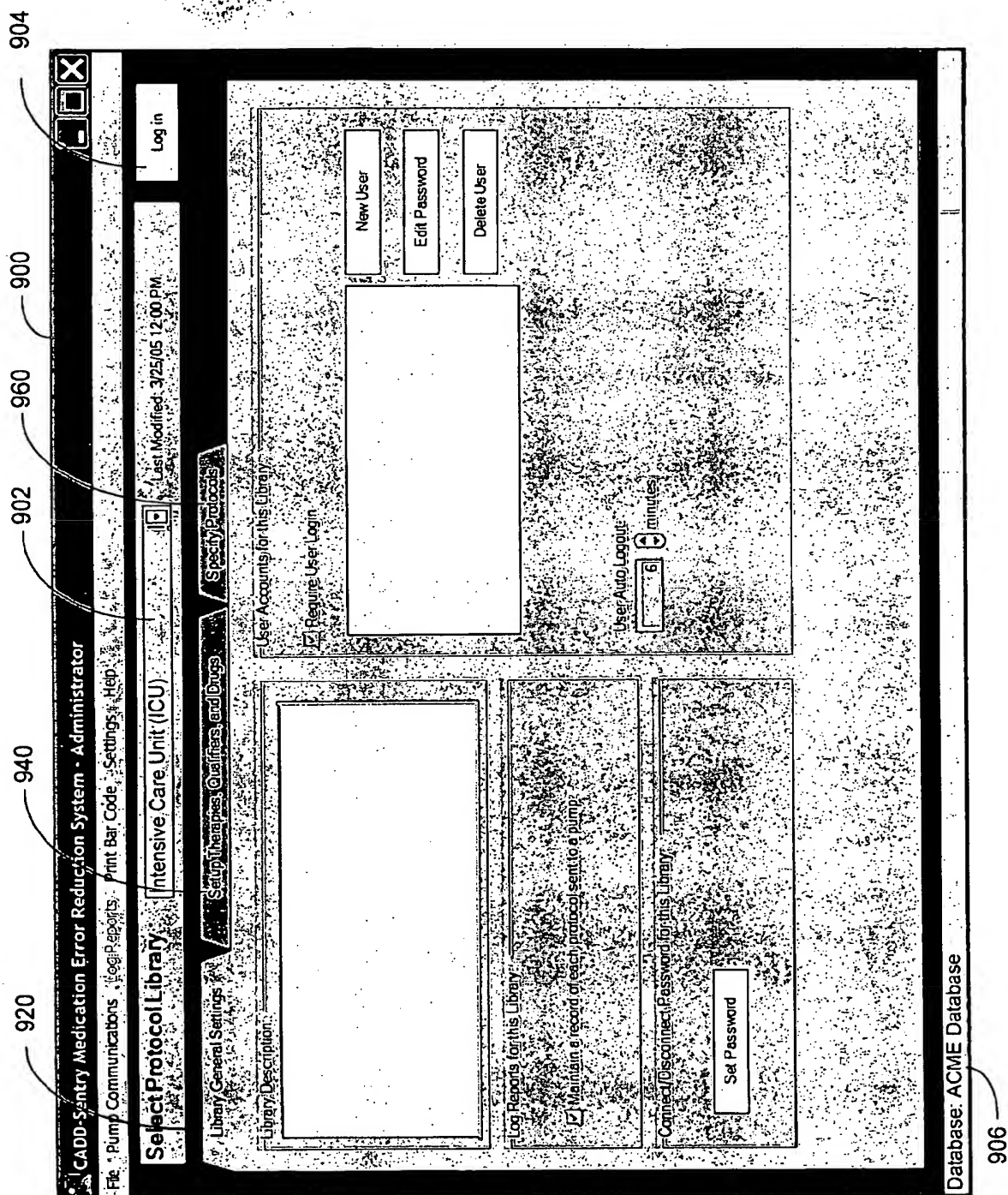


FIG. 10

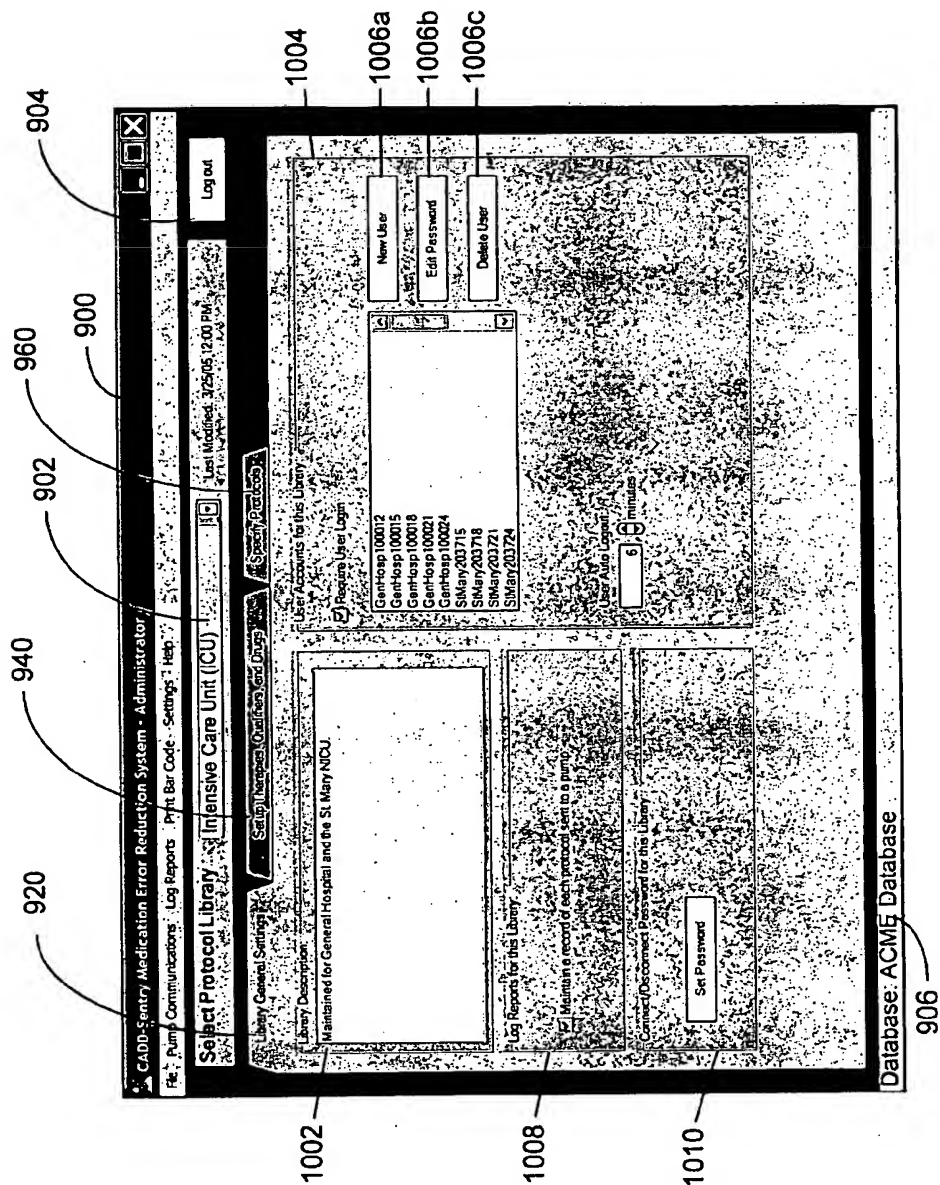


FIG. 11


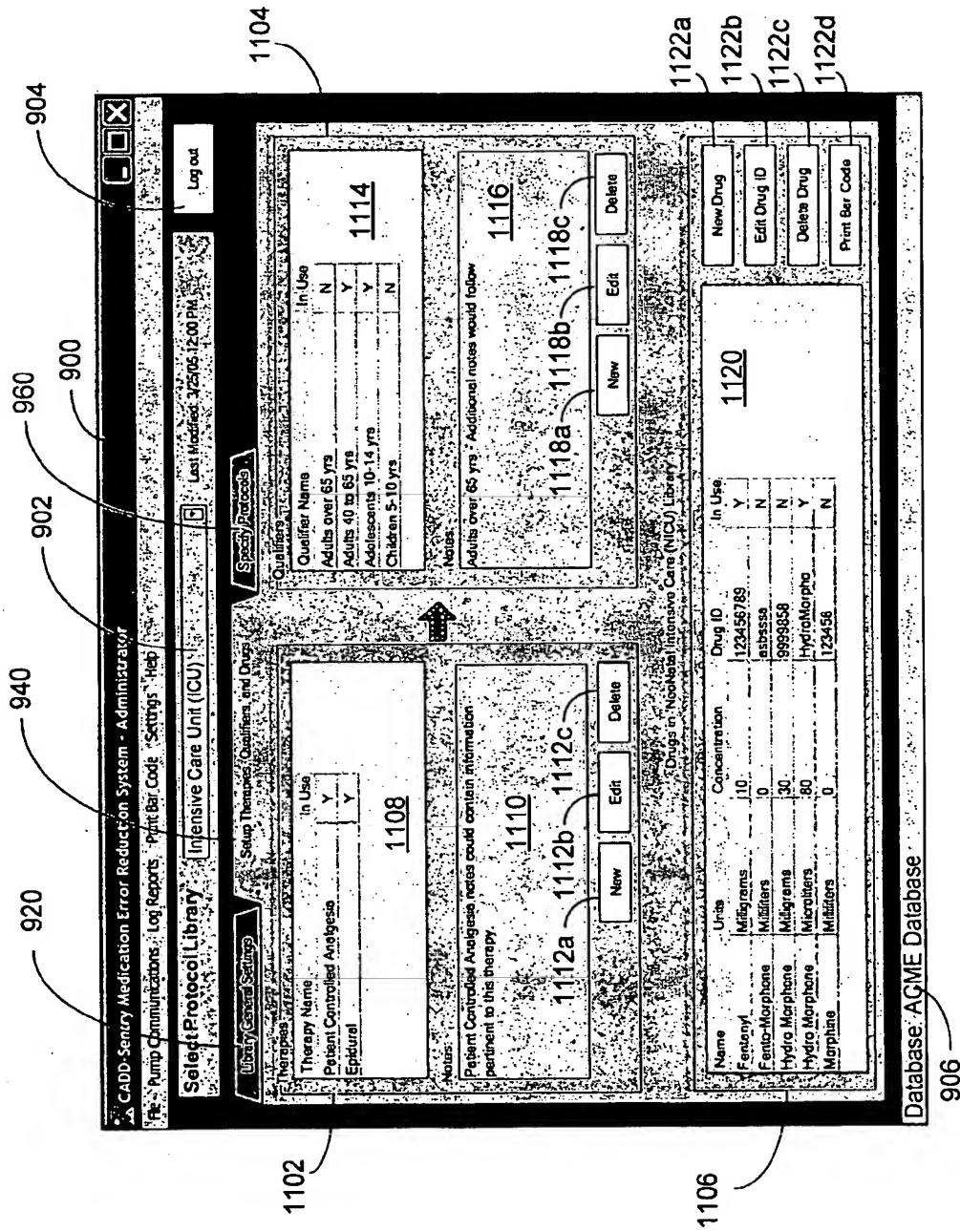


FIG. 12

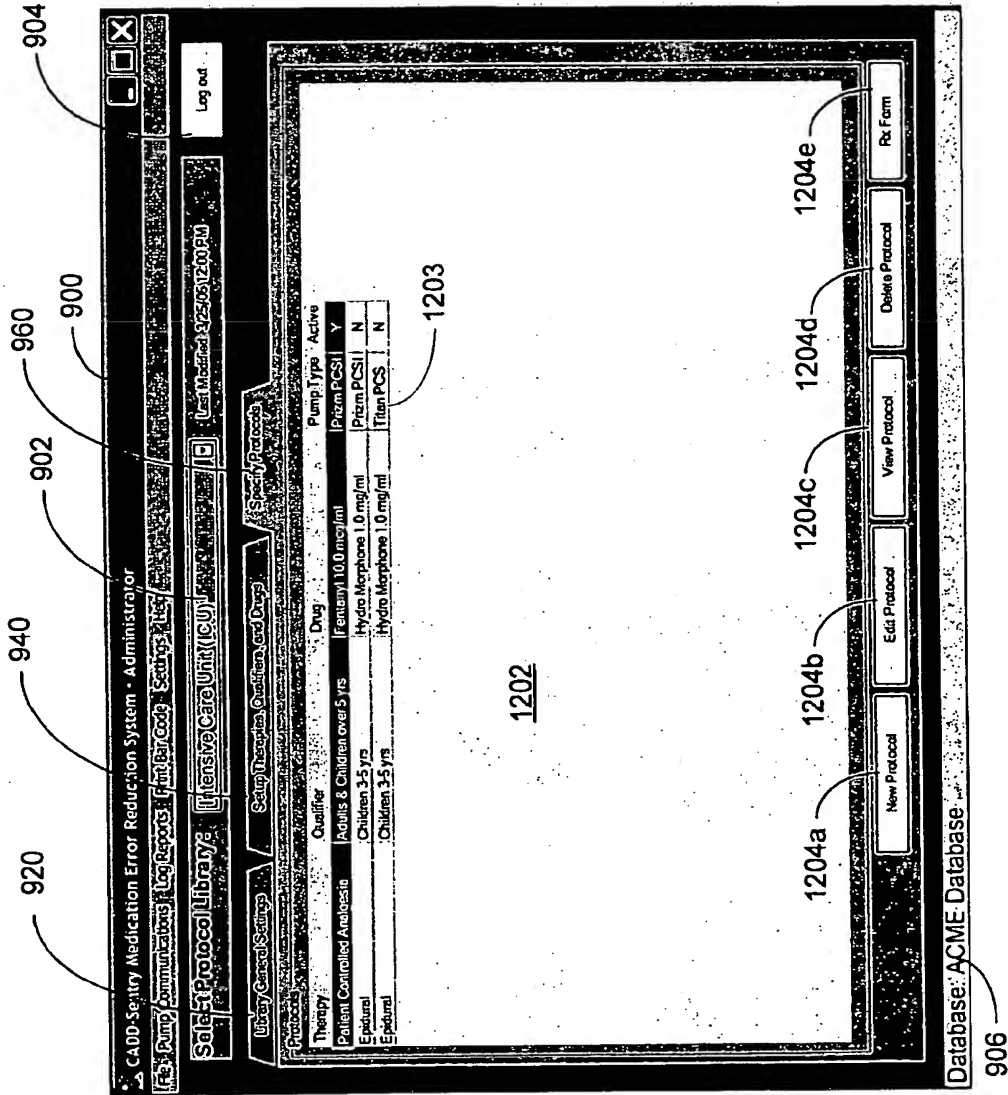

13/41

FIG. 13

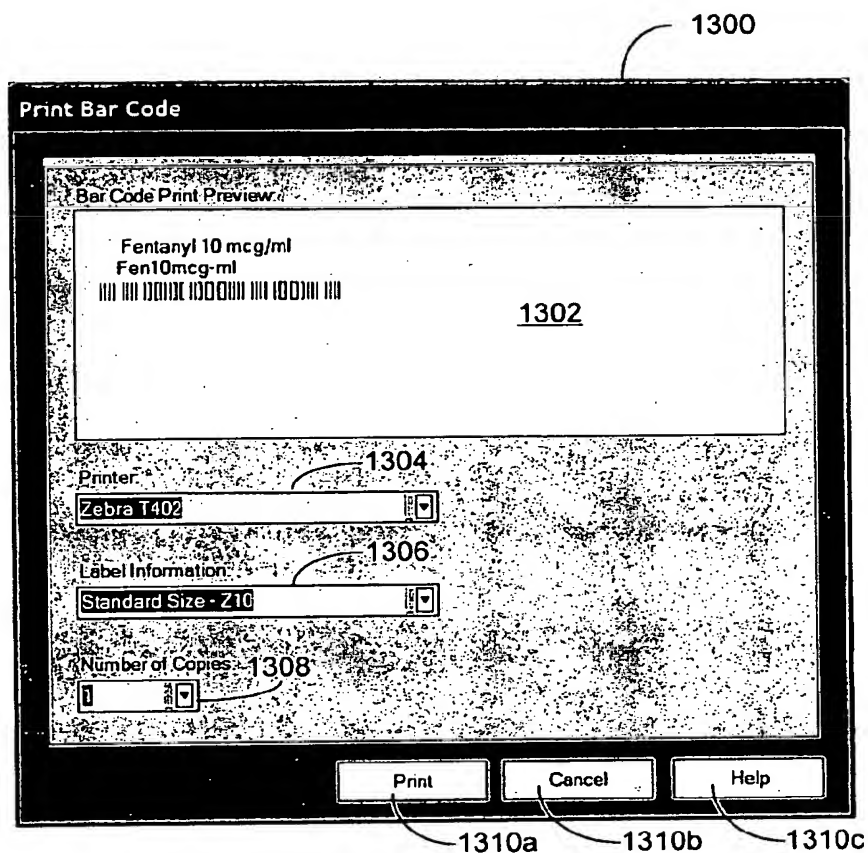


FIG. 14

1400

View Prescription Order Form

77%

10

1

1

1

1

1

1

1

Close

Print

Prescription Form:

Patient Controlled Analgesia
Adults & Children over 5 yrs

Date: _____
Time: _____

Therapy:	Patient Controlled Analgesia
Notes:	Patient controlled analgesia setup with morphine 300 and 5000 mg 3.1245.
Qualifier:	Adults & Children over 5 yrs
Notes:	Monitor and record respiratory rate and level of consciousness every two (2) hours times twelve (12) hours upon PCA initiation or dose change, then every four (4) hours. Record pain location and intensity every four (4) hours

Fentanyl 100.0 mcg/ml	
Concentration Rate mcg/hr	60 to 50 mcg/hr
Demand Dose mcg	5 to 20 mcg
Dose Lockout (min)	0.00 to 2.5 hr
4 Hour Lock Dose	60 to 100 mcg
Reservoir Volume ml	***
Prescribed Notes: Monitor for side effects such as: dizziness, decreased arterial and pulse, vomiting.	

Hydrocodone 1.0 mcg/ml	
Concentration Rate mcg/hr	60 to 50 mcg/hr
Demand Dose mcg	5 to 20 mcg
Dose Lockout (min)	0.00 to 2.5 hr
4 Hour Lock Dose	60 to 100 mcg
Reservoir Volume ml	***
Prescribed Notes: Consider reducing dose gradually. Monitor for side effects such as: dizziness, decreased arterial and pulse, decreased respiratory rate, nausea and vomiting.	

Morphine 1.0 mcg/ml	
Concentration Rate mcg/hr	60 to 50 mcg/hr
Demand Dose mcg	5 to 20 mcg
Dose Lockout (min)	0.00 to 2.5 hr
4 Hour Lock Dose	60 to 100 mcg
Reservoir Volume ml	

FIG. 15

1500

New Protocol

Therapy

1502

Patient Controlled Analgesia

Patient Controlled Analgesia

Epidural

Therapy Notes

Patient controlled analgesia

1504

Qualifier

1506

Drug

1510

Qualifier Notes

1508

Protocol Notes

1512

Next

Cancel

Help

FIG. 16

1500

New Protocol

Therapy	Qualifier	Drug	Protocol Notes
1502 Patient Controlled Analgesia	1506 Adults & Children over 5 yrs Children 3-5 yrs	1510	1512
Therapy Notes	Qualifier Notes		
1504 Patient controlled analgesia	1508		

Next Cancel Help

FIG. 17

New Protocol

Therapy	Qualifier	Drug
1502	1506	1510
Patient Controlled Analgesia Therapy Notes:	Adults & Children over 5 yrs Qualifier Notes:	Fentanyl 10.0 mcg/ml Hydro Morphine 1.0 mg/ml Morphine 1.0 mg/ml Morphine 10.0 mg/ml
Patient controlled analgesia Therapy Notes:	Monitor and record respiratory rate and level of consciousness every two (2) hours times twelve (12) hours upon PCA initiation or dose change, then every four (4) hours. Record pain location and intensity every four (4) hours.	1512
1504	1508	

Next
 Cancel
 Help

FIG. 18

Protocol Details

Protocol Will Be Active in Point of Care

Therapy: Patient Controlled Analgesia

Qualifier: Adults & Children over 5 years

Dilg: Fentanyl 100 mcg/ml

Delivery 1

Delivery 2

Alarms

Security

Display/Sound

Reports

Limits for Patient Specific Parameters

Continuous Rate

838a ☐ Demand Dose Lockout 1826

838b ☐ Display 1824

838c ☐ Allow Edit 1864b

838d ☐ Require Comment 1864c

838e ☐ Hard Maximum 1858

838f ☐ Soft Maximum 1862

838g ☐ Initial Value 1856

838h ☐ Soft Minimum 1862

838i ☐ Hard Minimum 1858

Demand Dose

1852a ☐ Demand Dose Lockout 1826

1852b ☐ Display 1824

1852c ☐ Allow Edit 1864b

1852d ☐ Require Comment 1864c

1852e ☐ Hard Maximum 1858

1852f ☐ Soft Maximum 1862

1852g ☐ Initial Value 1856

1852h ☐ Soft Minimum 1862

1852i ☐ Hard Minimum 1858

General Settings

1816 ☐ Require Drug Verification

1817 ☐ Require Second Nurse Verification

(User Login must be enabled to use this feature.)

1818 ☐ Set as a Weight-Based Protocol

Weight Limit: 0 kg

1806a Save

1806b Cancel

1806c Help

FIG. 19

Protocol Details

☐ Protocol Will Be Active in Point of Care

Therapy: Patient Controlled Analgesia

Qualifier: Adults & Children over 5 years

Drug: Fentanyl 100 mcg/ml

Delivery

☐ Delivery 2

Alarms

Security

Reports

Weight-Based Protocol

☐ Require Drug Verification

☐ Require Second Nurse Verification

(User Login must be enabled to use this feature)

Weight-Based

☒ Set as Weight-Based Protocol

Weight Limit: 100 kg

Limits for Patient Specific Parameters:

Continuous Rate mcg/kg/hr

☐ Display ☐ Allow Edit ☐ Require Comment

Hard Maximum: 0

Soft Maximum: 0

Initial Value: 0

Soft Minimum: 0

Hard Minimum: 0

Demand Dose mcg/kg

☐ Display ☐ Allow Edit ☐ Require Comment

Hard Maximum: 0

Soft Maximum: 0

Initial Value: 0

Soft Minimum: 0

Hard Minimum: 0

Demand Dose Lockout

☐ Display ☐ Allow Edit ☐ Require Comment

Hard Maximum: 00:00

Soft Maximum: 00:00

Initial Value: 00:00

Soft Minimum: 00:00

Hard Minimum: 00:00

Save **Cancel** **Help**

FIG. 20

1800

1802

1860

1850

1840

1830

1820

1810

1804

Protocol Details

☐ Protocol Will Be Active in Point of Care

Therapy: Patient Controlled Analgesia

Qualifier: Adults & Children over 5 years

Drug: Fentanyl 100 mcg/ml

Delivery 1: ☐ Delivery 2: ☐ Alarms: ☐ Security: ☐ Reports: ☐

Limits for Patient Specific Parameters:

Dosing Limit

☒ Delivery Limit

☐ No Limit

4 hours

☐ Doses per Hour

☐ Neither

4 Hour Delivery Limit

2026a ☐ Display

2026b ☐ Allow Edit

2026c ☐ Require Comment

Hard Maximum: 0

Soft Maximum: 0

Initial Value: 0

Soft Minimum: 0

Hard Minimum: 0

Maximum Demand Doses per Hour

2014 ☐ Display

2016 ☐ Allow Edit

2018 ☐ Require Comment

Hard Maximum: 0

Soft Maximum: 0

Initial Value: 0

Soft Minimum: 0

Hard Minimum: 0

Reservoir Volume Limit

2042 ☐ Display

2044 ☐ Allow Edit

2046a ☐ Display

2046b ☐ Allow Edit

Not in Use

Initial Value: 0

Titration

☐ No Titration

☐ Titration Limit

☐ No Limit

2006

Maximum Delivery Rate: 125 ml/hr

Maximum Crinician Bolus: 0 mcg

2008

2010

2002

2004

2012

2019

2020

2024

2028

2030

2032

2034

2036

2038

2040

2042

2044

2046

2048

2050

2052

2054

2056

2058

2060

2062

2064

2066

2068

2070

2072

2074

2076

2078

2080

2082

2084

2086

2088

2090

2092

2094

2096

2098

2100

2102

2104

2106

2108

2110

2112

2114

2116

2118

2120

2122

2124

2126

2128

2130

2132

2134

2136

2138

2140

2142

2144

2146

2148

2150

2152

2154

2156

2158

2160

2162

2164

2166

2168

2170

2172

2174

2176

2178

2180

2182

2184

2186

2188

2190

2192

2194

2196

2198

2200

2202

2204

2206

2208

2210

2212

2214

2216

2218

2220

2222

2224

2226

2228

2230

2232

2234

2236

2238

2240

2242

2244

2246

2248

2250

2252

2254

2256

2258

2260

2262

2264

2266

2268

2270

2272

2274

2276

2278

2280

2282

2284

2286

2288

2290

2292

2294

2296

2298

2300

2302

2304

2306

2308

2310

2312

2314

2316

2318

2320

2322

2324

2326

2328

2330

2332

2334

2336

2338

2340

2342

2344

2346

2348

2350

2352

2354

2356

2358

2360

2362

2364

2366

2368

2370

2372

2374

2376

2378

2380

2382

2384

2386

2388

2390

2392

2394

2396

2398

2400

2402

2404

2406

2408

2410

2412

2414

2416

2418

2420

2422

2424

2426

2428

2430

2432

2434

2436

2438

2440

2442

2444

2446

2448

2450

2452

2454

2456

2458

2460

2462

2464

2466

2468

2470

2472

2474

2476

2478

2480

2482

2484

2486

2488

2490

2492

2494

2496

2498

2500

2502

2504

2506

2508

2510

2512

2514

2516

2518

2520

2522

2524

2526

2528

2530

2532

2534

2536

2538

2540

2542

2544

2546

2548

2550

2552

2554

2556

2558

2560

2562

2564

2566

2568

2570

2572

2574

2576

2578

2580

2582

2584

2586

2588

2590

2592

2594

2596

2598

2600

2602

2604

2606

2608

2610

2612

2614

2616

2618

2620

2622

2624

2626

2628

2630

2632

2634

2636

2638

2640

2642

2644

2646

2648

2650

2652

2654

2656

2658

2660

2662

2664

2666

2668

2670

2672

2674

2676

2678

2680

2682

2684

2686

2688

2690

2692

2694

2696

2698

2700

2702

2704

2706

2708

2710

2712

2714

2716

2718

2720

2722

2724

2726

2728

2730

2732

2734

2736

2738

2740

2742

2744

2746

2748

2750

2752

2754

2756

2758

2760

2762

2764

2766

2768

2770

2772

2774

2776

2778

2780

2782

2784

2786

2788

2790

2792

2794

2796

2798

2800

2802

2804

2806

2808

2810

2812

2814

2816

2818

2820

2822

2824

2826

2828

2830

2832

2834

2836

2838

2840

2842

2844

2846

2848

2850

2852

2854

2856

2858

2860

2862

2864

2866

2868

2870

2872

2874

2876

2878

2880

2882

2884

2886

2888

2890

2892

2894

2896

2898

2900

2902

2904

2906

2908

2910

2912

2914

2916

2918

2920

2922

2924

2926

2928

2930

2932

2934

2936

2938

2940

2942

2944

2946

2948

2950

2952

2954

2956

2958

2960

2962

2964

2966

2968

2970

2972

2974

2976

2978

2980

2982

2984

2986

2988

2990

2992

2994

2996

2998

3000

3002

3004

3006

3008

3010

3012

3014

3016

3018

3020

3022

3024

3026

3028

3030

3032

3034

3036

3038

3040

3042

3044

3046

3048

3050

3052

3054

3056

3058

3060

3062

3064

3066

3068

3070

3072

3074

3076

3078

3080

3082

3084

3086

3088

3090

3092

3094

3096

3098

3100

3102

3104

3106

3108

3110

3112

3114

3116

3118

3120

3122

3124

3126

3128

3130

3132

3134

3136

3138

3140

3142

3144

3146

3148

3150

3152

3154

3156

3158

3160

3162

3164

3166

3168

3170

3172

3174

3176

3178

3180

3182

3184

3186

3188

3190

3192

3194

3196

3198

3200

3202

3204

3206

3208

3210

3212

3214

3216

3218

3220

3222

3224

3226

3228

3230

3232

3234

3236

3238

3240

3242

3244

3246

3248

3250

3252

3254

3256

3258

3260

3262

3264

3266

3268

3270

3272

3274

3276

3278

3280

3282

3284

3286

3288

3290

3292

3294

3296

3298

3300

3302

3304

3306

3308

3310

3312

3314

3316

3318

3320

3322

3324

3326

3328

3330

3332

3334

3336

3338

3340

3342

3344

3346

3348

3350

3352

3354

3356

3358

3360

3362

3364

3366

3368

3370

3372

3374

3376

3378

3380

3382

3384

3386

3388

3390

3392

3394

3396

3398

3400

3402

3404

3406

3408

3410

3412

3414

3416

3418

3420

3422

3424

3426

3428

3430

3432

3434

3436

3438

3440

3442

3444

3446

3448

3450

3452

3454

3456

3458

3460

3462

3464

3466

3468

3470

3472

3474

3476

3478

3480

3482

3484

3486

3488

3490

3492

3494

3496

3498

3500

3502

3504

3506

3508

3510

3512

3514

3516

3518

3520

3522

3524

3526

3528

3530

3532

3534

3536

3538

3540

3542

3544

3546

3548

3550

3552

3554

3556

3558

3560

3562

3564

3566

3568

3570

3572

3574

3576

3578

3580

3582

3584

3586

3588

3590

3592

3594

3596

3598

3600

3602

3604

3606

3608

3610

3612

3614

3616

3618

3620

3622

3624

3626

3628

3630

3632

3634

3636

3638

3640

3642

3644

3646

3648

3650

3652

3654

3656

3658

3660

3662

3664

3666

3668

3670

3672

3674

3676

3678

3680

3682

3684

3686

3688

3690

3692

3694

3696

3698

3700

3702

3704

3706

3708

3710

3712

3714

3716

3718

3720

3722

3724

3726

3728

3730

3732

3734

3736

3738

3740

3742

3744

3746

3748

3750

3752

3754

3756

3758

3760

3762

3764

3766

3768

3770

3772

3774

3776

3778

3780

3782

3784

3786

3788

3790

3792

3794

3796

3798

3800

3802

3804

3806

3808

3810

3812

3814

3816

3818

3820

3822

3824

3826

3828

3830

3832

3834

3836

3838

3840

3842

3844

3846

3848

3850

3852

3854

3856

3858

3860

3862

3864

3866

3868

3870

3872

3874

3876

3878

3880

3882

3884

3886

3888

3890

3892

3894

3896

3898

3900

3902

3904

3906

3908

3910

3912

3914

3916

3918

3920

3922

3924

3926

3928

3930

3932

3934

3936

3938

3940

3942

3944

3946

3948

3950

3952

3954

3956

3958

3960

3962

3964

3966

3968

3970

3972

3974

3976

3978

3980

3982

3984

3986

3988

3990

3992

3994

3996

3998

4000

4002

4004

4006

4008

4010

4012

4014

4016

4018

4020

4022

4024

4026

4028

4030

4032

4034

4036

4038

4040

4042

4044

4046

4048

4050

4052

4054

4056

4058

4060

4062

4064

4066

4068

4070

4072

4074

4076

4078

4080

4082

4084

4086

4088

4090

4092

4094

4096

4098

4100

4102

4104

4106

4108

4110

4112

4114

4116

4118

4120

4122

4124

4126

4128

4130

4132

4134

4136

4138

4140

4142

4144

4146

4148

4150

4152

4154

4156

4158

4160

4162

4164

4166

4168

4170

4172

4174

4176

4178

4180

4182

4184

4186

4188

4190

4192

4194

4196

4198

4200

4202

4204

4206

4208

4210

4212

4214

4216

4218

4220

4222

4224

4226

4228

4230

4232

4234

4236

4238

4240

4242

4244

4246

4248

4250

4252

4254

4256

4258

4260

4262

4264

4266

4268

4270

4272

4274

4276

4278

4280

4282

4284

4286

4288

4290

4292

4294

4296

4298

FIG. 21

1800

1802

1860

1850

1840

1830

1820

1810

1804

Protocol Details

☐ Protocol Will Be Active in Point of Care

Therapy: Patient Controlled Analgesia

Qualifier: Adults & Children over 5 years

Drug: Fentanyl 100 mcg/ml

Delivery Alarms Security Display/Sound Reports

Reservoir Low Alarm 2102

☐ Standard

☐ Trip Point: 0.6 ml

Reservoir Empty Alarm Type 2104

☐ Single

☐ Insistent

Pump Stopped Alarm 2106

☐ Beeps

☐ Two Tone

Preventive Maintenance 2108

☐ Reminder On

☐ Reset Reminder

Interval: 0.6 months

Due: 03/02/2005

Pump Hardware 2110

☐ Air Detector Required

☐ Air Detector On

☐ Upstream Sensor On

Save 1806a Cancel 1806b Help 1806c

FIG. 22

1800

1802

1860

1850

1840

1830

1820

1810

1804

Protocol Details

Therapy: Patient Controlled Analgesia

Qualifier: Adults & Children over 5 years

Drug: Fentanyl 10.0 mg/ml

Delivery 1

Delivery 2

Alarms

Security

Display/Sound

Reports

Auto Lock: Not in Use

Lock Code: 61

Clinician Code: 997

Initial Lock Level: 2

Epidural Mode: 2204

Epidural Mode On

Save

Cancel

Help

1806a

1806b

1806c

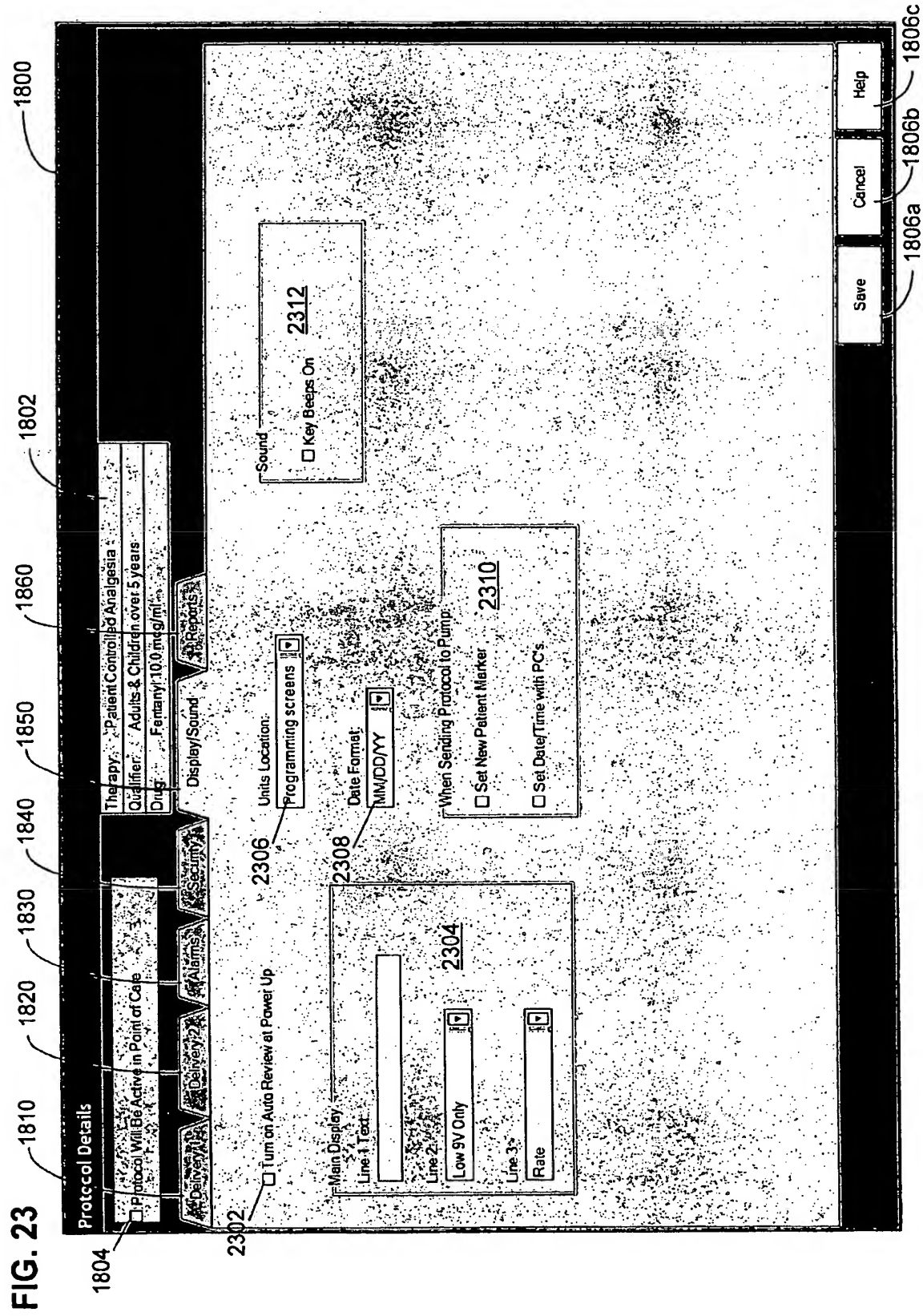


FIG. 24

1800

1802

1860

1850

1840

1830

1820

1810

1804

Protocol Details

☐ Protocol Will Be Active in Point of Care

Therapy: Patient Controlled Analgesia

Qualifier: Adults & Children over 5 years

Drug: Fentanyl 10.0 mg/ml

Delivery

Alarms

Security

Display Sound

Reports

New Patient Marker

2402

☐ Reports/No Clear

☐ Power Up/No Clear

☐ Reports/Clear

☐ Power Up/Clear

Custom Reports

2404

☐ Dose Counters

☐ Given

☐ Doses Hour by Hour

☐ Patient Review

☐ Pain Scale Entry

☐ Pain Scale Log

☐ Delivery Log

☐ Event Log

☐ New Patient Marker

Save

Cancel

Help

1806a

1806b

1806c

25/41

FIG. 25

2500

Export Library

Select Protocol Libraries

☐ Neonatal Intensive Care Unit (NICU) 2502

☐ Post Anesthesia Care Unit (PACU)

☐ Intensive Care Unit (ICU)

☐ Our Library Base Template

☐ Library One

☐ Library Two

Select All Clear All

Export Option

☒ Read Only for Use at Point of Care 2504

☐ Read-Write Template for Use in Administrator


Create New Passwords

(This password will be used when importing the selected libraries.)

OK Cancel Help

26/41

FIG. 26

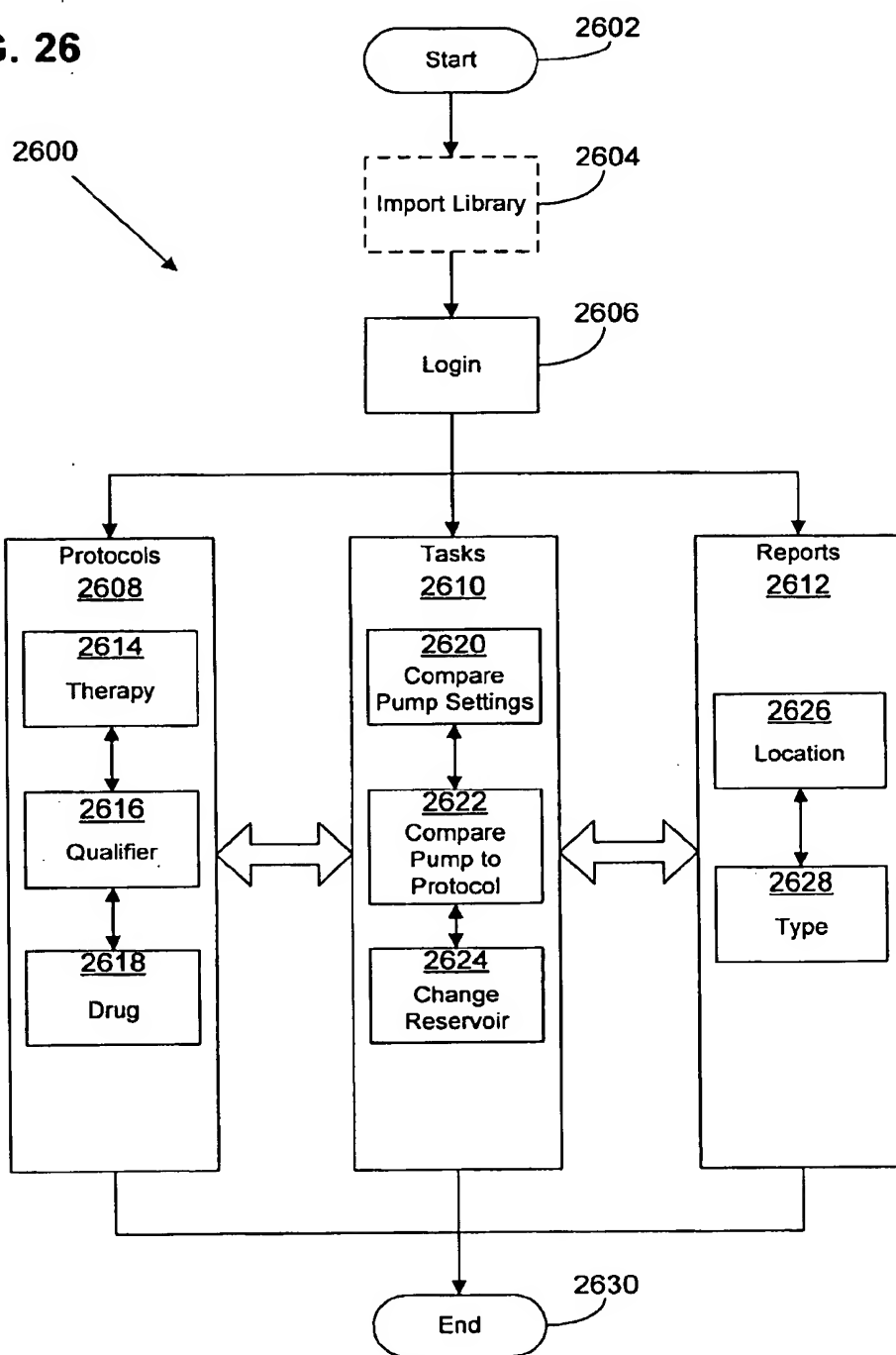

2700

Import Library to Local Computer

Help

2702

Import From Location

ExportFromAdmin.mdb

Browse

Library

2704

Name

Description

Last Modified On

Enter Library Password

2708a

OK

2708b


Cancel


2706

FIG. 27

28/41

FIG. 28

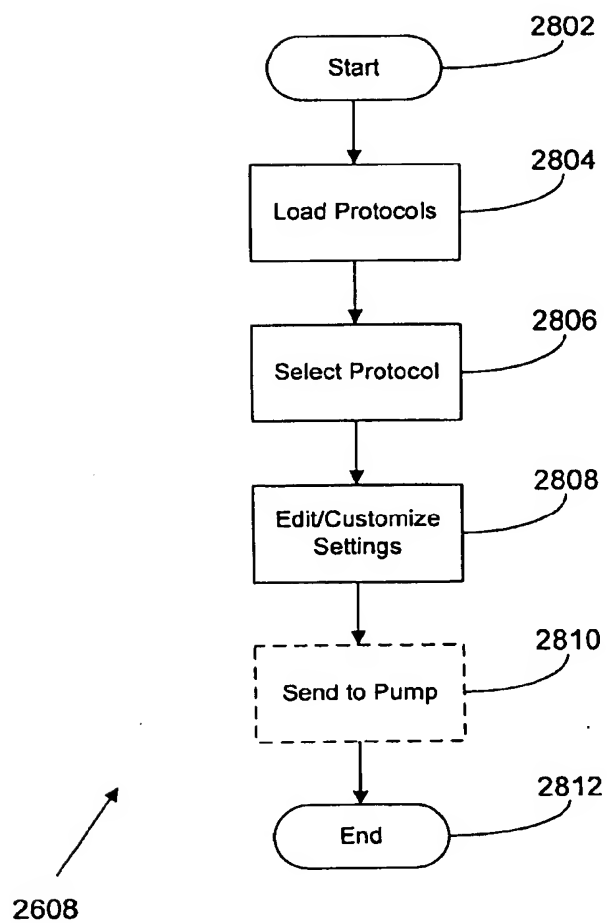

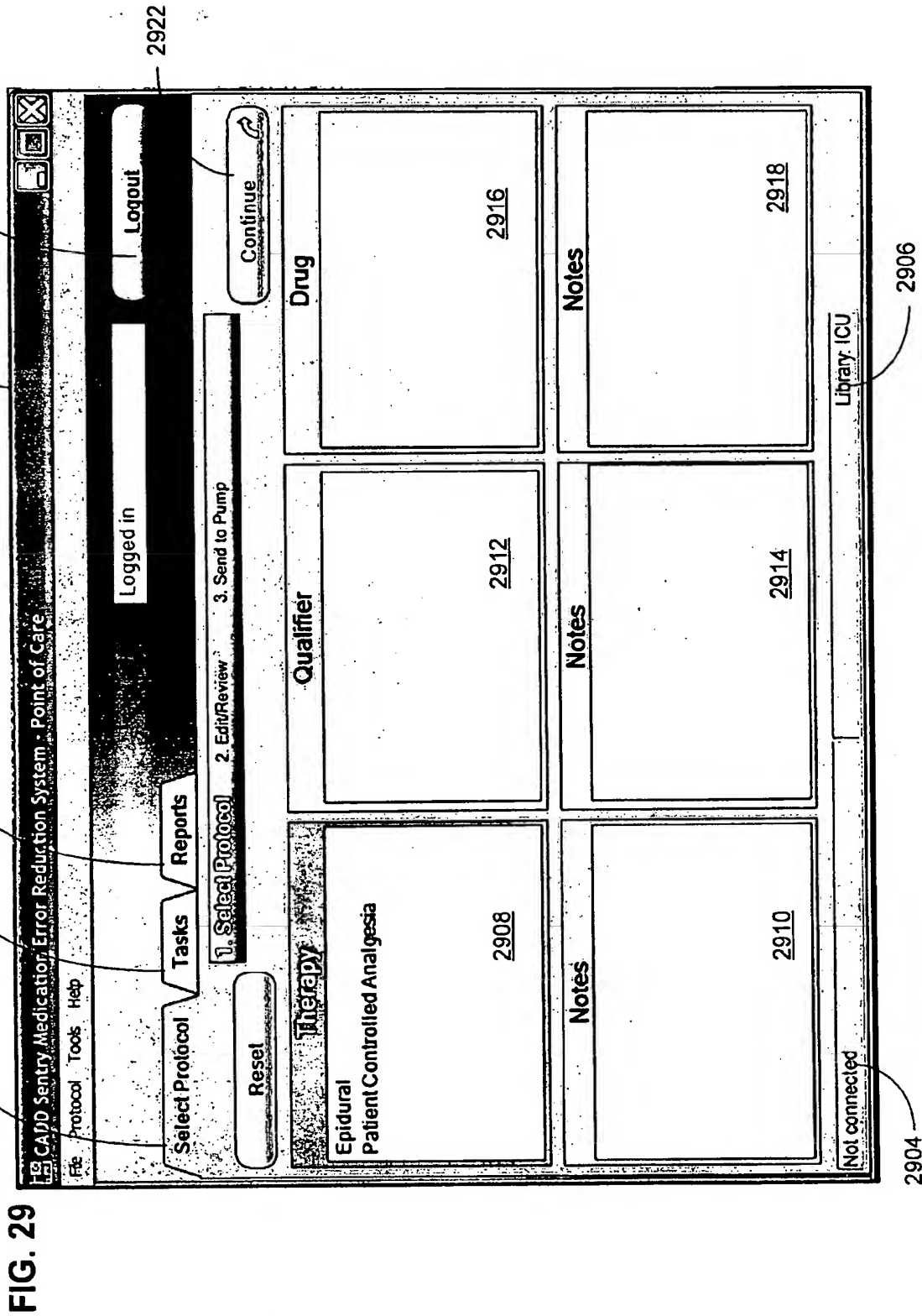

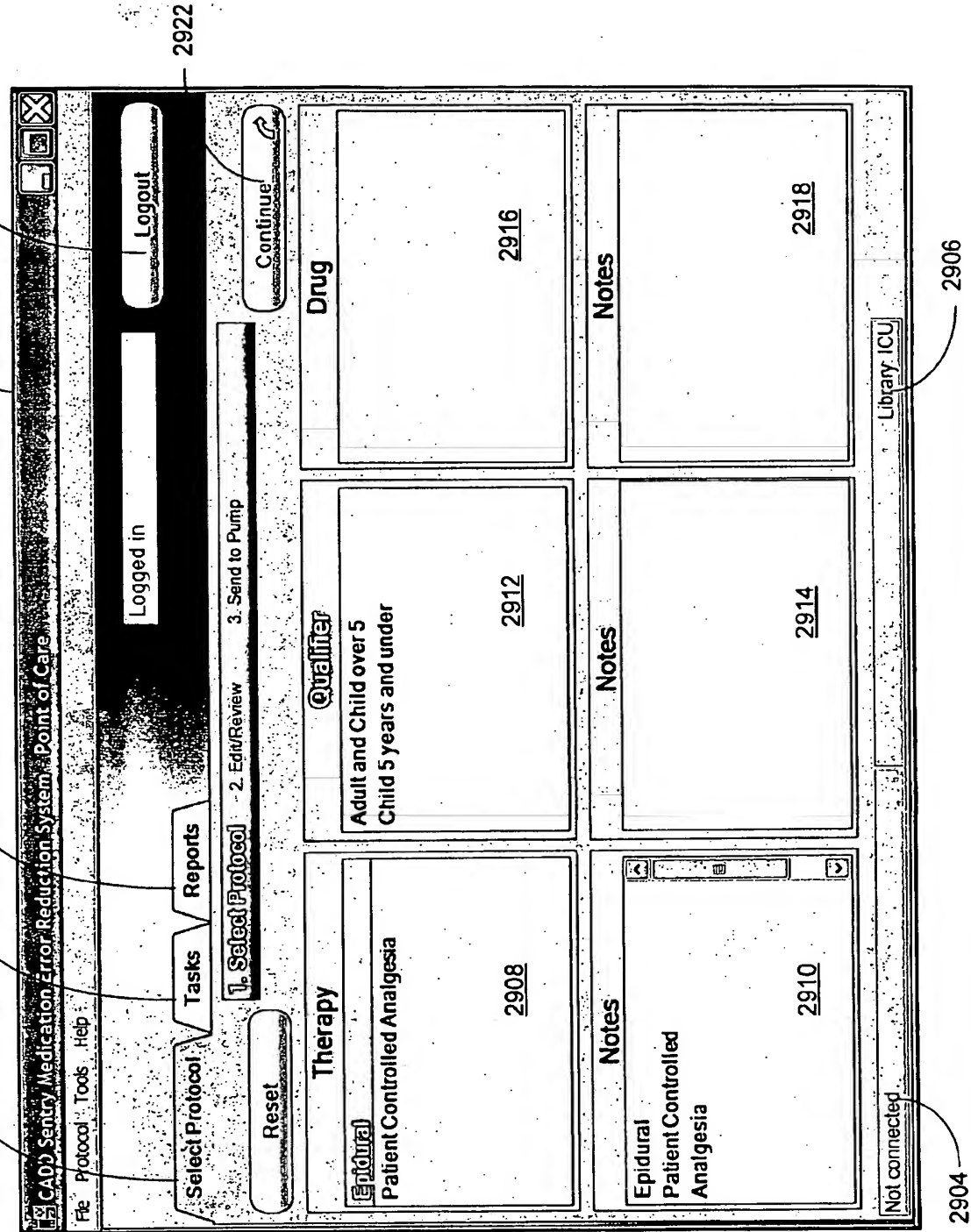


FIG. 30

FIG. 31

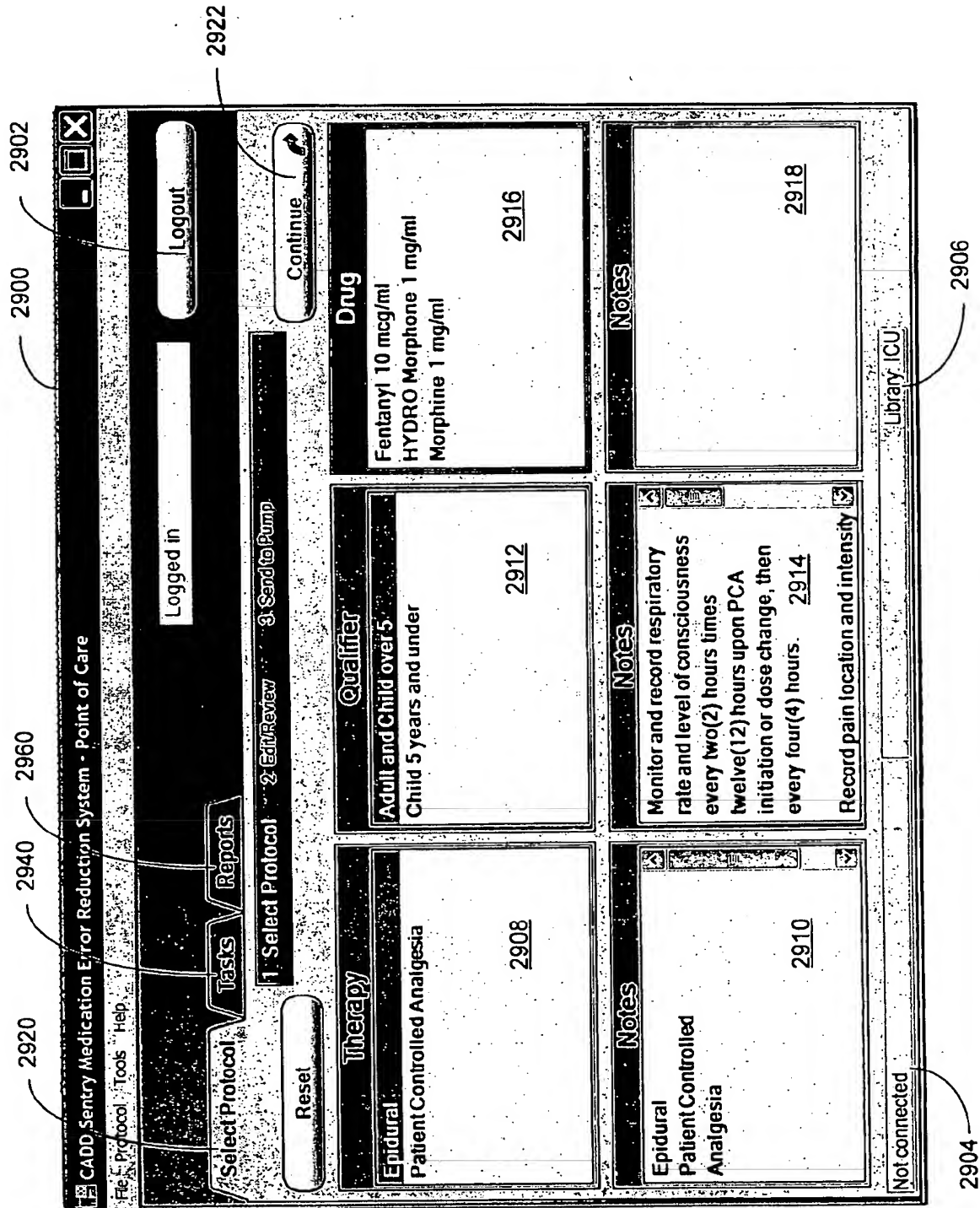


FIG. 32


The screenshot displays the CADD Sentry Medication Error Reduction System interface. At the top, a menu bar includes 'File', 'Protocol', 'Tools', and 'Help'. Below the menu, a status bar shows 'Logged in' and a 'Logout' button. The main interface is divided into several sections:


- Navigation Bar:** Contains buttons for 'Select Protocol', 'Tasks', and 'Reports'. Below these, a progress indicator shows '1. Select Protocol', '2. Edit Review', and '3. Send to Pump'.
- Buttons:** A 'Reset' button is located on the left, and a 'Continue' button with a right-pointing arrow is on the right.
- Therapy Section:** A box labeled 'Therapy' contains the text 'Epidural Patient Controlled Analgesia' and the number '2908'.
- Qualifier Section:** A box labeled 'Qualifier' contains the text 'Adult and Child over 5', 'Adults over 65 years', and 'Child 5 years and under', along with the number '2912'.
- Drug Section:** A box labeled 'Drug' contains the text 'Fentanyl 10 mcg/ml', 'HYDRO Morphine 1 mg/ml', and 'Morphine 1 mg/ml', along with the number '2916'.
- Notes Section:** A box labeled 'Notes' contains the text 'Patient Controlled Analgesia' and the number '2910'.
- Notes Section (Right):** A box labeled 'Notes' contains the text 'Consider reducing dose gradually. Monitor for side effects such as: dizziness, SOB, anxiety, confusion, decreased urine output, hallucinations, constipation, dry mouth, nausea and vomiting. Monitor for symptoms of overdose: cold clammy' and the number '2918'.
- Status Bar:** At the bottom, it shows 'Not connected' and 'Library ICU'.

Reference numerals are placed around the interface: 2920 (top left), 2940 (top center), 2960 (top right), 2900 (middle left), 2902 (middle right), 2922 (bottom left), 2904 (bottom center), and 2906 (bottom right).


33/41

FIG. 33

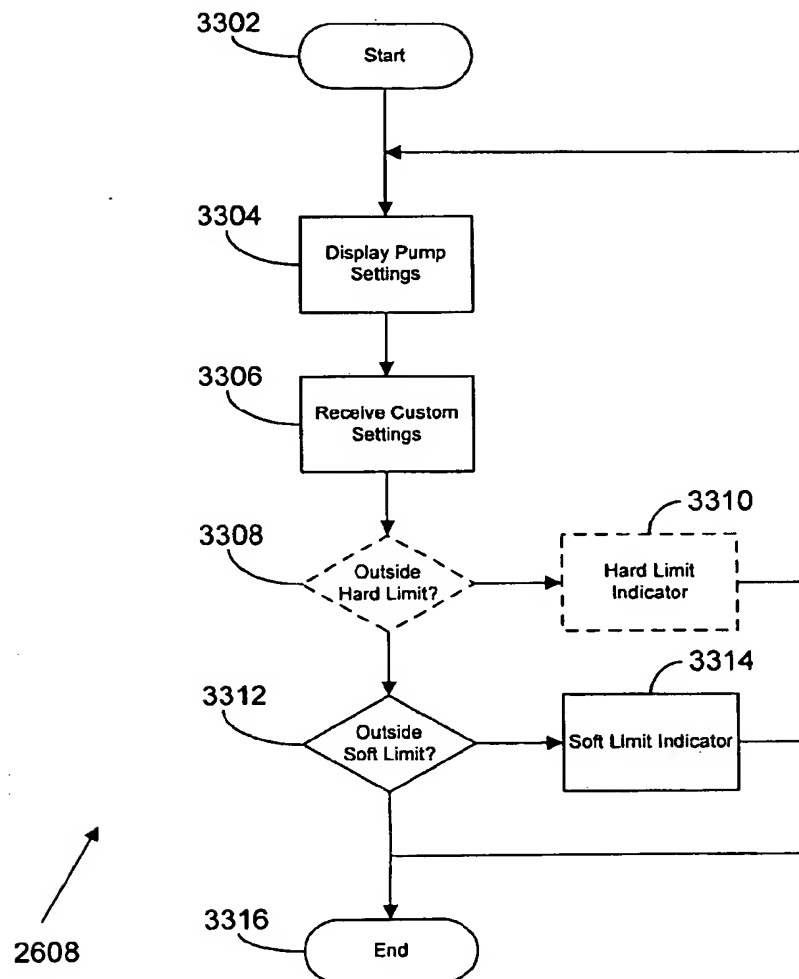

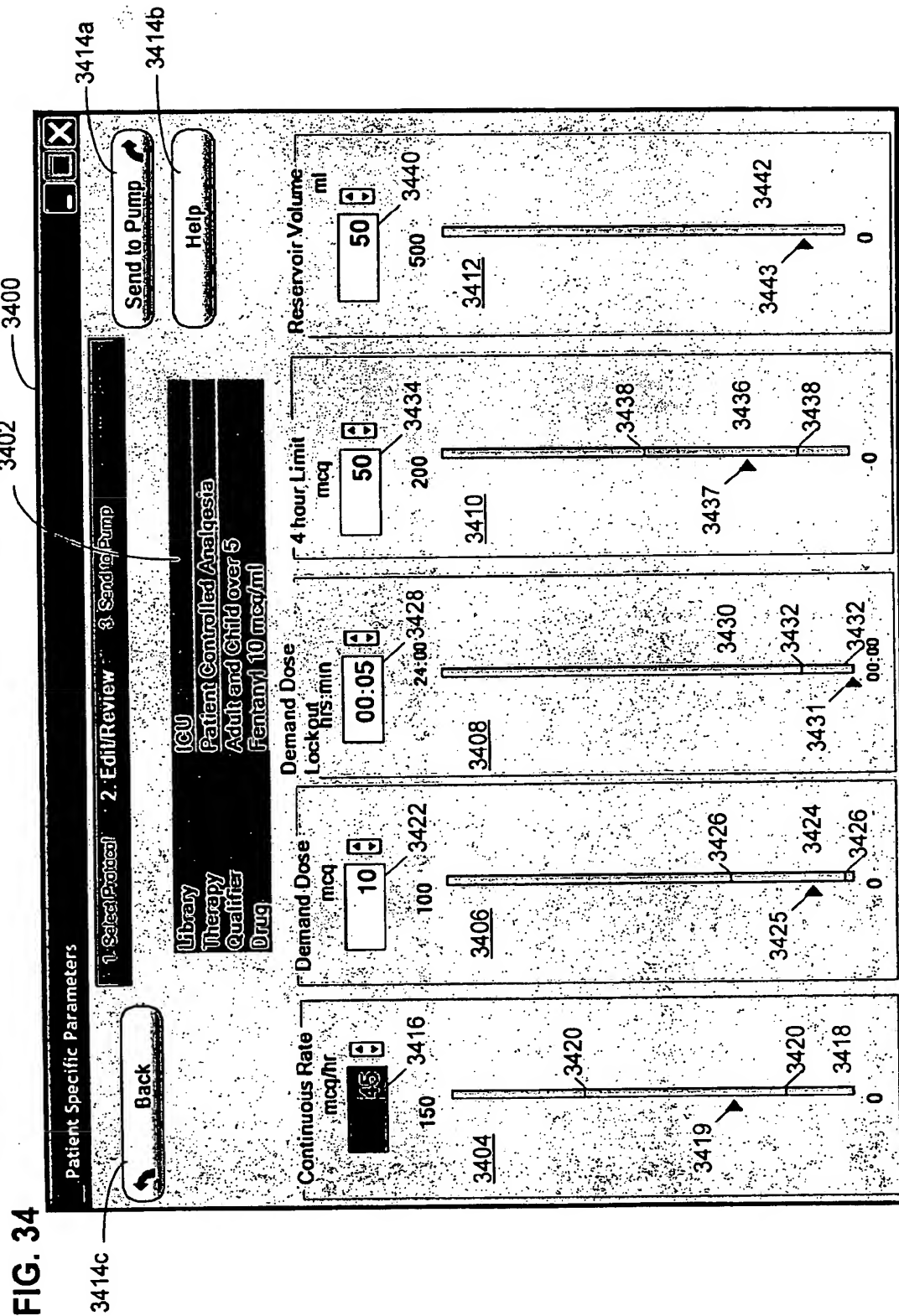

35/41

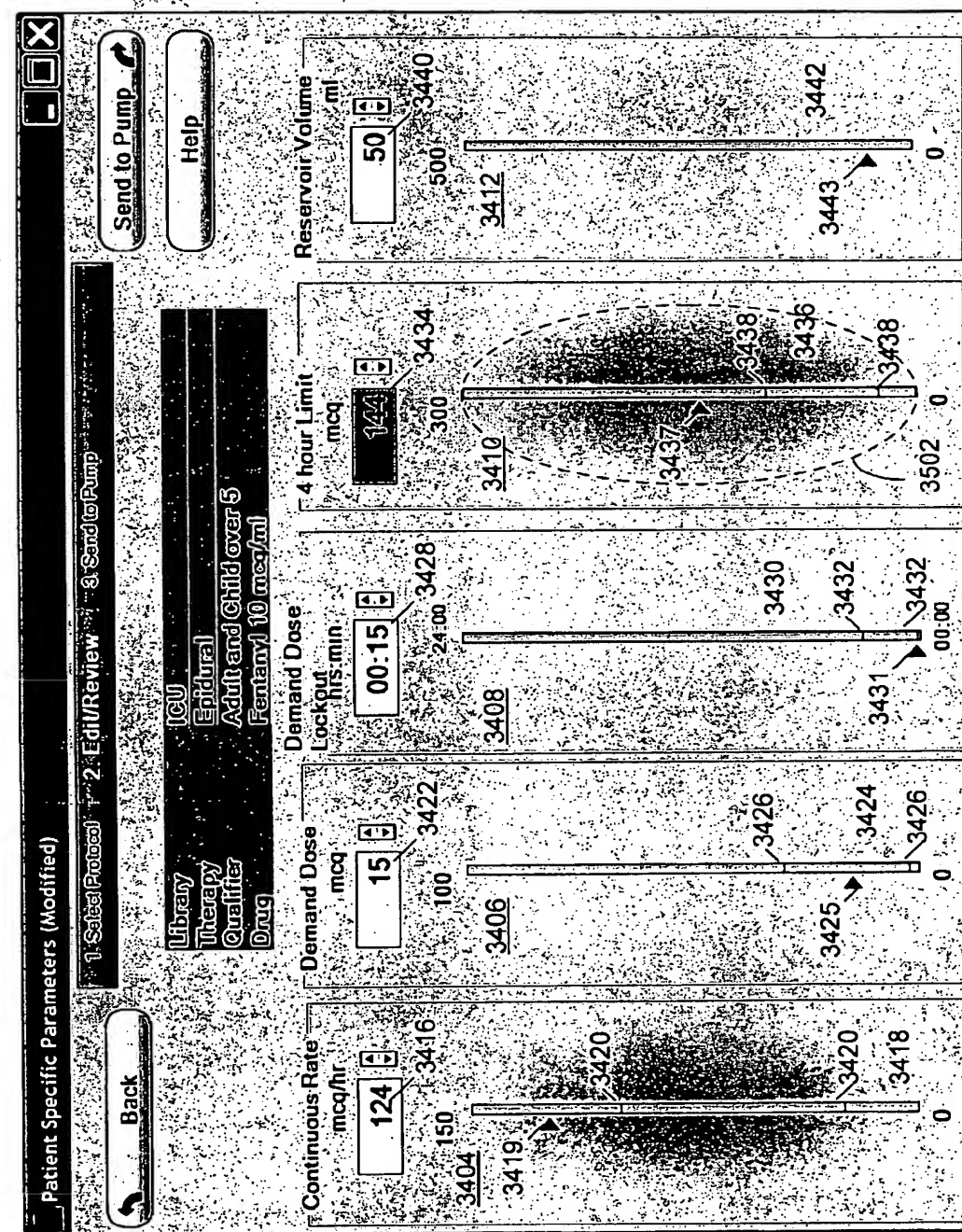

Send to Pump

Back

1. Select Protocol2. Edit/Review3. Send to Pump

LibraryTherapyQualifierDrug

ICUPatient Controlled AnalgesiaAdult and Child over 5Fentanyl 10 mcg/ml

Continuous Rate45 mcg/hr3602a

DemandDose0.5 mcg3602b

4 hour Limit50 hrs min3602c

Reservoir Volume50 ml3602e

Not ConnectedLibrary ICU

CloseHelp

Send

Please connect a pump. If a pump is connected you may need to change the COM setting.

FIG. 36

37/41

FIG. 37

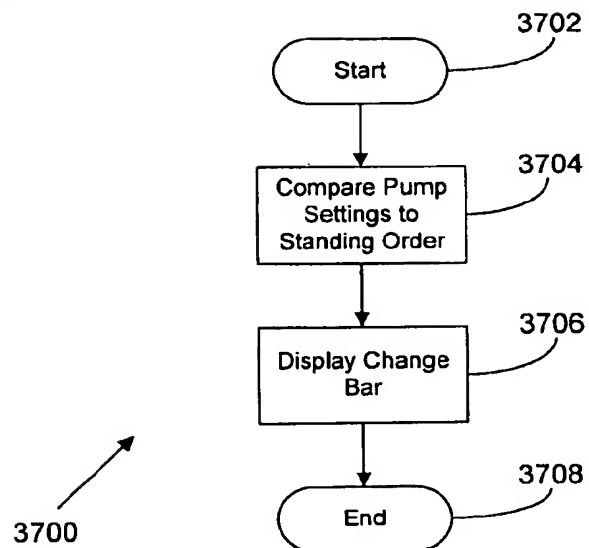
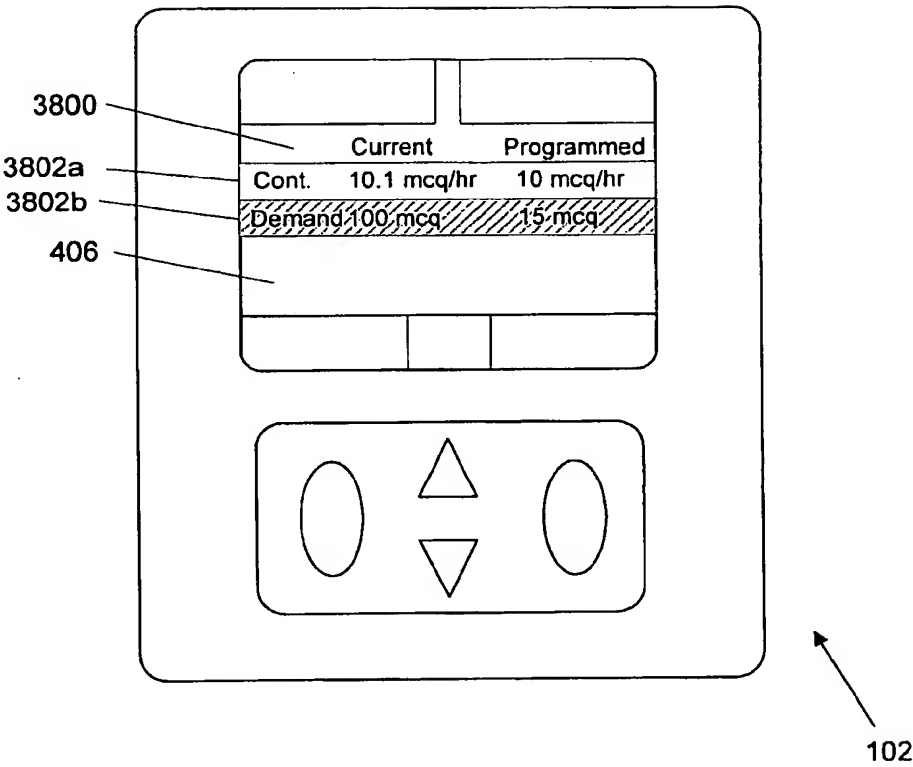


FIG. 38




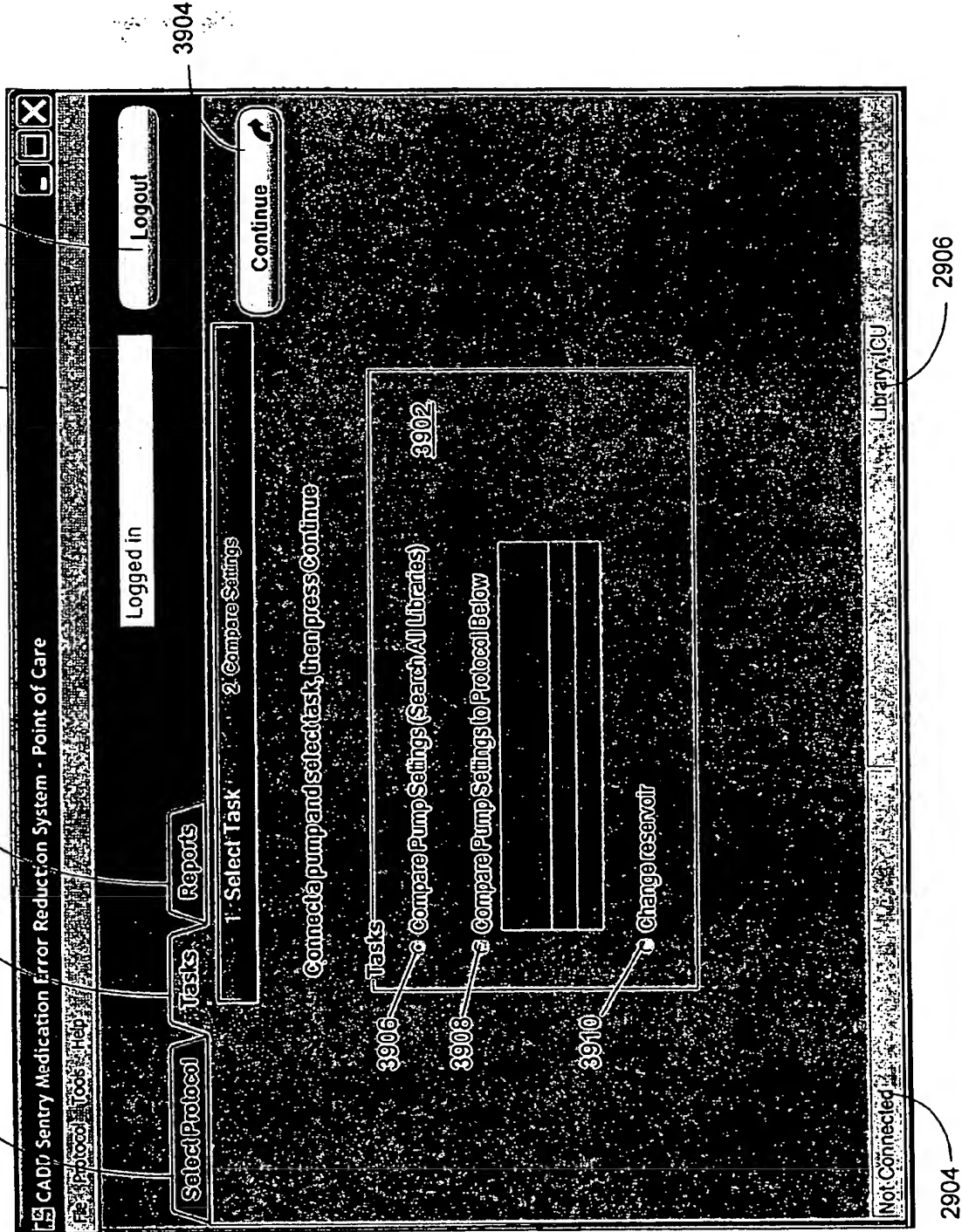


FIG. 39

FIG. 40

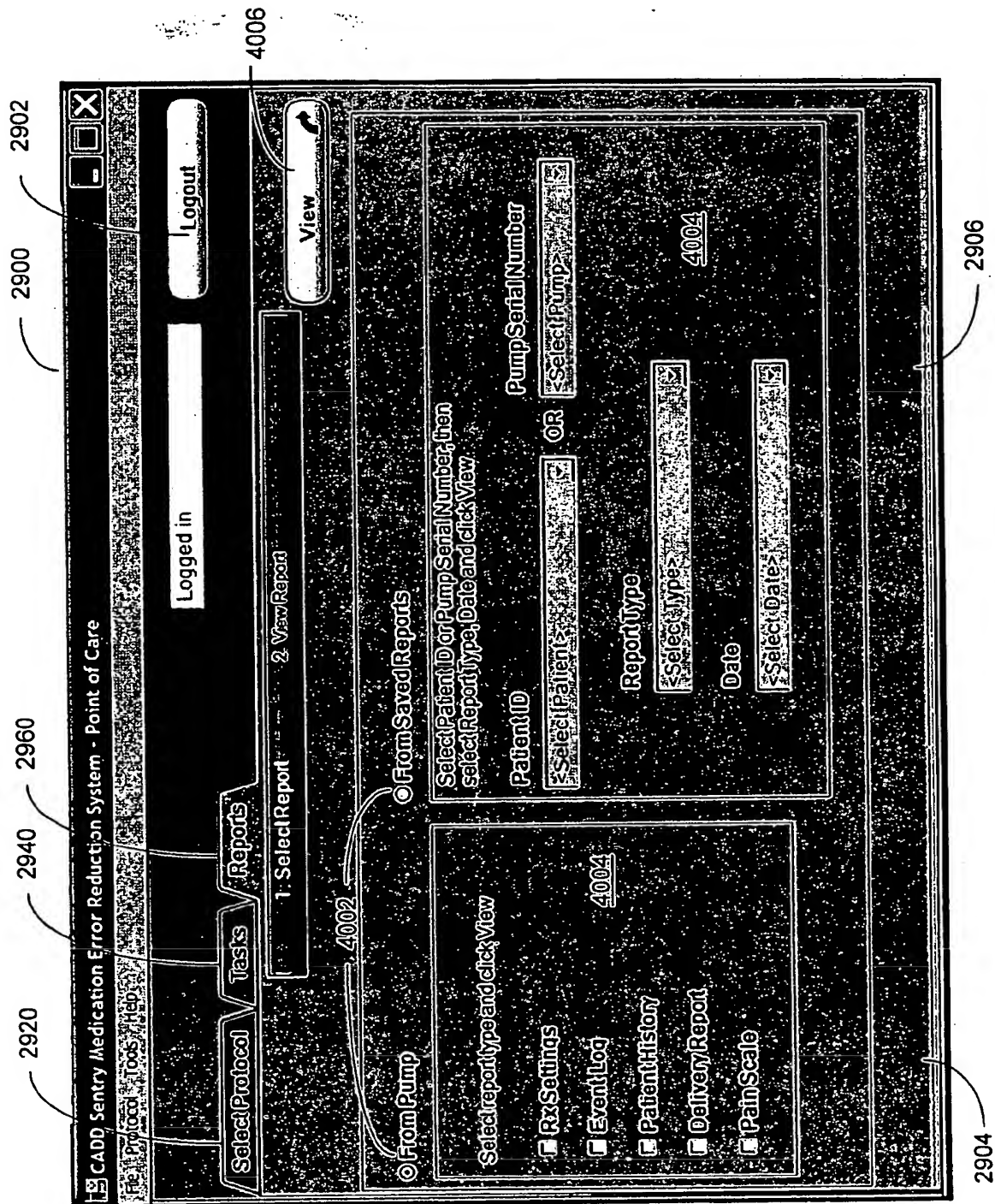

41/41

4100

FIG. 41

Report Viewer

1. Select Report 2. View Report

Back Print Save

001> <<<<<DELIVERY REPORT>>>>>

002>

003> Patient ID: _____

004>

005> Patient's Name: _____

006>

007> CADP-Prizm+ (TM) Model 6101 SN 00611188

008> Date Printed: 01/12/05 Time: 11:41

009> Delivery Mode: PCS 6215-

010>

011> Drug: _____

012>

013>

014>

015>

Units	Current Rx SETTINGS
016> Milligrams	100.0 mg/ml
017> Concentration	0.00 mg/hr
018> Continuous Rate	0.00 mg
019> Demand Dose	1.0 ml
020> Reservoir Volume	0.00 mg since 01/06/05 at 09:34
021> Milligrams Given	LL2
022> Lock Level	On
023> Upstream Sensor Enable	
024>	
025>	
026>	

DELIVERY SINCE NEW PATIENT